

## Oasis Sample Extraction Products: Enhanced Sensitivity Using $\mu$ Elution Plate

Waters Corporation



This is an Application Brief and does not contain a detailed Experimental section.

---

## Abstract

This application brief describes about Oasis  $\mu$ Elution plates in combination with Oasis sorbents that are an ideal solution for small sample volumes up to 375 mL, providing up to a 15X increase in concentration.

---

## Introduction

Waters Oasis  $\mu$ Elution plates combine patented plate design and proven Oasis chemistries to provide high analyte recovery and clean extracts in elution volumes as low as 25  $\mu$ L. This can enhance a bioanalytical assay's sensitivity and allow scientists to easily achieve their required LLOQs.

- Ideal for small sample volumes (10-25  $\mu$ L)
- Eliminates evaporation and reconstitution
- Offers up to a 15X increase in concentration without evaporation
- Excellent recovery for a diverse range of therapeutic peptides


---

## Results and Discussion

### Superior Sensitivity with $\mu$ Elution format

The Oasis  $\mu$ Elution format offers sample enrichment of up to 15X, providing superior sensitivity to bioanalytical assays.

### Risperidone response: 15X increase in concentration

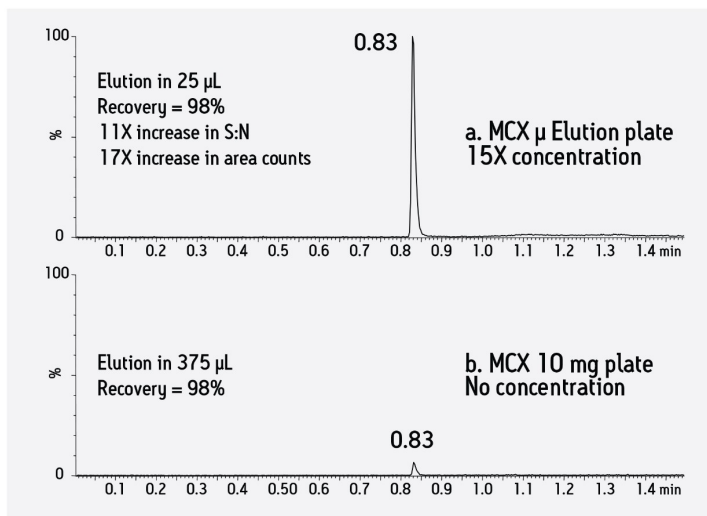


*375  $\mu$ L of 0.5 ng/mL risperidone in human plasma sample extracted with Oasis MCX  $\mu$ Elution plate a) MCX  $\mu$ Elution and b) MCX 10 mg.*

### Highest SPE Recoveries

Oasis sorbents are water-wettable which allows them to maintain high retention and capacity for a variety of analytes even if the well or cartridge runs dry. With Oasis sorbents, the analytical scientist can be assured of reliable, robust performance and the highest degree of recovery possible.

The efficient therapeutic peptide solid-phase extraction (SPE) screening strategy, based on a single protocol and two of the Oasis mixed-mode sorbents in  $\mu$ Elution format, simplifies method development for the extraction of peptides from human plasma and provide excellent recoveries for diverse set of peptides.

## Excellent SPE Recovery for 12 Peptides

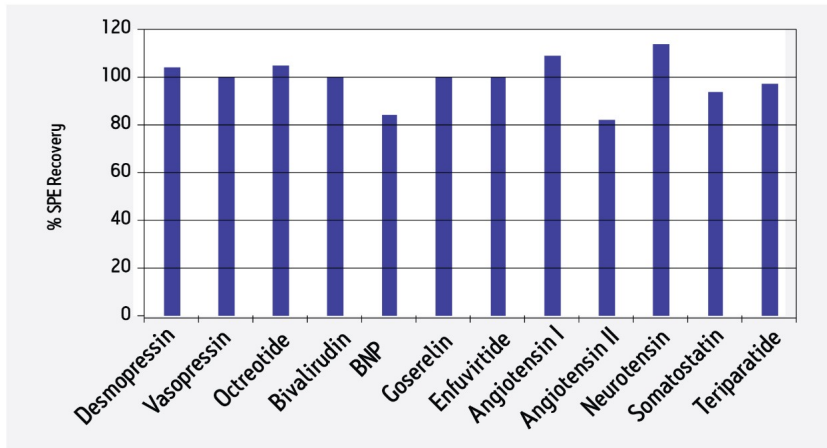


Excellent recovery for 12 diverse peptides using Oasis Mixed-Mode SPE in  $\mu$ Elution plate format.

## Lowest Matrix Effects

The combination of Oasis mixed-mode SPE and the  $\mu$ Elution format provide the cleanest SPE extract possible, facilitating the meeting of regulatory bioanalytical requirements.

## Matrix Effects for 1 ng/mL of Imipramine, Nonafluoropentanoic Acid, Ibuprofen, and Valethamate in Human Plasma

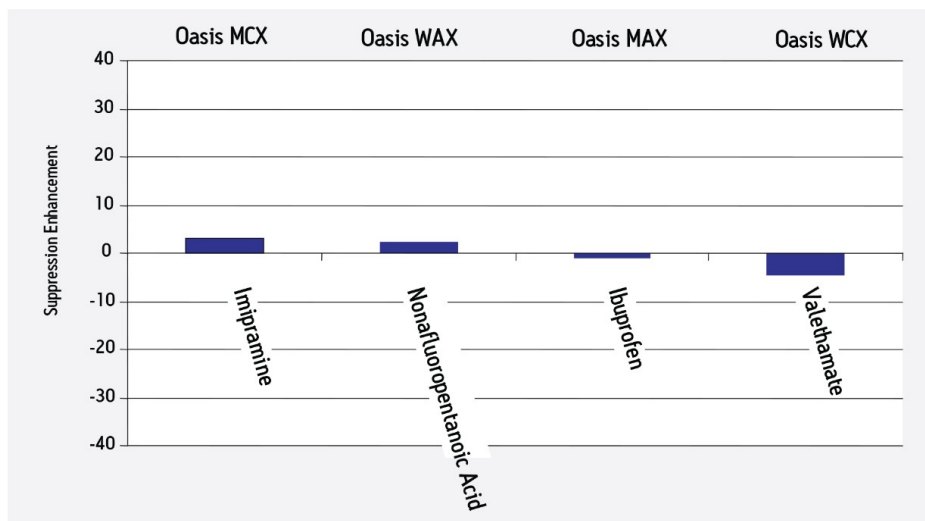


Matrix effect results of 4 diverse small molecules on 4 mixed-mode Oasis chemistries in (1 ng/mL in human plasma).

## $\mu$ Elution Plate Loading Capacity


The physicochemical properties of Oasis sorbents are designed to provide exceptionally high loading capacity, even though each well in a Waters Oasis  $\mu$ Elution plate contains only 2 mg of Oasis sorbent.

To determine the Oasis  $\mu$ Elution plate capacity, increasing volumes of plasma and urine samples (from 50  $\mu$ L to 350  $\mu$ L in 50  $\mu$ L increments) were spiked with 200 ng/mL Imipramine (non-polar base) and 200 ng/mL Atenolol (polar base). The plasma aliquots were diluted 1:1 with 4% aqueous  $H_3PO_4$  and the urine aliquots were diluted 1:1 with  $H_2O$  and then loaded onto the  $\mu$ Elution plate. SPE recovery was calculated and plotted for each loading level.

### SPE Recovery for 200 ng/mL Imipramine and 200 ng/mL Atenolol on Oasis MCX $\mu$ Elution Plate

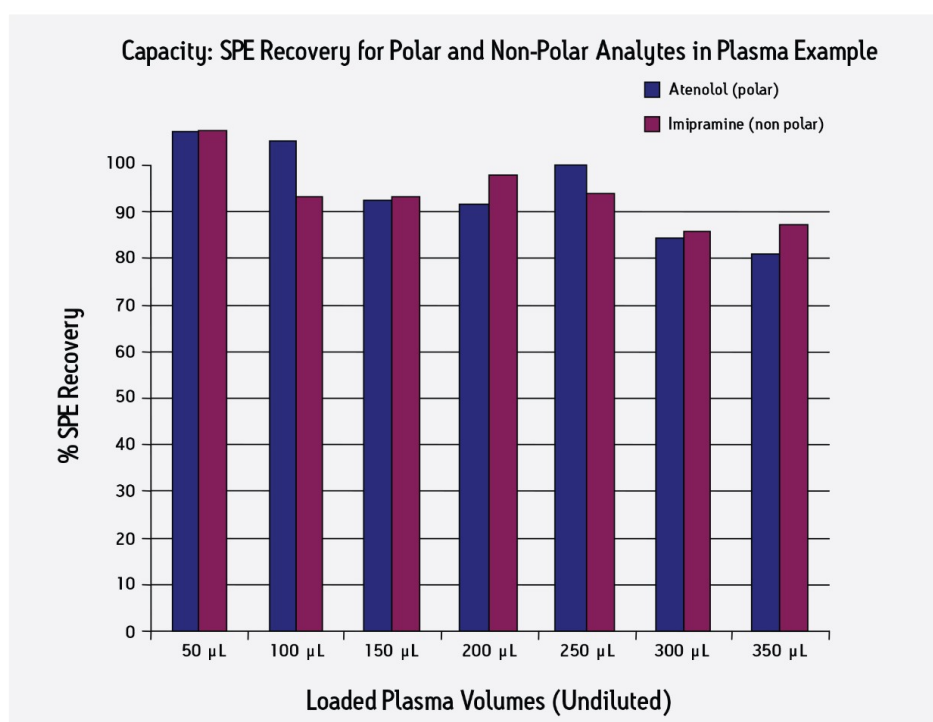

### Capacity: SPE Recovery for Polar and Non-Polar Analytes in Urine Example

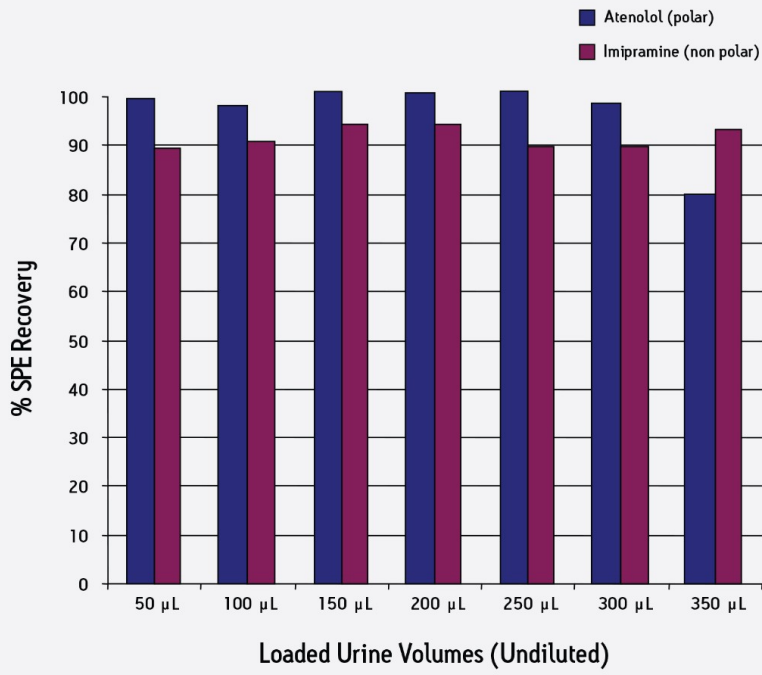

---

## Conclusion

Using a unique plate design in combination with Oasis sorbents, the Oasis  $\mu$ Elution plates are an ideal solution for small sample volumes up to 375 mL, providing up to a 15X increase in concentration. This concentration, along with excellent recoveries, equates to unparalleled sensitivity for bioanalytical samples.

720003680, September 2010

©2019 Waters Corporation. All Rights Reserved.