



Pheniramine

Waters Corporation




This is an Application Brief and does not contain a detailed Experimental section.

Abstract

This application brief highlights the analysis of pheniramine using XTerra RP₁₈ columns.

Introduction

Pheniramine has been analyzed in this study.

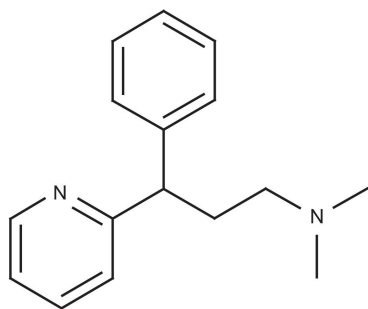

Pheniramine

Experimental

HPLC Conditions

Column:	XTerra RP ₁₈ 4.6 x 150 mm, 5 µm (p/n: 186000492)
Mobile phase:	At pH 3.0: H ₂ O/ACN/100 mM, NH ₄ COOH, pH 3.0, 55:35:10 At pH 7.0: H ₂ O/ACN/100 mM, NH ₄ HCO ₃ , pH 7.0, 30:60:10 At pH 10.0: H ₂ O/ACN/100 mM, NH ₄ HCO ₃ , pH 10.0, 60:30:10
Flow rate:	1.0 mL/min

Injection volume: 5 μ L of 250 μ g/mL


Temperature: 30 $^{\circ}$ C

Detection: UV @ 265 nm

Instrument: Alliance 2695, 2996 PDA

Mobile phase pH	USP Tailing
3.0	1.19
7.0	1.30
10.0	1.23

Results and Discussion

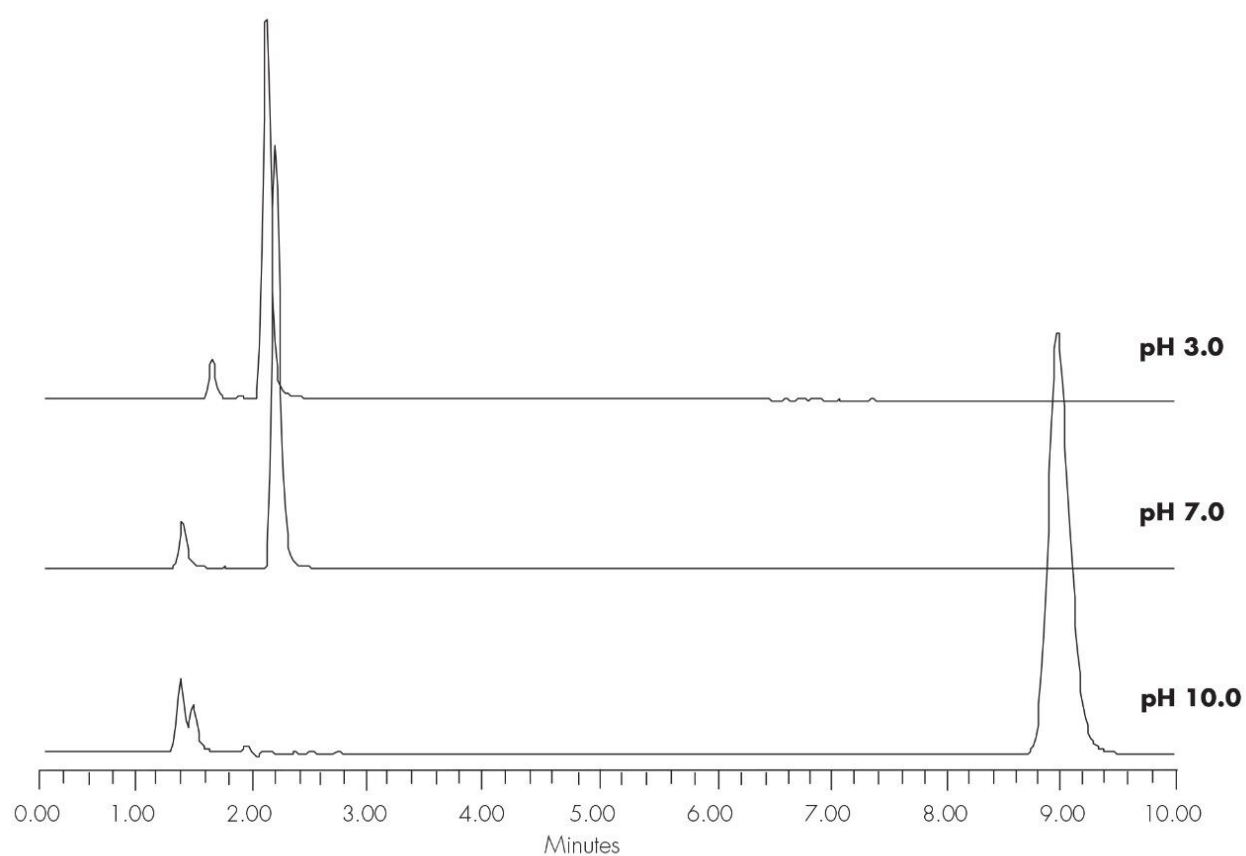

Featured Products

Alliance HPLC <<https://www.waters.com/514248>>

WA20738.087, June 2002

