

## QC of Synthetic Oligonucleotides

---

Waters Corporation



This is an Application Brief and does not contain a detailed Experimental section.

---

### Abstract

This application brief demonstrates QC of synthetic oligonucleotides.

---

# Introduction

The compound analyzed in this study is synthetic oligonucleotides.

---

## Experimental

### Conditions

Column:	Xterra MS C <sub>18</sub> 4.6 x 50 mm, 2.5 µm
Part number:	186000602
Mobile phase A:	5% ACN, 95% 100 mM Triethylammonium acetate, pH 7.0
Mobile phase B:	20% ACN, 80% 100 mM Triethylammonium acetate, pH 7.0
Flow rate:	1.0 mL/min
Injection volume:	0.6 µL
Temperature:	60 °C
Detection:	UV @ 260 nm
Instrument:	Alliance 2695, 2996 PDA


### Gradient

Time(min)	Profile
-----------	---------

	%A	%B
0.0	90.7	9.3
20.8	56.0	44.0

---

## Results and Discussion

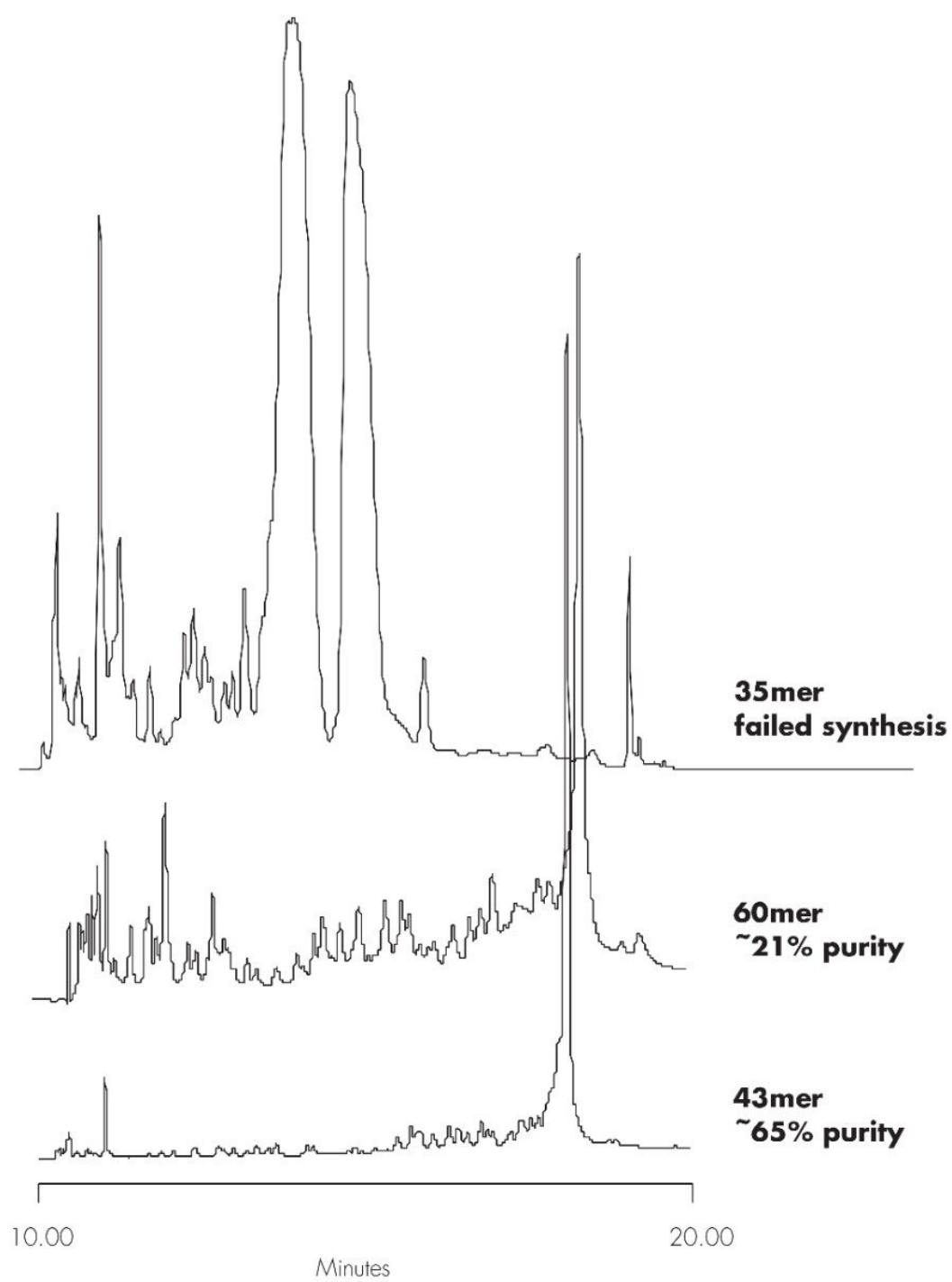

---

Featured Products

Alliance HPLC System <<https://www.waters.com/534293>>

WA20738.095, June 2002

© 2021 Waters Corporation. All Rights Reserved.