

Salbutamol in Rat Plasma

Waters Corporation




This is an Application Brief and does not contain a detailed Experimental section.

Abstract

This application brief demonstrates analysis of salbutamol in rat plasma.

Introduction

The compound analyzed in this study is salbutamol in rat plasma.

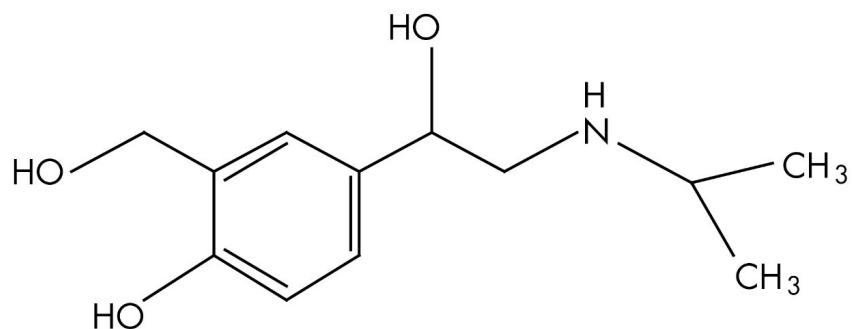

Salbutamol


Experimental

Conditions

Column:	Xterra MS C ₁₈ 2.1 x 30 mm, 3.5 μm
Part number:	186000398
Mobile phase A:	100 mM NH ₄ COOH
Mobile phase B:	ACN
Isocratic mobile phase composition:	70% A; 30% B
Flow rate:	0.2 mL/min
Injection volume:	50 μL

Detection:	MS ESI+
Instrument:	Alliance 2790, Micromass Quattro Ultima
Ion source:	ESI+
Source temperature:	150 °C
Gas cell:	1.5e-3 mbar, 15eV
Desolvation temperature:	350 °C
Cone gas flow:	150 L/hr
Drying gas flow:	600 L/hr
Cone voltage:	20V

Oasis® MCX Extraction Method
Oasis® MCX Extraction Plate, 10 mg/96-well
Part Number 186000259

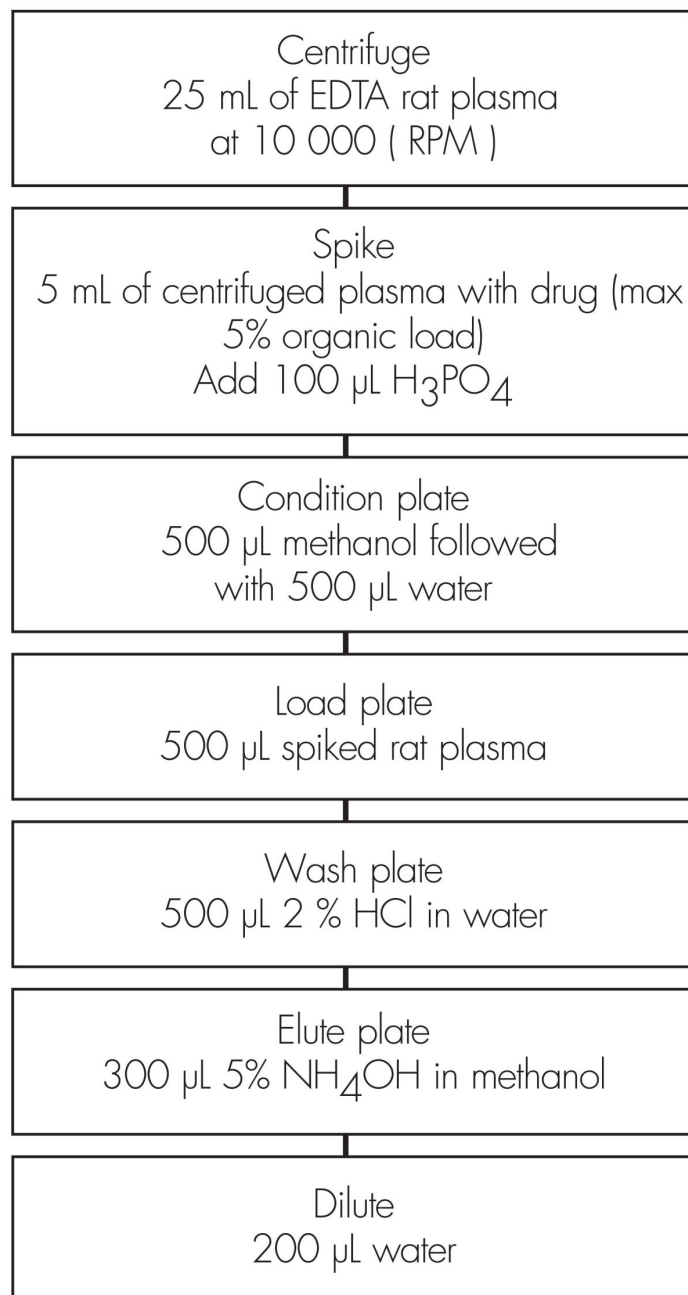


Results and Discussion

CID mass spectra

esi cid 15
3 1 (1.005)

Daughters of 240ES+
2.35e9

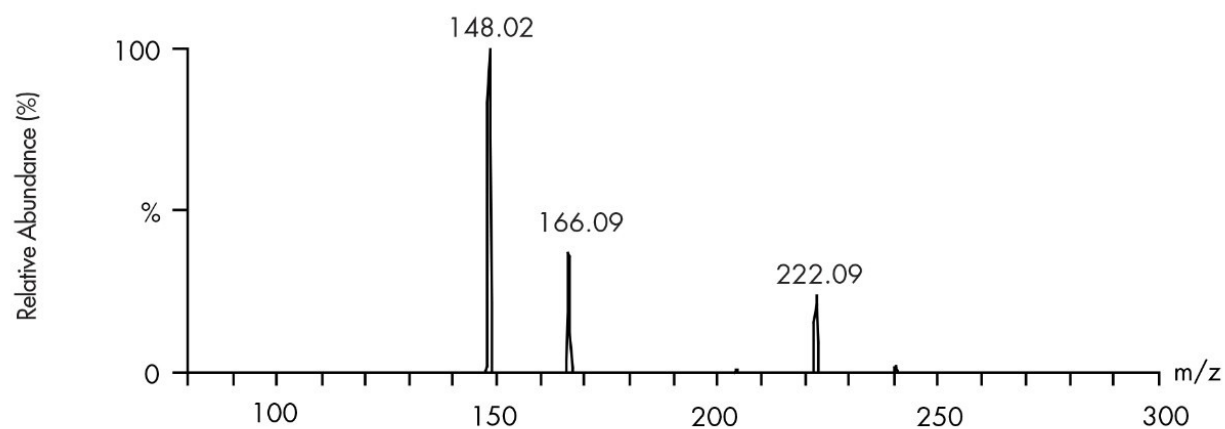

Featured Products

· [Alliance HPLC System <https://www.waters.com/534293>](https://www.waters.com/534293)

WA20738.101, June 2002



© 2021 Waters Corporation. All Rights Reserved.