

Note d'application

## Promethazine – Impurities Profile

---

Waters Corporation



This is an Application Brief and does not contain a detailed Experimental section.

---


### Abstract

This application brief highlights the impurities profile of Promethazine using SymmetryShield Columns.

---

## Introduction

Promethazine HCL impurities are studied in this application brief.

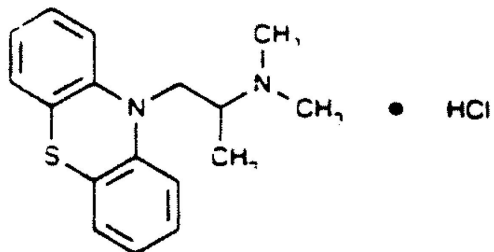

### 1. Promethazine HCL

---

## Experimental

### HPLC Method


Column:	SymmetryShield RP <sub>18</sub> , 3.0 x 150 mm, 5 μm (p/n: 186000692)
Mobile phase A:	20 mM potassium phosphate, pH 7.0
Mobile phase B:	Acetonitrile
Flow rate:	1.0 mL/min
Injection volume:	10 μL of 0.5 mg/mL
Detection:	UV @ 254 nm

Gradient

Time (min)	Profile	
	%A	%B
0	95	5
30	20	80

---

Results and Discussion

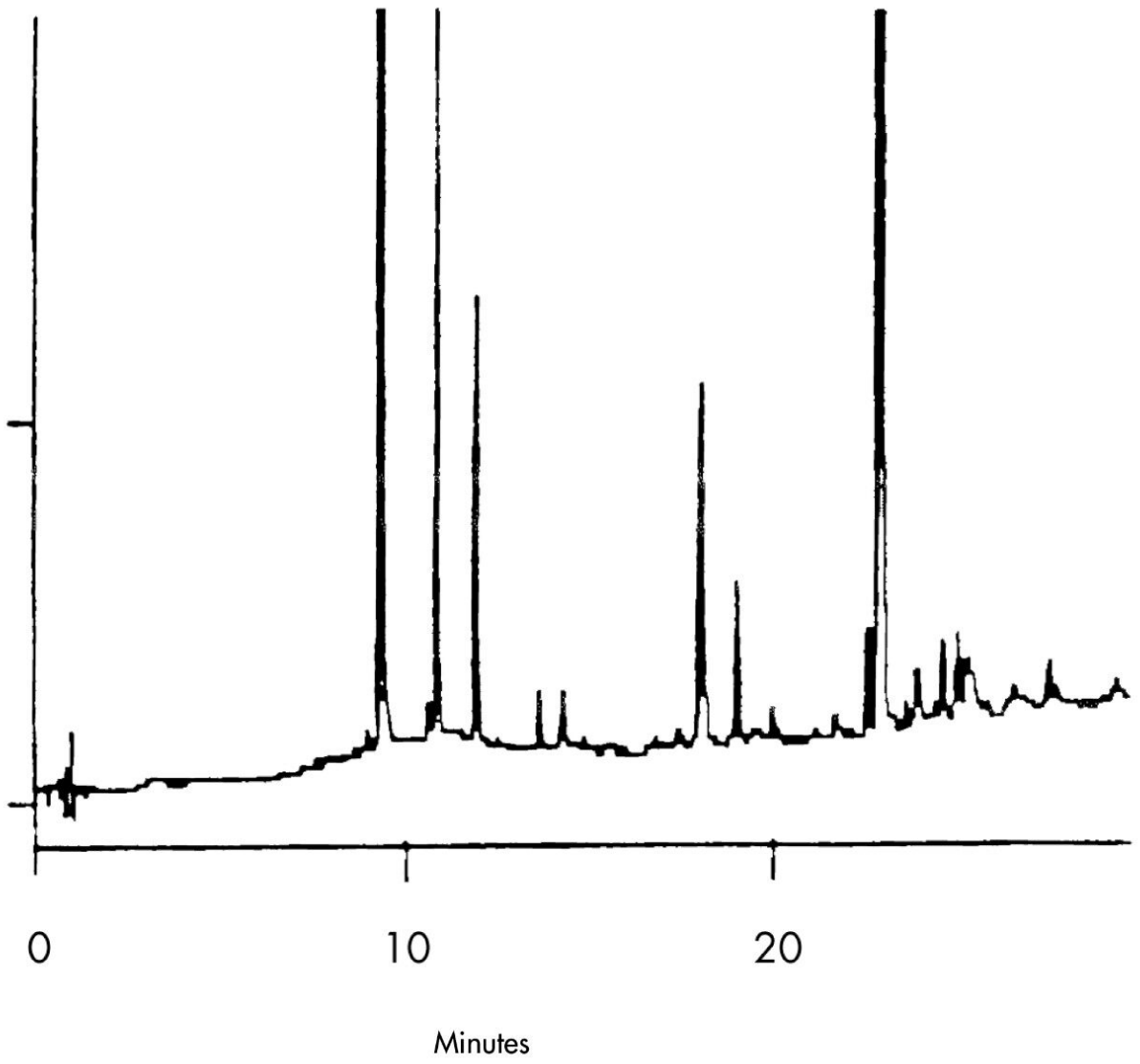

---

## Featured Products

WA31763.133, June 2003

© 2021 Waters Corporation. All Rights Reserved.