

Nalidixic Acid Antibiotics by LC-MS, 2.0 Minute Gradient – 2.1 x 20 mm Intelligent Speed Separation

Waters Corporation




This is an Application Brief and does not contain a detailed Experimental section.

Abstract


This application brief details on the analysis of nalidixic acid antibiotics by LC-MS.

Introduction


| Compound | MW |
|-------------------|-------|
| 1. Cinoxacin | 262.2 |
| 2. Oxolinic Acid | 261.2 |
| 3. Nalidixic Acid | 232.2 |

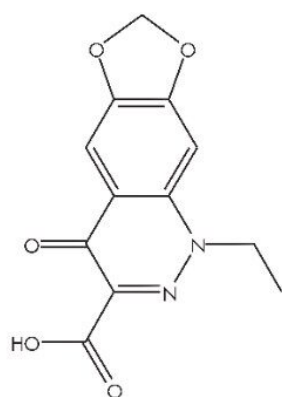

Cinoxacin

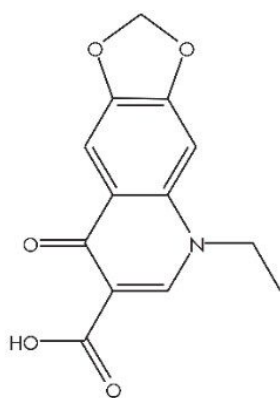

Oxolinic acid

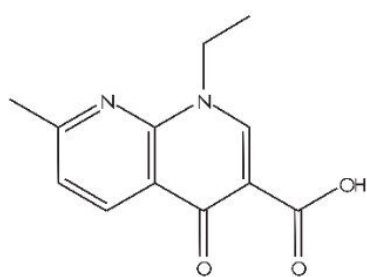

Nalidixic acid

Experimental

Conditions

| | |
|-----------------------|--|
| Column: | Atlantis dC ₁₈ , 2.1 x 20 mm IS, 3.0 µm, (P/N: 186002058) |
| Mobile phase A: | Water |
| Mobile phase B: | Methanol |
| Mobile phase C: | 1% HCOOH in Water |
| Flow rate: | 0.4 mL/min |
| Injection volume: | 2 µL |
| Sample concentration: | 10 µg/mL |
| Temperature: | 30°C |
| Instrument: | Alliance 2795 and Waters ZQ |

Gradient

| Time (min) | Profile | | |
|------------|---------|----|----|
| | %A | %B | %C |
| 0.0 | 60 | 30 | 10 |
| 2.0 | 40 | 50 | 10 |

MS Conditions

Waters ZQ


ES+

| | |
|-----------------------------|-----|
| Capillary (kV): | 3.5 |
| Cone (V): | 5.0 |
| Extractor: | 3.0 |
| RF lens: | 0.1 |
| Source temp. (°C): | 150 |
| Desolvation temp. (°C): | 400 |
| Cone gas flow (L/Hr): | 50 |
| Desolvation gas fLow(L/Hr): | 500 |
| LM resolution: | 15 |
| HM resolution: | 15 |
| Ion energy: | 0.5 |
| Multiplier (V): | 650 |

Results and Discussion

The top figure is the total ion current, followed by the extracted ion signals for each of the three analytes.

*The "extra" peak in the cinoxacin panel is the isotope from oxolinic acid.

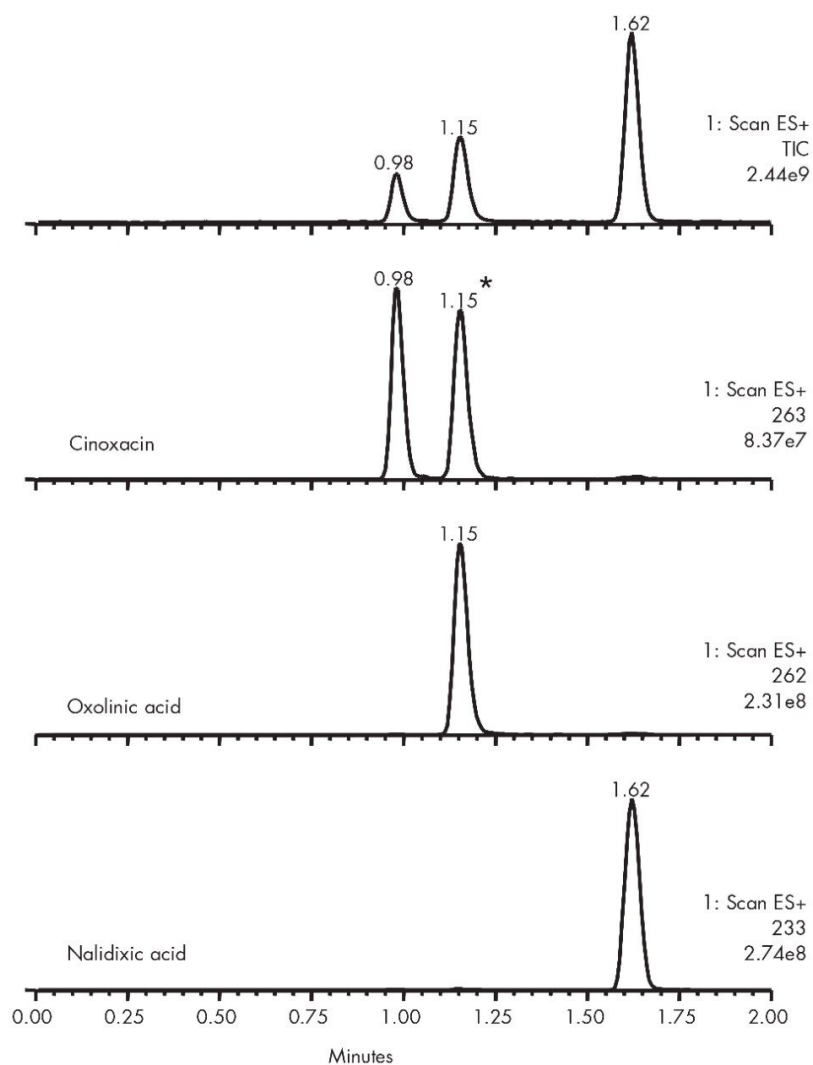

Featured Products

Alliance HPLC System <<https://www.waters.com/534293>>

WA31787.17, June 2003

