


Isocratic Analysis of Acyclovir on Atlantis T3 Columns

Waters Corporation

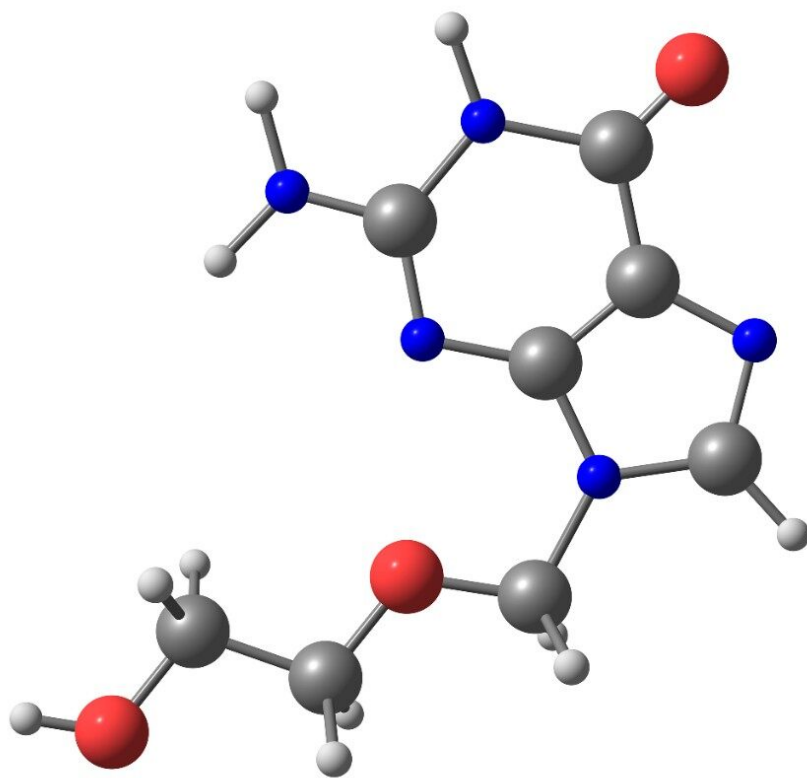


This is an Application Brief and does not contain a detailed Experimental section.

Abstract

This application brief demonstrates isocratic analysis of acyclovir on Atlantis T3 Columns.

Introduction

The compound used in this study is Acyclovir.

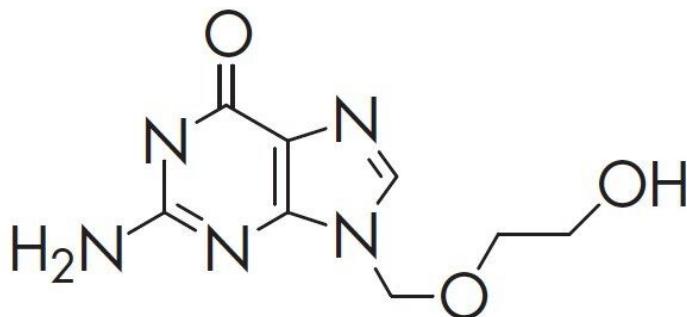

Acyclovir


Experimental

Test Conditions

Column:	Atlantis T3, 4.6 x 50 mm, 3 μ m
Part number:	186003726
Mobile phase A:	H ₂ O
Mobile phase B:	ACN
Mobile phase C:	100 mM HCOONH ₄ , pH 3.0
Flow rate:	1.0 mL/min
Isocratic mobile phase composition:	87% A; 3% B; 10% C

Injection volume:	5 μ L
Sample:	Acyclovir (10 μ g/mL) in H ₂ O
Column temperature:	30 $^{\circ}$ C
Sample temperature:	20 $^{\circ}$ C
Detection:	UV @ 254 nm
Sampling rate:	5 point/sec
Filter response:	0.2
Instrument:	Waters Alliance 2695 with 2996 PDA

Results and Discussion

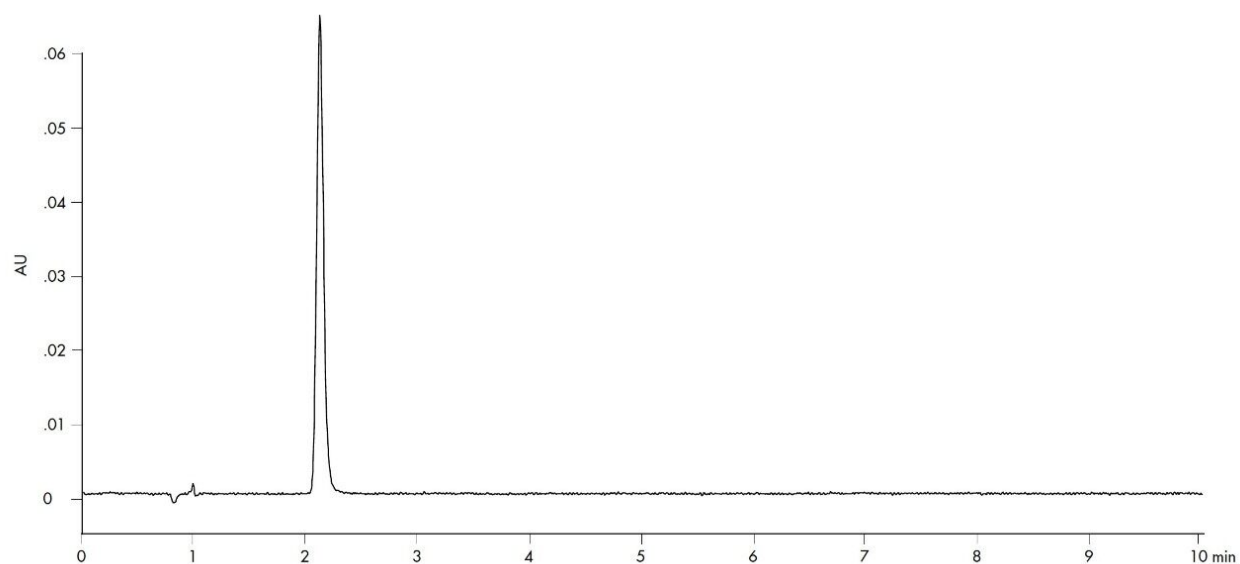

Featured Products

Alliance HPLC System <<https://www.waters.com/534293>>

WA60213, January 2009

© 2021 Waters Corporation. All Rights Reserved.