

Analysis of Aromatic Amino Acids Using Atlantis T3

Waters Corporation



This is an Application Brief and does not contain a detailed Experimental section.

Abstract

This application brief highlights the analysis of aromatic amino acids.

Introduction

Compounds analysed in this application brief :

1. Tyrosine
2. Phenylalanine
3. Tryptophan

Experimental

Test Conditions


Column:	Atlantis T3, 4.6 x 50 mm, 3 μ m
Part Number:	186003726
Mobile Phase A:	H ₂ O
Mobile Phase B:	MeOH
Mobile Phase C:	100 mM HCOONH ₄ , pH 3.0
Flow Rate:	1.0 mL/min
Injection Volume:	10 μ L
Sample:	Tyrosine (150 μ g/ mL), Phenylalanine (200 μ g/mL), Tryptophan (30 μ g/mL) in H ₂ O

Column Temperature:	30 °C
Sample Temperature:	20 °C
Detection:	UV @ 254 nm
Sampling Rate:	5 point/sec
Filter Response:	0.2
Instrument:	Waters Alliance 2695 with 2996 PDA

Gradient

Time (min)	Profile
	%A
0.00	90
10.00	10
11.00	90
15.00	90

Results and Discussion

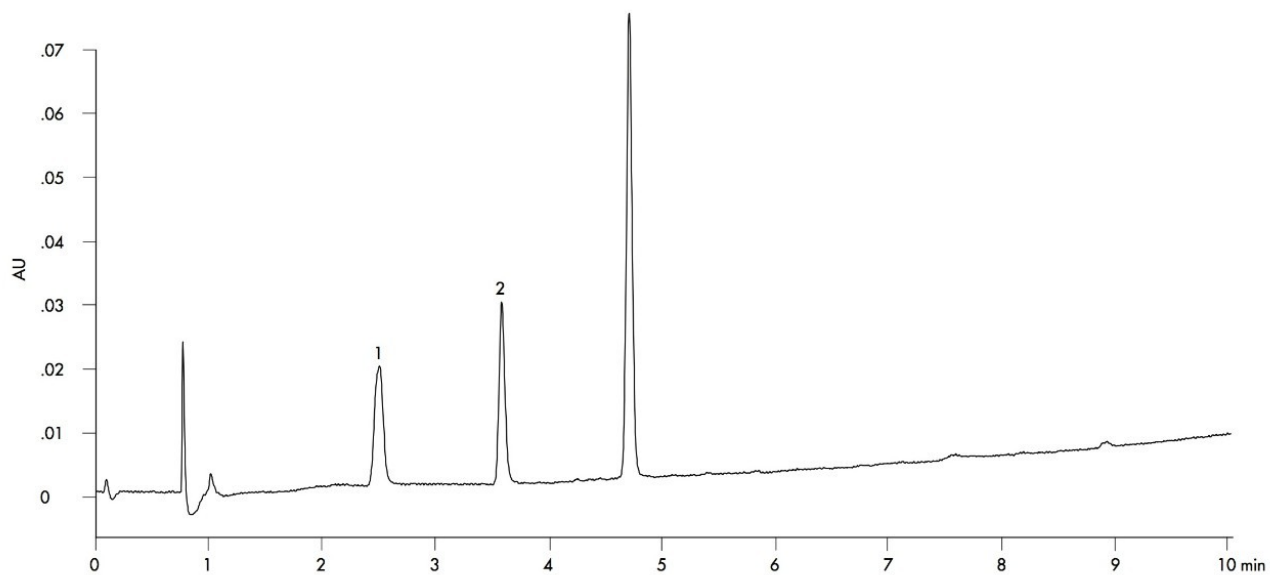

1. Tyrosine, 2. Phenylalanine, 3. Tryptophan

Featured Products

- [Alliance HPLC System <https://www.waters.com/534293>](https://www.waters.com/534293)
- [2998 Photodiode Array \(PDA\) Detector <https://www.waters.com/1001362>](https://www.waters.com/1001362)

WA60221, January 2009



©2019 Waters Corporation. All Rights Reserved.