

응용 자료

## Multi-Residue Analysis of Pesticides in Grapes Using AOAC QuEChERS Method by UPLC-MS/MS

---

Waters Corporation



This is an Application Brief and does not contain a detailed Experimental section.

---

## Abstract

This application brief demonstrates multi-residue analysis of pesticides in grapes using AOAC QuEChERS method by UPLC-MS/MS.

---

## Experimental

### Test Conditions

#### LC Conditions

LC System:	Waters ACQUITY UPLC System
Column:	ACQUITY UPLC BEH C <sub>18</sub> , 2.1 x 100 mm, 1.7 µm
Column Temp:	40 °C
Sample Temp:	4 °C
Flow Rate:	0.3 mL/min.
Mobile Phase A:	Water + 0.1% formic acid
Mobile Phase B:	Methanol + 0.1% formic acid
Injection Volume:	15 µL, Partial loop injection

### Gradient:


Time	Flow Rate	A%	B%
0.00	0.3	75	25
0.25	0.3	75	25
7.75	0.3	5	100
8.50	0.3	0	100
8.51	0.5	75	25
10.50	0.5	75	25
11.0	0.3	75	25

## MS Conditions

Instrument:	Waters ACQUITY TQ Detector
Ionization:	Positive electrospray (ESI+)
Acquisition:	Multiple reaction monitoring (MRM)

---

## Results and Discussion

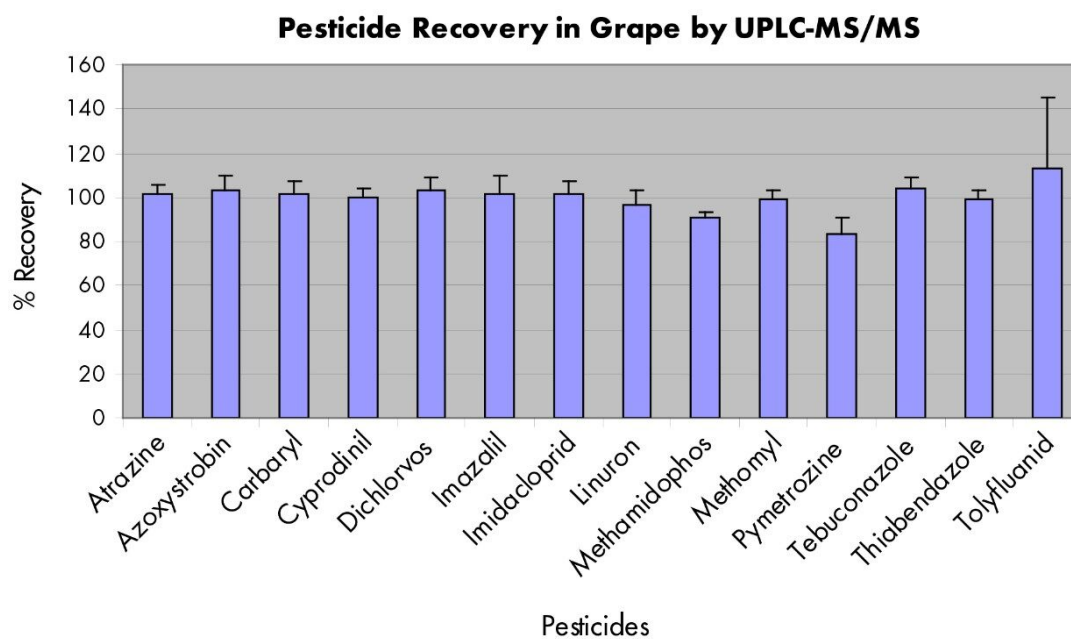

---

## Featured Products

ACQUITY UPLC System <<https://www.waters.com/514207>>

720003369, March 2010