



## Multi-Residue Analysis of Pesticides in Avocado Using AOAC QuEChERS Method by UPLC-MS/MS

---

Waters Corporation



This is an Application Brief and does not contain a detailed Experimental section.

---

## Abstract

This application brief demonstrates multi-residue analysis of pesticides in avocado using AOAC QuEChERS method by UPLC-MS/MS.

---

## Experimental

### Extraction Procedure

1. Add 15 mL 1% acetic acid in acetonitrile into the 50 mL DisQuE extraction tube.
2. Add 15 g of homogenized sample into the 50 mL tube.
3. Add any internal standards and standard mixture.
4. Shake vigorously for 1 minute and centrifuge > 1500 rcf for 5 minute.
5. Transfer 1 mL of the acetonitrile extract into the 2 mL clean-up tube containing 50 mg PSA, 150 mg MgSO<sub>4</sub>, and 50 mg C<sub>18</sub>.
6. Shake for 30 seconds and centrifuge >1500 rcf for 1 minute.
7. Transfer 100 µL of final extract into a 1.5 mL centrifuge tube.
8. Add any post-extraction internal standards.
9. Dilute as needed with an appropriate buffer or solvent.
10. Centrifuge > 16000 rcf for 5 minutes.
11. Transfer to autosampler vial.

### Test Conditions

#### LC Conditions

LC System: Waters ACQUITY UPLC

## LC Conditions

### System

Column:	ACQUITY UPLC BEH C <sub>18</sub> , 2.1 x 100 mm, 1.7 µm
Column Temp:	40 °C
Sample Temp:	4 °C
Flow Rate:	0.3 mL/min.
Mobile Phase A:	Water + 0.1% formic acid
Mobile Phase B:	Methanol + 0.1% formic acid
Injection Volume:	15 µL, Partial loop injection

### Gradient:

Time	Flow Rate	A%	B%
0.00	0.3	75	25
0.25	0.3	75	25
7.75	0.3	5	100
8.50	0.3	0	100
8.51	0.5	75	25
10.50	0.5	75	25

Time	Flow Rate	A%	B%
11.0	0.3	75	25


## MS Conditions

Instrument: Waters ACQUITY TQ Detector

Ionization: Positive electrospray (ESI+)

Acquisition: Multiple reaction monitoring  
(MRM)

## Results and Discussion

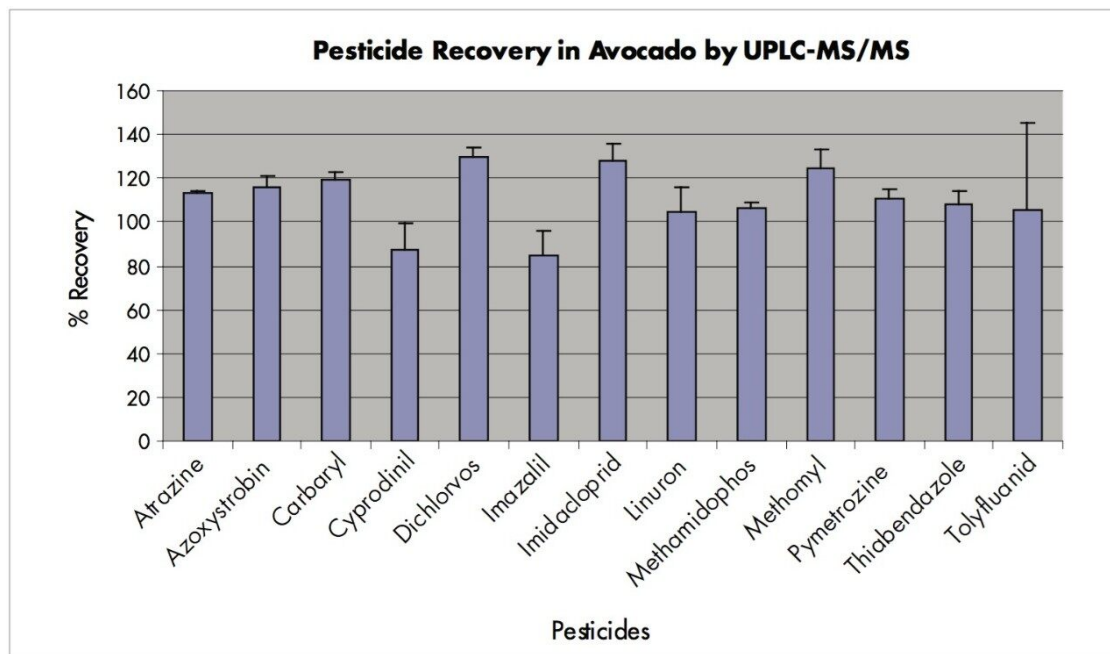

*Pesticides in Avocados by UPLC-MS/MS*

---

## Featured Products

ACQUITY UPLC System <<https://www.waters.com/514207>>

720003373, March 2010

©2019 Waters Corporation. All Rights Reserved.