

Nota applicativa

Analysis of Herbal Medicine on an ACQUITY UPLC BEH Amide Column

Waters Corporation


This is an Application Brief and does not contain a detailed
Experimental section.

Abstract

This application brief demonstrates the analysis of herbal medicine on an ACQUITY UPLC BEH amide
column.

Introduction

Compound

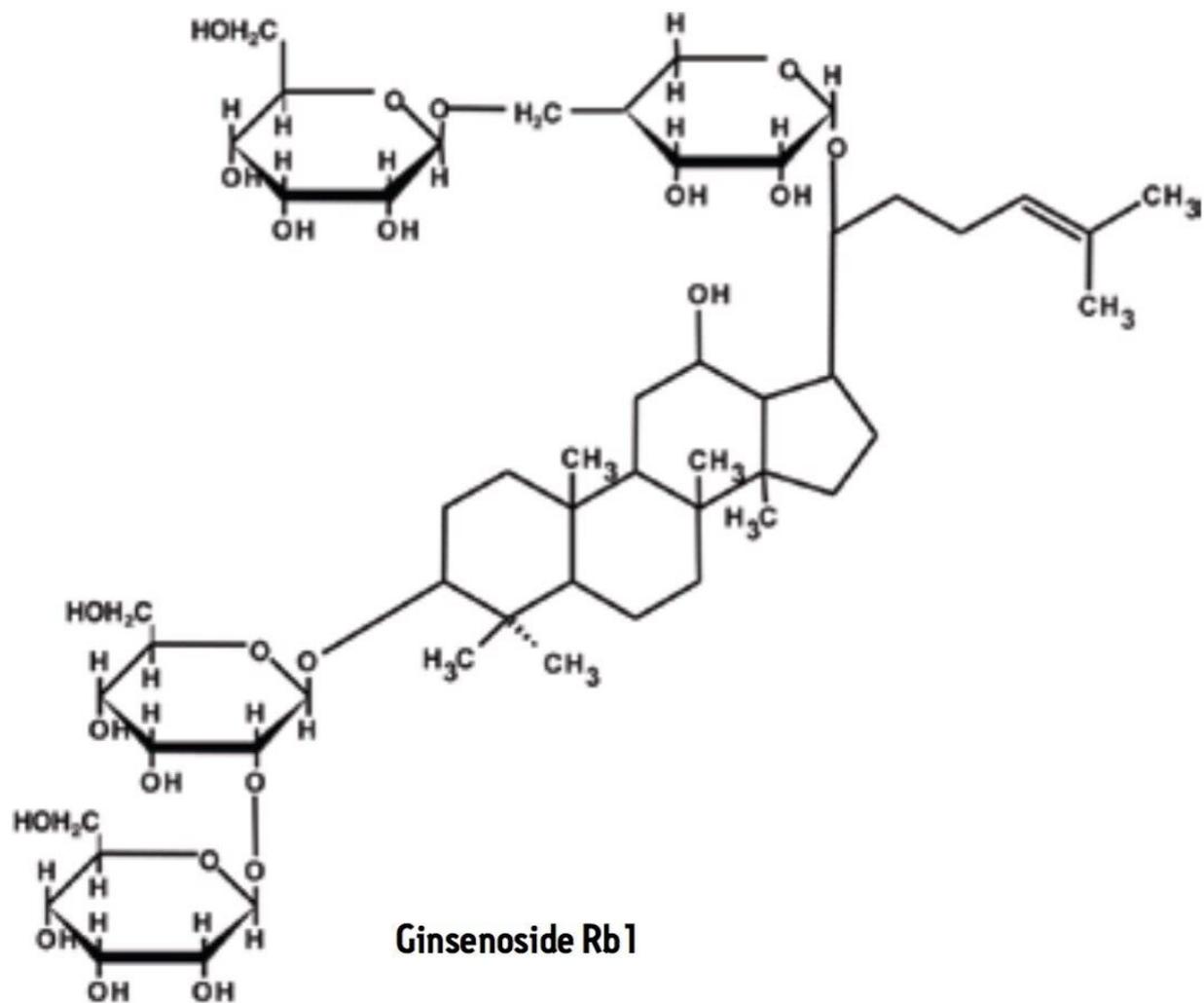

Experimental

UPLC Conditions

Column: ACQUITY UPLC BEH Amide,
1.7 μm , 2.1 x 100 mm

Part Number: 186004801

| | |
|----------------------------|---------------------------------------|
| Mobile Phase: | 80:20 MeCN:H ₂ O |
| Isocratic Flow Rate: | 0.6 mL/min |
| Column Temp.: | 60 °C |
| Sample Temp.: | 10 °C |
| Injection Vol.: | 1.7 µL; PLNO on 10 µL loop |
| Strong & Weak Needle Wash: | 95:5 MeCN:H ₂ O |
| Seal Wash: | 10:90 MeOH:H ₂ O |
| UV: | 203 nm |
| Sampling Rate: | 20 Hz |
| Filter Time Constant: | 0.2 sec |
| Total Run Time: | 2.5 min |
| Instrument: | ACQUITY UPLC with ACQUITY UPLC PDA |

Pretreatment

1. Weigh 2 g of herbal medicine powder into a centrifuge tube.
2. Add 30 mL of 60% MeOH/40% H₂O.
3. Shake for 15 min.
4. Centrifuge at 4,000 rpm for 10 min.
5. Obtain the supernatant.
6. Repeat steps 2-5 with the residue using 15 mL of 60% MeOH/40% H₂O.

7. Combine the supernatant, and make exactly 50 mL by adding 60% MeOH/40% H₂O.
8. Take 10 mL of this solution and add 3 mL of NaOH test solution (1 mol/L).
9. Let stand for 30 min.
10. Add 3 mL of HCl test solution (1 mol/L).
11. Add 60% MeOH/40% H₂O to make exactly 20 mL.

Solid-Phase Extraction

SPE Device: Sep-Pak Plus C₁₈
cartridge 360 mg
(55-105 µm)

Part Number: WAT020515

1. Condition with 2 mL MeOH.
2. Equilibrate with 2 mL of 30% MeOH/70% H₂O just before loading.
3. Load 5 mL of the solution from step 11 in the pretreatment stage.
4. Wash with 2 mL of 30% MeOH/70% H₂O.
5. Wash with 1 mL of Na₂CO test solution (1 mol/L).
6. Wash with 10 mL of 30% MeOH/70% H₂O.
7. Elute with 5 mL MeOH (this is the injection solution).

Featured Products

ACQUITY UPLC System <<https://www.waters.com/514207>>

ACQUITY UPLC PDA Detector <<https://www.waters.com/514225>>

WA64085, February 2010

©2019 Waters Corporation. All Rights Reserved.