

Application Note

Separation of Free Fatty Acids by Mixed-Mode Anion Exchange Chromatography

Kerri M. Smith, Paul D. Rainville

Waters Corporation



This is an Application Brief and does not contain a detailed Experimental section.

Abstract

This application brief demonstrates the separation of free fatty acids from C3:0–C18:0.

Benefits

The Atlantis Premier BEH C₁₈ AX Column provides excellent separation of free fatty acids

Introduction

Fatty acids are present in foods, as part of biological systems and microbiota fermentation of complex fiber, and represent an important reservoir for energy and cell wall components for mammals.¹ Their analysis spans multiple disciplines and can be performed in several ways, including gas chromatography (following derivatization of the carboxylic acid moiety),² reversed-phase liquid chromatography,³ and supercritical fluid chromatography,⁴ among others. Here, we show the separation of free fatty acids using a mixed-mode reversed-phase/anion exchange column, the Waters Atlantis Premier BEH C₁₈ AX.

Results and Discussion

The Atlantis Premier BEH C₁₈ AX is a mixed-mode reversed-phase/anion exchange column having C₁₈ ligand as well as tertiary alkylamine moieties bonded to a ethylene-bridged hybrid (BEH) silica particle.⁵ At pH below 8.0, the surface of the stationary phase is positively charged and can undergo ion exchange with negatively charged molecules, such as acids. A free fatty acid standard mix from Nu-Chek Prep, Inc. (GLC-412) was prepared together with 10 µL of concentrated propionic acid (Sigma Aldrich) and brought to 1 mL with 3:1 methanol/chloroform. The sample was further diluted 100x with methanol and placed in an autosampler vial for analysis.

The mobile phase consisted of aqueous 1 mM ammonium formate, pH 3.0 in A, and 50:50 acetonitrile/isopropanol containing 1 mM ammonium formate, pH 3.0, in B. The gradient was programmed as 0–100% B over six minutes using a gradient curve of five. The flow rate was 0.6 mL/min, column temperature was set to 60 °C, and injection volume was 3 µL. Figure 1 shows the separation of free fatty acids from C3:0 propionic acid to C18:0 steric acid. The k' of propionic acid is 2.36 and increases in order of chain length. With the pH at 3.0, lower than the pKa of

the carboxylic moiety, the aliphatic acids are uncharged and the retention is achieved by interaction with the C₁₈ ligand. However, the pore size of the Atlantis BEH C₁₈ AX material is 95 Å, smaller than the 130 Å of other BEH particles.⁵ The smaller pore size particles have greater surface area, which leads to greater retention. This is illustrated in Figure 2, where retention of fatty acids is increased in the Atlantis Premier BEH C₁₈ AX Column with surface area of 270 m²/g compared to the BEH C₁₈ Column with surface area of 185 m²/g. Additionally, the peak shape for propionic acid improved. To test injection reproducibility, retention time RSD were calculated for five replicate injections. In Figure 3, the injections are stacked and retention time RSD ranged from 0.02–0.14% for all the analytes.

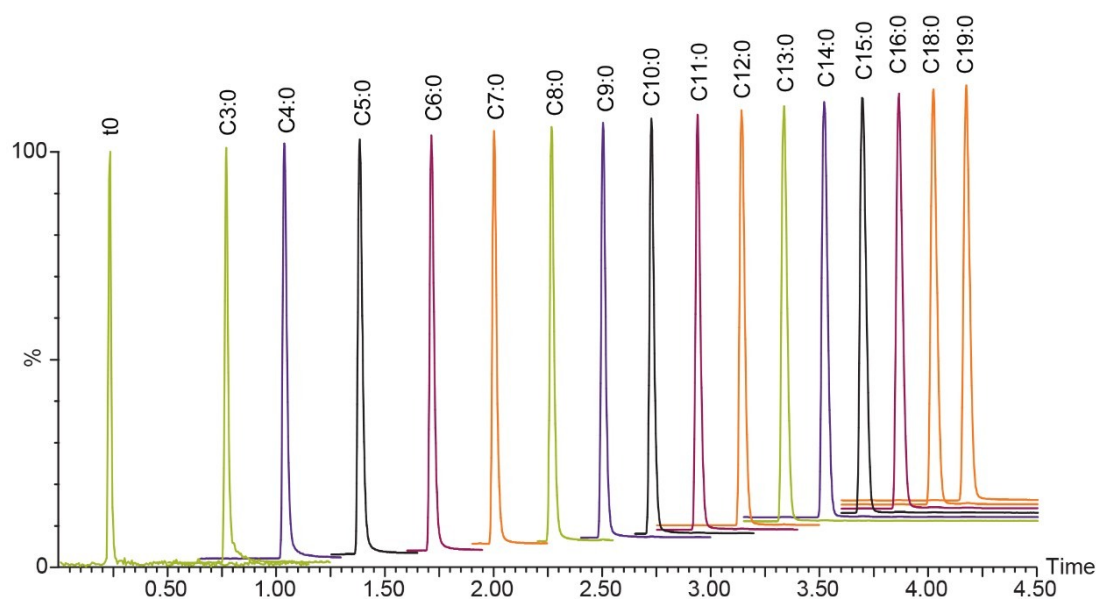


Figure 1. Separation of free fatty acids in a standard mix.

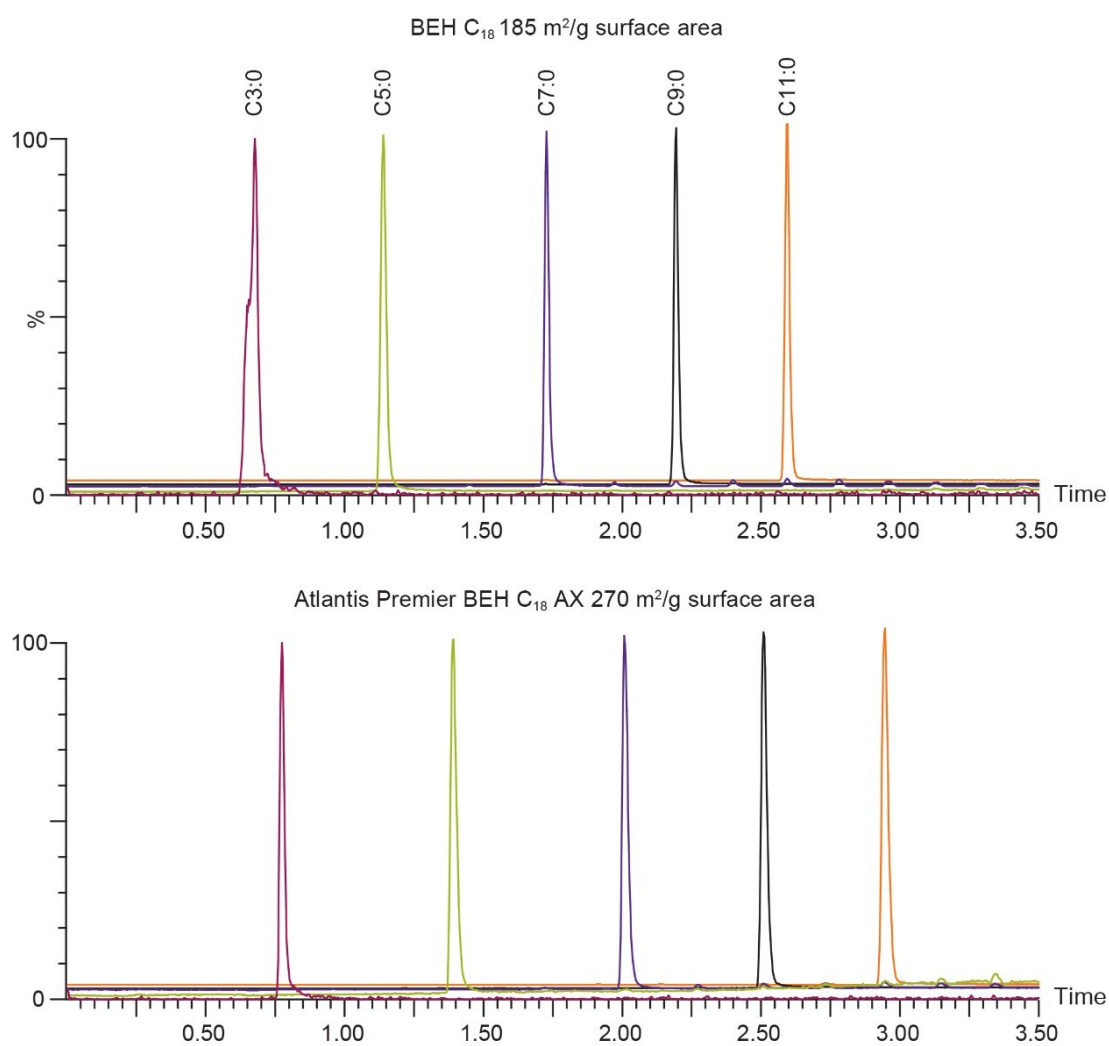


Figure 2. Effect of surface area on retention of free fatty acids. Comparison of the BEH C₁₈, 2.1 x 50 mm, 1.7 μ m column (top) with the Atlantis Premier BEH C₁₈ AX 2.1 x 50 mm, 1.7 μ m column (bottom). Retention is increased with higher surface area.

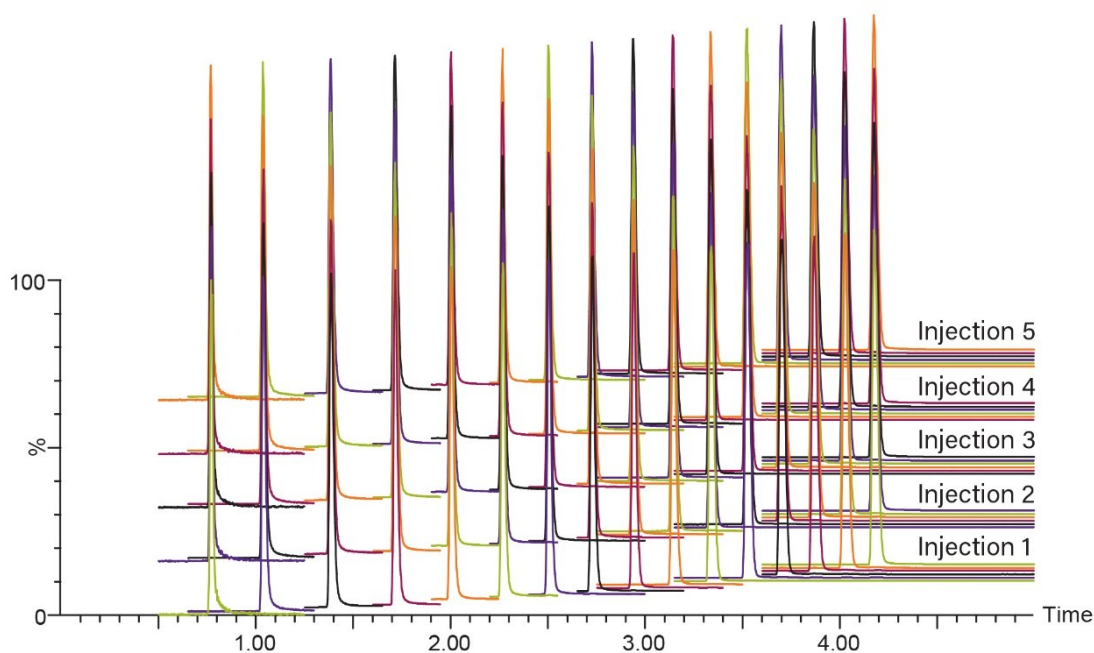


Figure 3. Injection reproducibility of free fatty acids.

Conclusion

The Atlantis Premier BEH C₁₈ AX Column provides good retention and separation of free fatty acids with excellent injection reproducibility.

References

1. Murray R.K.; Bender D.A.; *et al.* Harper's *Illustrated Biochemistry*, 28th ed. New York: McGraw-Hill, 2009; pp 121-142.
2. Ichihara K. *et al.* Preparation of Fatty Acid Methyl Esters for Gas-Liquid Chromatography. *J. Lipid Res.* 2010, 51, 635–640.
3. Nikolova-Damyanova B. *Reversed-Phase High- Performance Liquid Chromatography: General Principles and Application to the Analysis of Fatty Acids and Triacylglycerols in Advances in Lipid Methodology*, Volume 4.

Woodhead: Philadelphia, 1997, pp 193–251.

4. Yang Y. *et al.* Advances of Supercritical Fluid Chromatography in Lipid Profiling *J. Pharm. Anal.* 2019, 9, 1–8.
5. Walter, T.H. *et al.* A New Mixed-Mode Reversed Phase/Anion Exchange Stationary Phase Based on Hybrid Particles. Waters Application Note, 720006742 <<https://www.waters.com/nextgen/in/en/library/application-notes/2020/a-new-mixed-mode-reversed-phase-anion-exchange-stationary-phase-based-on-hybrid-particles.html>> , 2020.

Featured Products

720006747, January 2020