

Fast and Simple Bioanalytical Plasma Sample Extraction Using Oasis® PRiME HLB SPE for High Analyte Recovery

Mary Trudeau, Margot Lee

Waters Corporation

这是一份应用简报，不包含详细的实验部分。

Abstract

This application describes a fast and simple method for bioanalytical sample preparation from plasma. Utilizing Oasis PRiME HLB solid-phase extraction (SPE), the protocol requires no method development and achieves high recovery and low matrix effects for several small molecule pharmaceuticals.

Experimental

Oasis PRiME HLB SPE Protocol
96-well plate, 3 mg sorbent per well,
P/N 186008052



Sample pretreatment:

200 μ L 4% H_3PO_4 to 200 μ L
plasma sample

Load:

400 μ L pre-treated sample

Wash:

200 μ L 95:5 Water: Methanol

Elute:

2x 25 μ L 90:10
Acetonitrile/Methanol

Dilute:

50 μ L H_2O

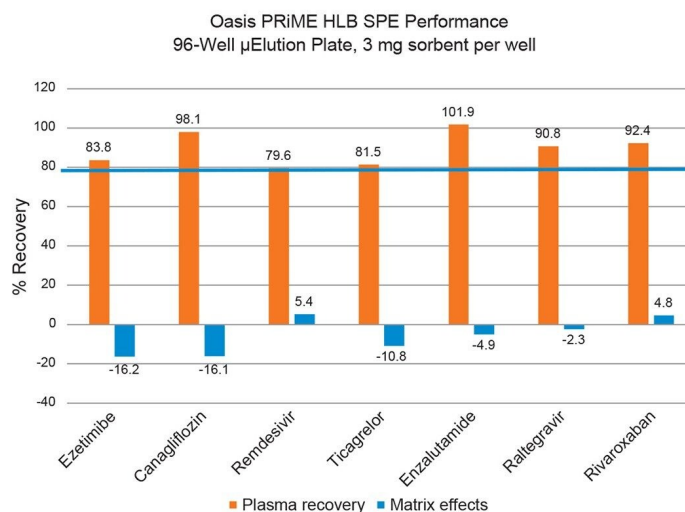


Figure 1. Demonstration of Oasis PRiME HLB SPE bioanalytical performance, requiring no optimization of the starting protocol, with high recovery (> 79%) and low matrix effects (-16.2 to 5.4%) for several pharmaceuticals extracted from plasma.

Results and Discussion

High Analyte Recovery with low matrix from plasma is achieved Using Oasis PRiME HLB SPE for top-selling small molecule pharmaceuticals in 2023.¹

References

1. Williams, RE and Leatherwood HM [Top 200 Small Molecule Drugs by Retail Sales in 2023](https://bpb-us-e2.wpmucdn.com/sites.arizona.edu/dist/9/130/files/2024/05/2023Top200SmallMoleculePosterV5.pdf) <<https://bpb-us-e2.wpmucdn.com/sites.arizona.edu/dist/9/130/files/2024/05/2023Top200SmallMoleculePosterV5.pdf>> . (accessed 21 June 2024).

Featured Products

720008423, July 2024



© 2024 Waters Corporation. All Rights Reserved.

[使用条款](#) [隐私策略](#) [商标](#) [招聘](#) [法律和隐私声明](#) [危险化学品生产经营许可证](#) [Cookie](#) [Cookie 设置](#)

[沪ICP备06003546号-2](#) [京公网安备 31011502007476号](#)