


Note d'application

## Atorvastatin - XTerra RP<sub>18</sub>

---

Waters Corporation

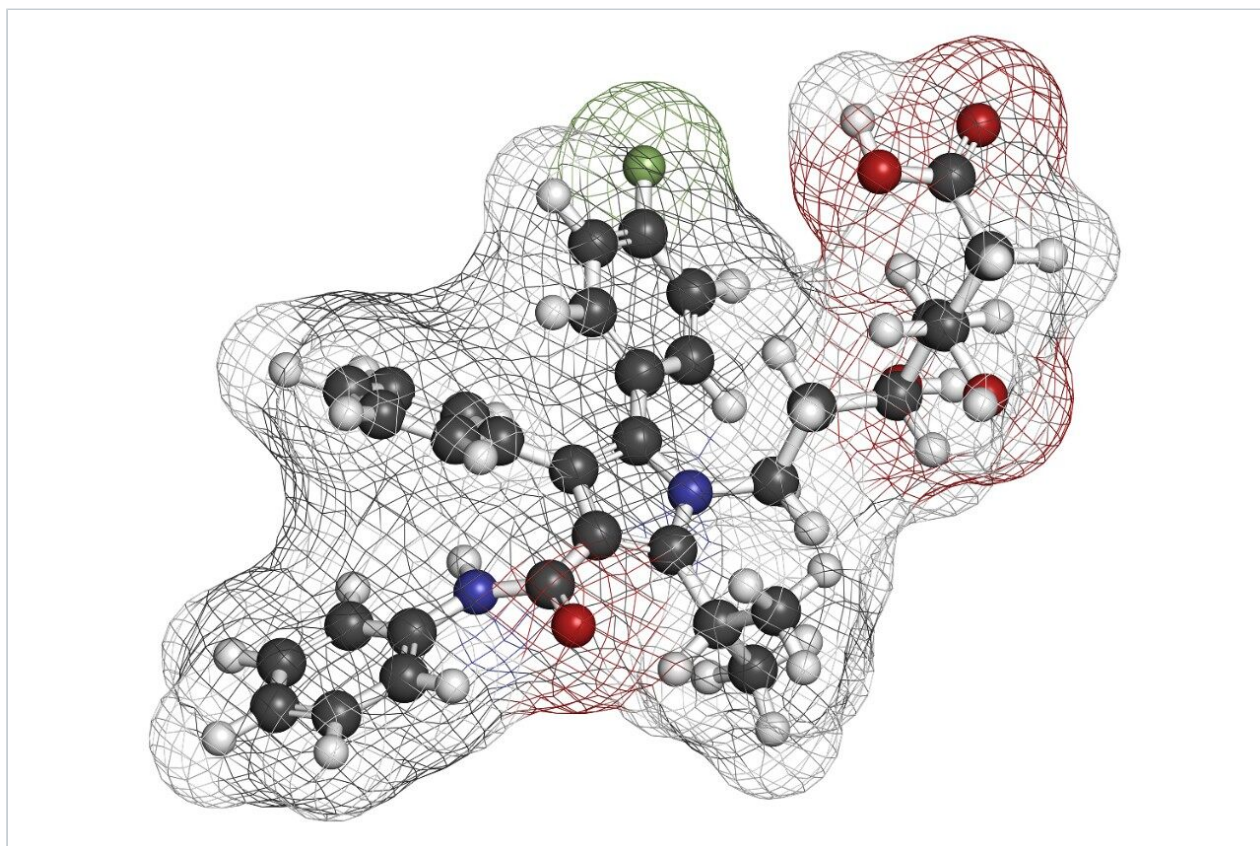

This is an Application Brief and does not contain a detailed Experimental section.

---


### Abstract

This application brief demonstrates the analysis of atorvastatin using XTerra RP<sub>18</sub> Columns.

---

## Introduction

In this study the compound analyzed is atorvastatin.

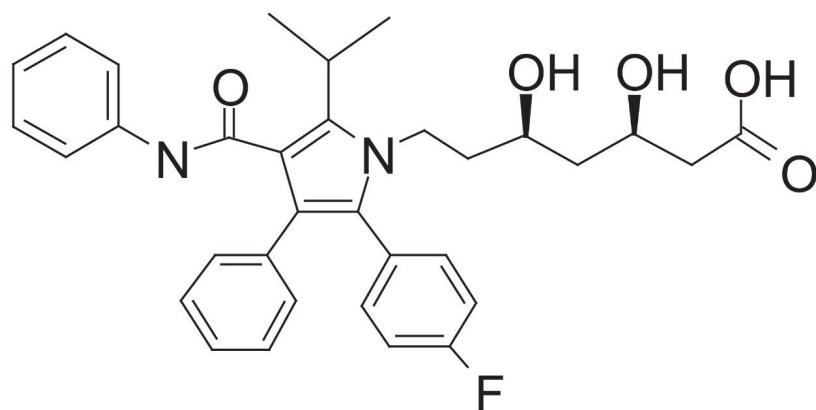

Atorvastatin

---

## Experimental

### Conditions

Column:	XTerra RP <sub>18</sub> 4.6 x 100 mm, 3.5 μm
Part number:	186000438
Mobile phase A:	10 mM NH <sub>4</sub> COOH, pH 3.0
Mobile phase A:	10 mM NH <sub>4</sub> HCO <sub>3</sub> , pH 7.0
Mobile phase A:	10 mM NH <sub>4</sub> HCO <sub>3</sub> , pH 10.0
Mobile phase B:	ACN


### Isocratic Conditions as Indicated

Flow rate:	1.4 mL/min
Injection volume:	2 $\mu$ L
Sample Concentration:	0.25 mg/mL
Temperature:	30 °C
Detection:	UV @ 254 nm
Instrument:	Alliance 2795 HT, 2996 PDA

Mobile Phase pH	USP Tailing
3.0	1.04
7.0	1.19
10.0	1.11

---

## Results and Discussion

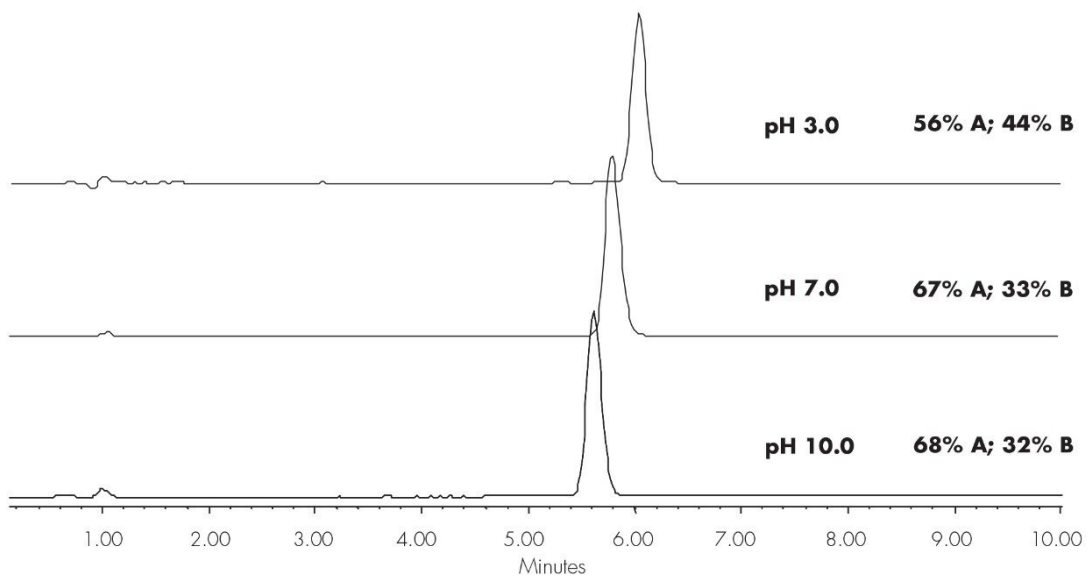

---

## Featured Products

WA20738.014, June 2002