


Atorvastatin - XTerra Phenyl

Waters Corporation

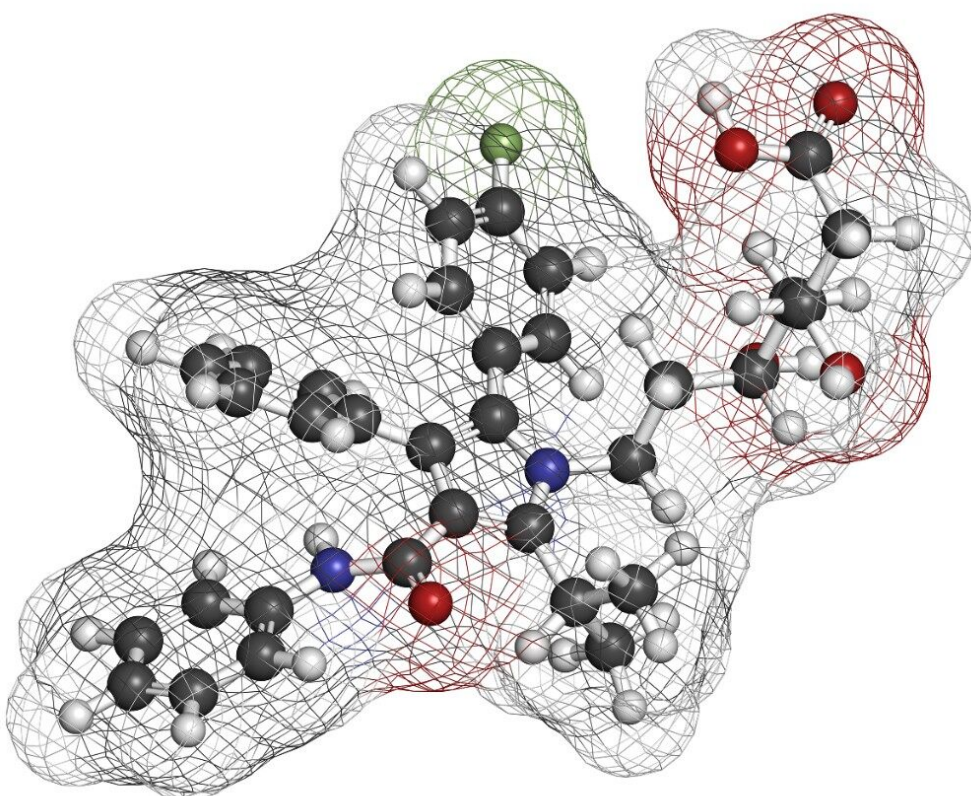


This is an Application Brief and does not contain a detailed Experimental section.

Abstract

This application brief demonstrates the analysis of atorvastatin using XTerra Columns.

Introduction

In this study the compound analyzed is atorvastatin.

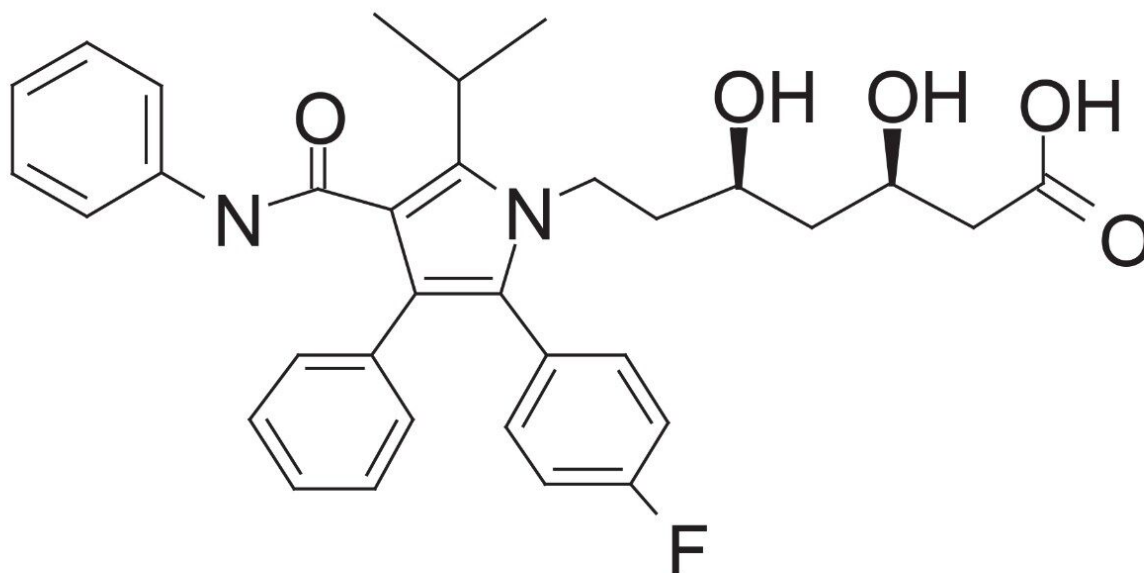

Atorvastatin

Experimental

Conditions

Column: XTerra Phenyl, 4.6 x 100 mm,
3.5 μ m

Part number: 186001139

Mobile phase A: 10 mM NH_4COOH , pH 3.0

Mobile phase A:	10 mM NH ₄ HCO ₃ , pH 7.0
-----------------	---

Mobile phase A:	10 mM NH ₄ HCO ₃ , pH 10.0
-----------------	--

Mobile phase B:	ACN
-----------------	-----

Isocratic Conditions as
Indicated

Flow rate:	1.4 mL/min
------------	------------

Injection volume:	2 µL
-------------------	------

Sample Concentration:	0.25 mg/mL
-----------------------	------------

Temperature:	30 °C
--------------	-------

Detection:	UV @ 254 nm
------------	-------------

Instrument:	Alliance 2695, 2996 PDA
-------------	-------------------------


Mobile Phase pH	USP Tailing
-----------------	-------------

3.0	0.95
-----	------

7.0	1.11
-----	------

10.0	1.06
------	------

Results and Discussion

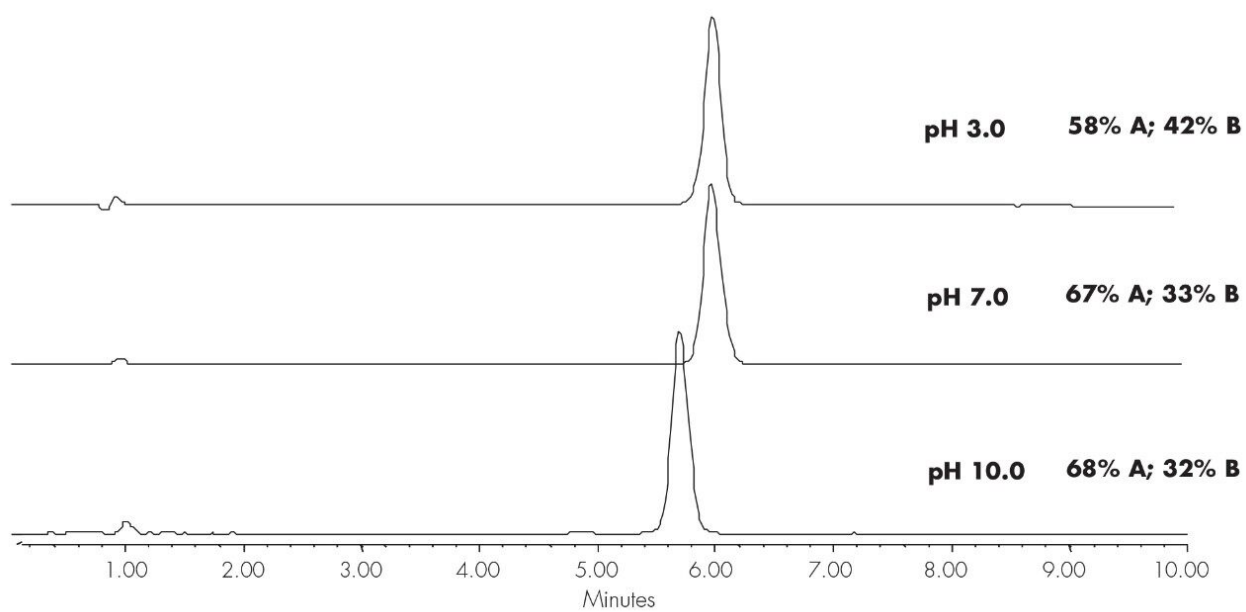

Featured Products

WA20738.015, June 2002