


Nota applicativa

## Cytochrome C Tryptic Digest – pH 10.0, XTerra RP<sub>18</sub>

---

Waters Corporation

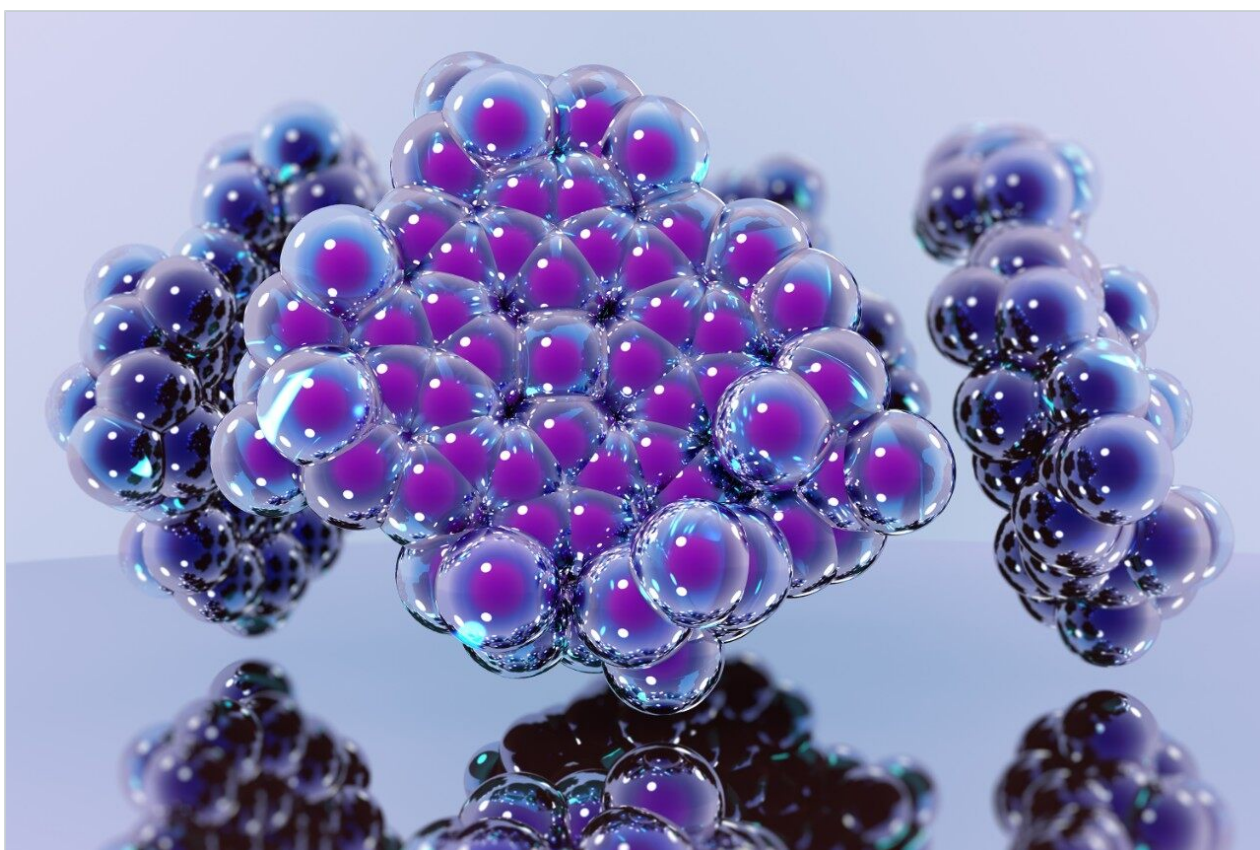

This is an Application Brief and does not contain a detailed Experimental section.

---

### Abstract

This application brief highlights the analysis of cytochrome C using XTerra RP<sub>18</sub> columns.

---

## Introduction

Cytochrome C has been studied in this application brief.

---

## Experimental

### HPLC Conditions

Column:	XTerra RP <sub>18</sub> 4.6 x 50 mm, 3.5 µm (p/n: 186000434)
Mobile phase A:	H <sub>2</sub> O
Mobile phase B:	ACN
Mobile phase C:	100 mM HCOOH, pH 10.0
Flow rate:	0.75 mL/min, 0.2 mL/min to MS
Injection volume:	50 µL (25 µg)
Temperature:	Ambient
Detection:	MS ESI+
Instrument:	Alliance 2790, Micromass Quattro Ultima

## Gradient


Time (min)	Profile		
	%A	%B	%C
0.0	90	0	10
45.0	50	40	10
45.0	10	80	10
50.0	10	80	10

## MS Conditions

MS system:	Micromass ZQ
Ion source:	ESI+
Capillary:	3.5 kV
Cone:	30 V
Extractor:	3.0
Source temp.:	150 °C
Desolvation temp.:	400 °C
Cone gas flow:	60 L/hr
Desolvation gas flow:	550 L/hr

---

## Results and Discussion

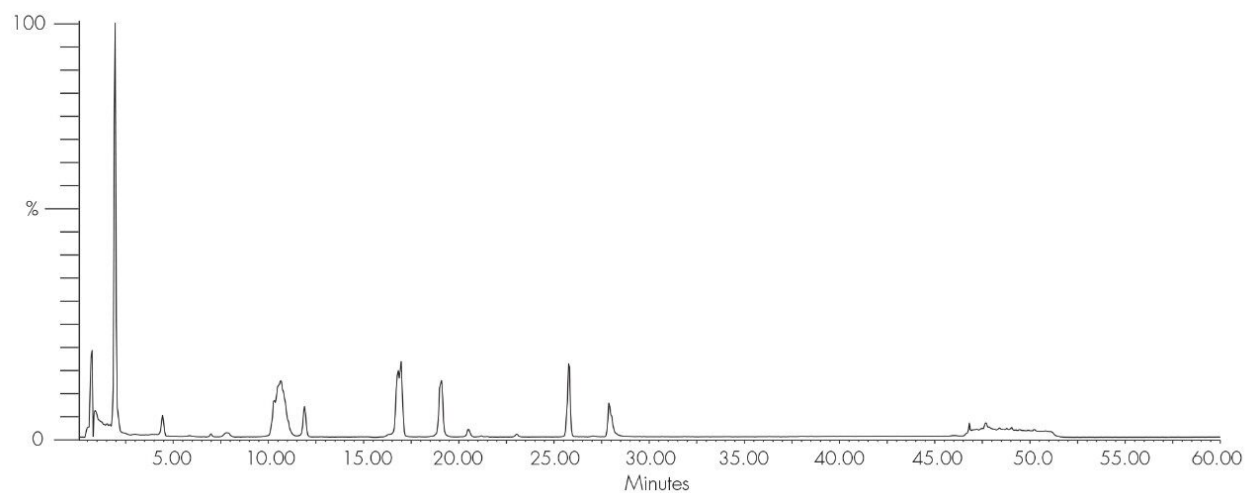

---

## Featured Products

Alliance HPLC <<https://www.waters.com/514248>>

WA20738.030, June 2002