

## Dextromethorphan

---

Waters Corporation



This is an Application Brief and does not contain a detailed Experimental section.

---


### Abstract

This application brief highlights the analysis of dextromethorphan using XTerra RP<sub>18</sub> columns.

---

## Introduction

Dextromethorphan has been analyzed in this study.

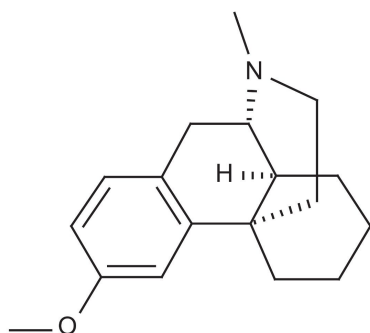

Dextromethorphan

---

## Experimental

### HPLC Conditions


Column:	XTerra RP <sub>18</sub> 4.6 x 150 mm, 5 µm (p/n: 186000492)
Mobile phase:	At pH 3.0: H <sub>2</sub> O/ACN/100 Mm NH <sub>4</sub> COOH, pH 3.0, 40:50:10 At pH 7.0: H <sub>2</sub> O/ACN/100 mM NH <sub>4</sub> HCO <sub>3</sub> , pH 7.0, 50:40:10 At pH 10.0: H <sub>2</sub> O/ACN/100 mM NH <sub>4</sub> HCO <sub>3</sub> , pH 10.0, 22:68:10
Flow rate:	1.0 mL/min
Injection volume:	5 µL of 250 µg/mL
Temperature:	30 °C

Detection:UV @ 228 nm

Instrument:Alliance 2695, 2996 PDA

Mobile Phase pH	USP Tailing
3.0	1.25
7.0	1.39
10.0	1.39

Results and Discussion

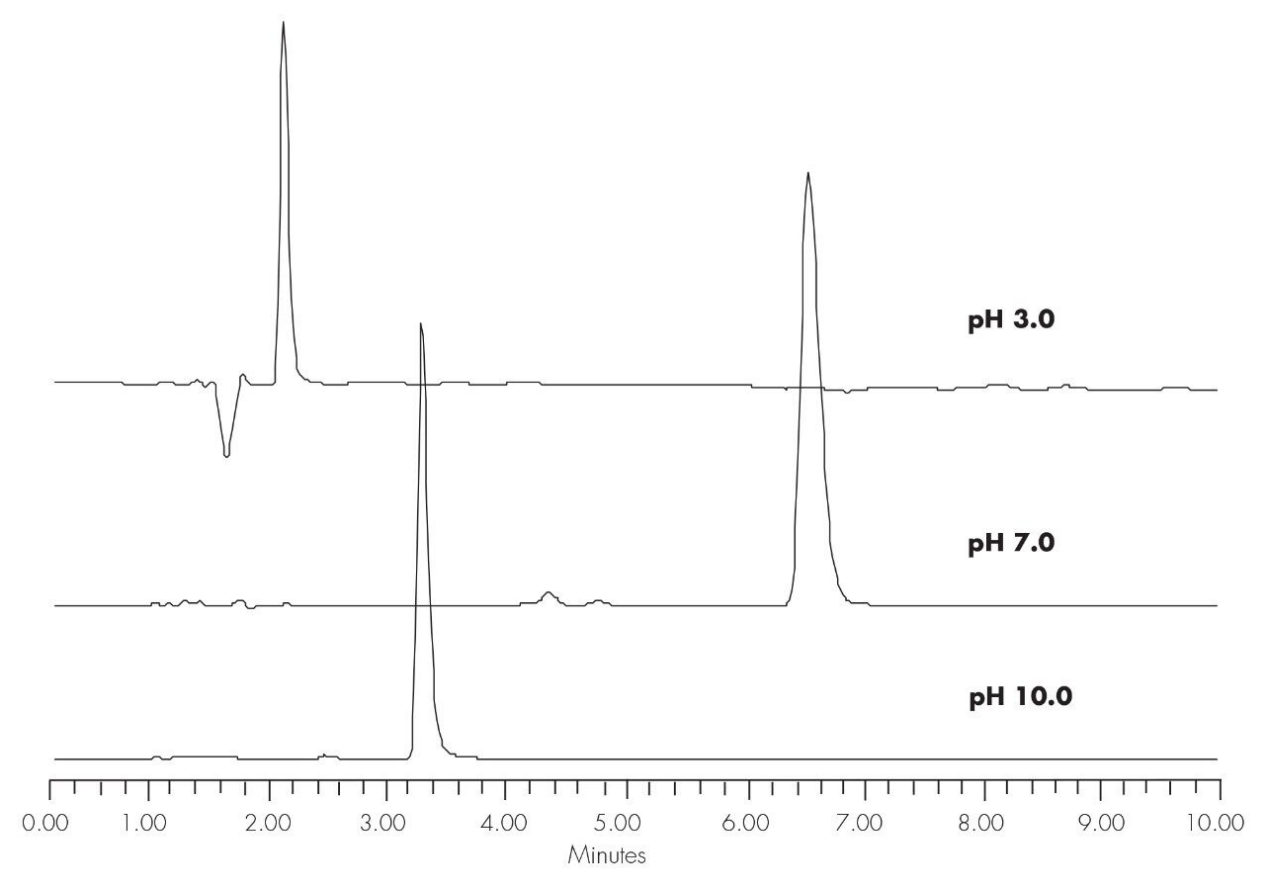

---

## Featured Products

Alliance HPLC <<https://www.waters.com/514248>>

WA20738.033, June 2002

© 2021 Waters Corporation. All Rights Reserved.