

Metoprolol - LC/UV

Waters Corporation




This is an Application Brief and does not contain a detailed Experimental section.

Abstract

This application brief highlights the analysis of metoprolol using XTerra RP₁₈ columns.

Introduction

Metoprolol has been analyzed in this study.

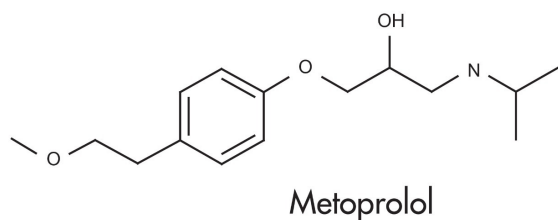

Experimental

HPLC Conditions


Column:	XTerra RP ₁₈ 4.6 x 150 mm, 5 µm (p/n: 186000492)
Mobile phase:	At pH 3.0: H ₂ O/ACN/100 mM, NH ₄ COOH, pH 3.0, 55:35:10 At pH 7.0: H ₂ O/ACN/100 mM, NH ₄ HCO ₃ , pH 7.0, 35:55:10 At pH 10.0 H ₂ O/ACN/100 mM, NH ₄ HCO ₃ , pH 10.0, 60:30:10
Flow rate:	1.0 mL/min
Injection volume:	5 µL of 250 µg/mL
Temperature:	30 °C
Detection:	UV @ 270 nm

Instrument:

Alliance 2695, 2996 PDA

Mobile Phase pH	USP Tailing
3.0	1.17
7.0	1.33
10.0	1.14

Results and Discussion

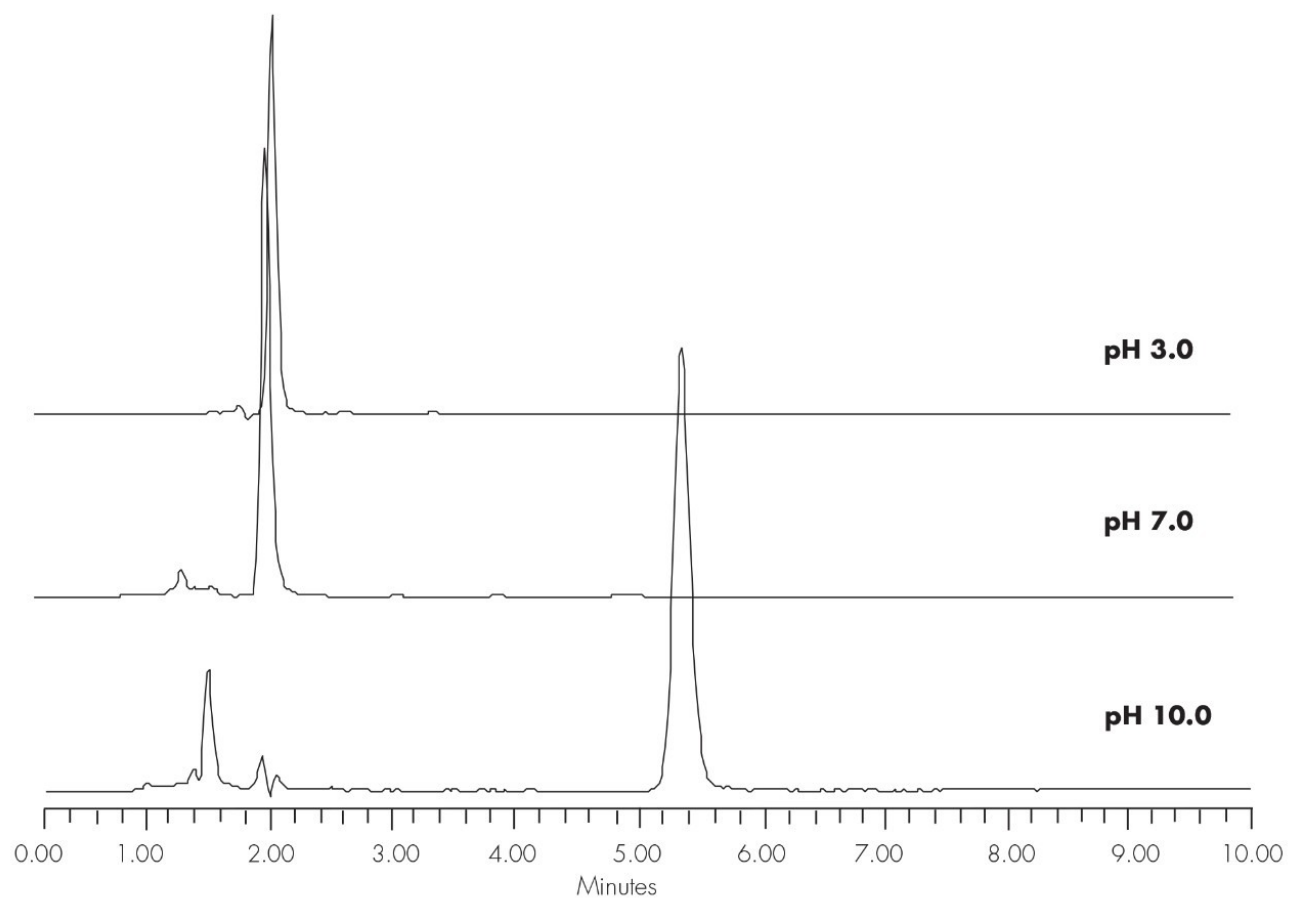

Featured Products

· [Alliance HPLC <https://www.waters.com/514248>](https://www.waters.com/514248)

WA20738.071, June 2002



© 2021 Waters Corporation. All Rights Reserved.