

Procaine

Waters Corporation




This is an Application Brief and does not contain a detailed Experimental section.

Abstract

This application brief highlights the analysis of procaine using XTerra RP₁₈ columns.

Introduction

Procaine has been analyzed in this application brief.

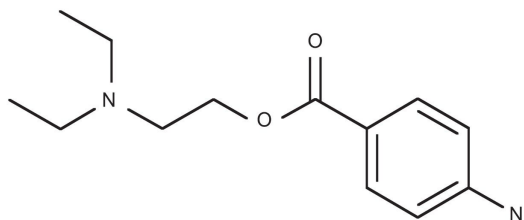

Procaine

Experimental

HPLC Method

Column:	XTerra RP ₁₈ 4.6 x 150 mm, 5 µm (p/n: 186000492)
Mobile phase:	At pH 3.0: H ₂ O/ACN/100 mM NH ₄ COOH, pH 3.0, 65:25:10 At pH 7.0: H ₂ O/ACN/100 mM NH ₄ HCO ₃ , pH 7.0, 35:55:10 At pH 10.0: H ₂ O/ACN/100 mM NH ₄ HCO ₃ , pH 10.0, 60:30:10
Flow rate:	1.0 mL/min
Injection volume:	5 µL of 250 µg/mL


Temperature: 30 °C

Detection: UV @ 280 nm

Instrument: Alliance 2695, 2996 PDA

Mobile Phase pH	USP Tailing
3.0	1.28
7.0	1.23
10.0	1.08

Results and Discussion

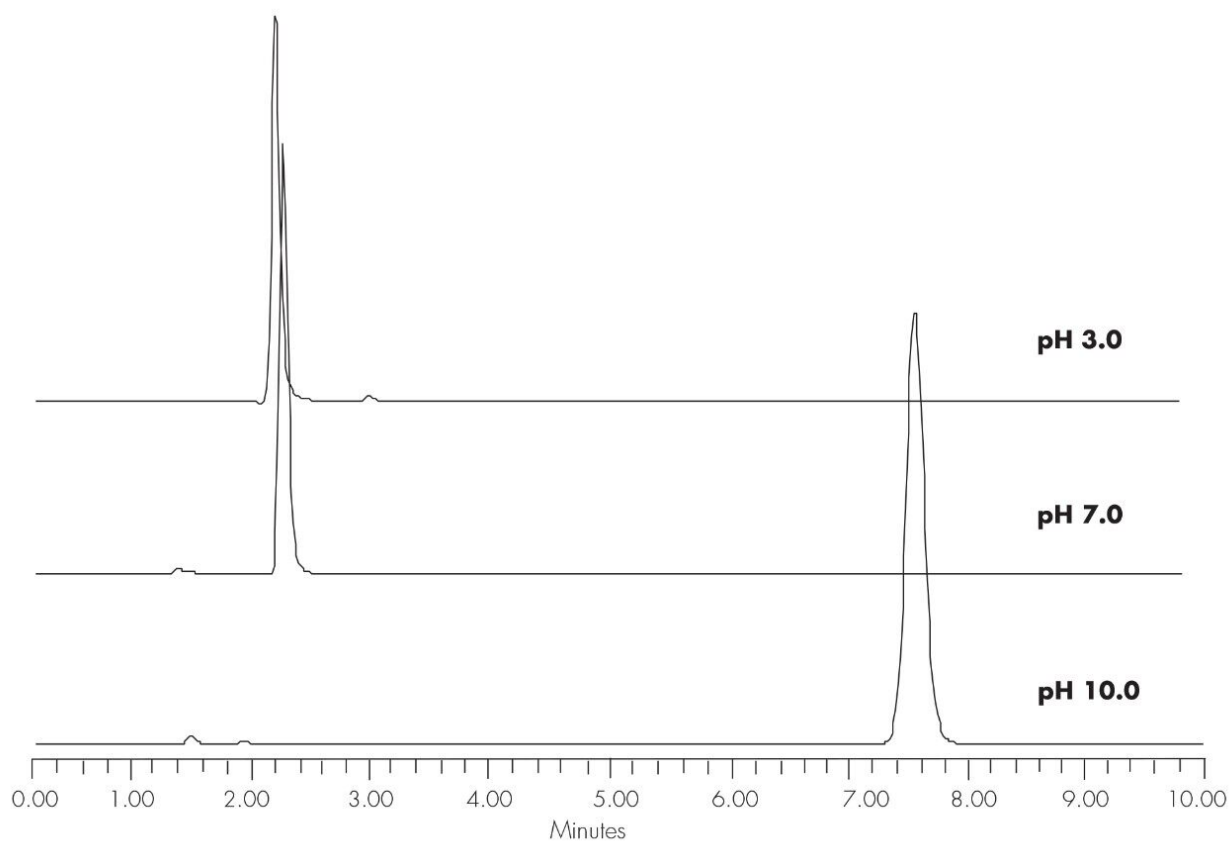

Featured Products

Alliance HPLC <<https://www.waters.com/514248>>

WA20738.091, June 2002

© 2021 Waters Corporation. All Rights Reserved.