


Ranitidine in Rat Plasma

Waters Corporation

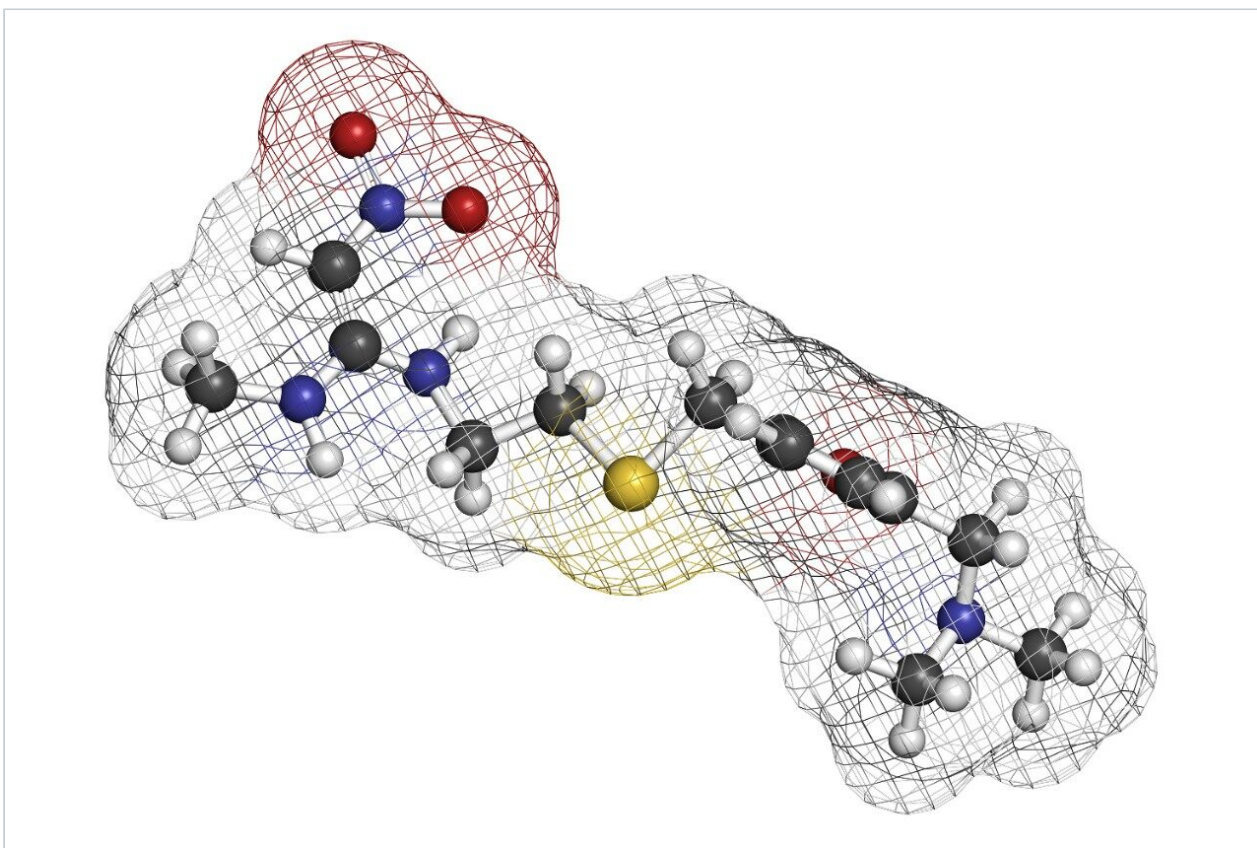


This is an Application Brief and does not contain a detailed Experimental section.

Abstract

This application brief demonstrates analysis of ranitidine in rat plasma.

Detection:	MS ESI+
Instrument:	Alliance 2790, Micromass Quattro Ultima
Ion source:	ESI+
Source temperature:	150 °C
Gas cell:	1.5e-3 mbar, 20eV
Desolvation temperature:	350 °C
Cone gas flow:	150 L/hr
Drying gas flow:	600 L/hr
Cone voltage:	20V

Oasis® MCX Extraction Method
Oasis® MCX Extraction Plate, 10 mg/96-well
Part Number 186000259

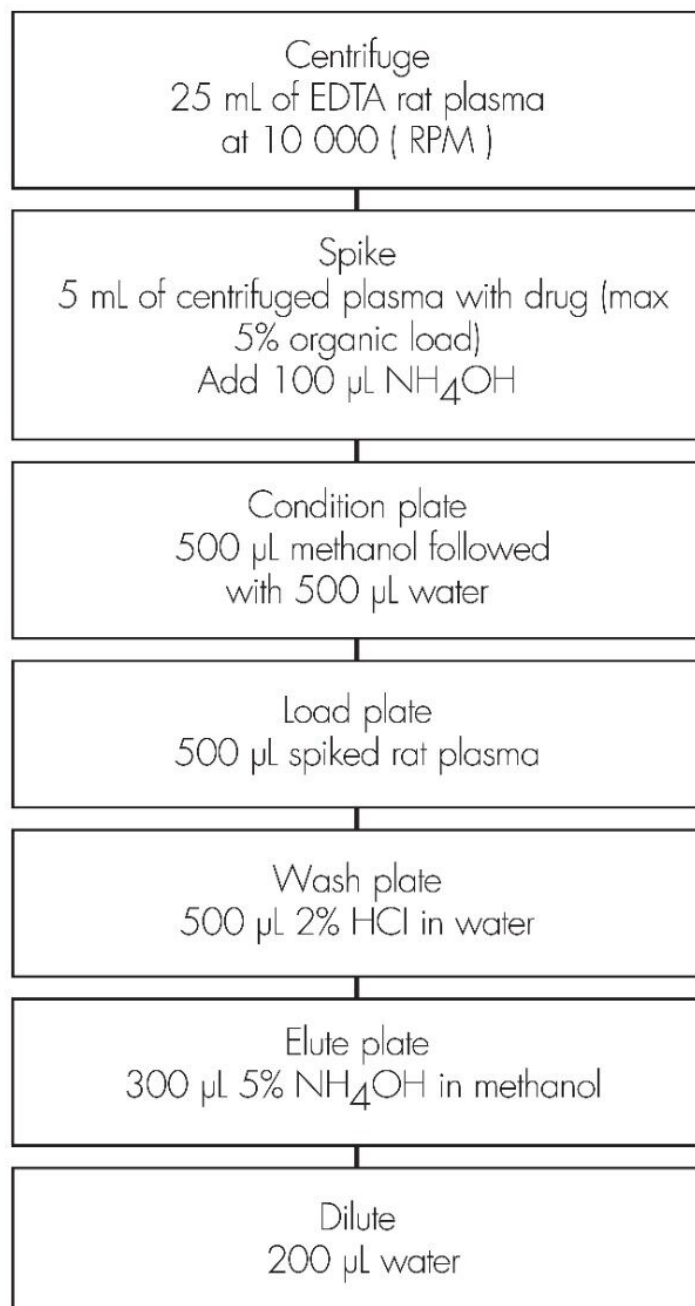


Results and Discussion

CID mass spectra


cid 15
S3 1 (1.005)

Daughters of 315ES+
2.99e8

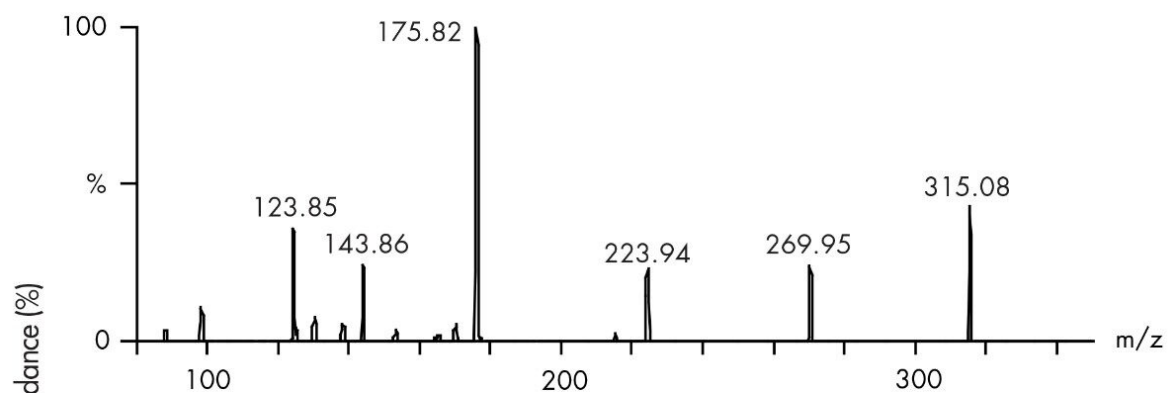

IS daughter of 271.0
2 88 (1.639) Cm (87:90)

Daughters of 271ES+
6.87e4

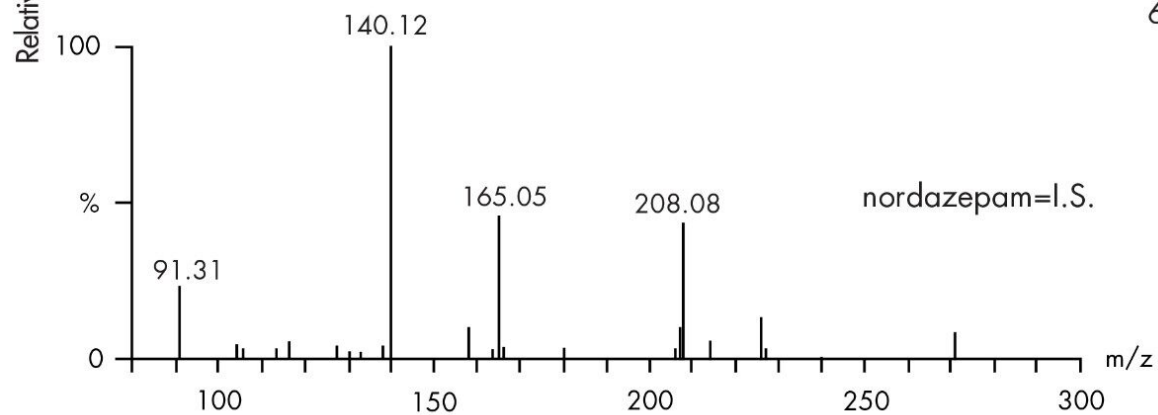

Featured Products

Alliance HPLC System <<https://www.waters.com/534293>>

WA20738.096, June 2002

