

Nota applicativa

## Metoprolol Tartrate and Oxprenolol

---

Waters Corporation



This is an Application Brief and does not contain a detailed Experimental section.

---

### Abstract

This application brief demonstrates analysis of Metoprolol Tartrate and Oxprenolol.

---


## Introduction

The compounds analyzed are - Metoprolol Tartrate and Oxprenolol.


USP Tailing Factors

1. 1.29

2. 1.98

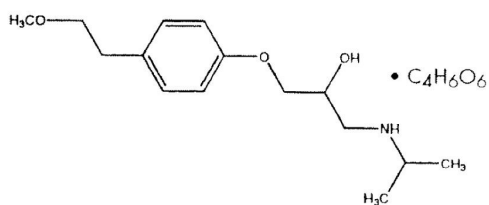

**1. Metoprolol tartrate**

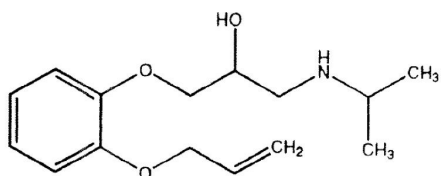

**2. Oxprenolol**

---

## Experimental

### HPLC Method


Column:

Symmetry Shield RP<sub>18</sub>, 4.6 x 250 mm, 5 μm

Part number:	186000112
Mobile phase:	1.92 g of 1-pentansulfonic acid sodium salt monohydrate, 0.27 g of sodium acetate trihydrate and 1.14 mL glacial acetic acid in 940 mL of water; 54% methanol/46% buffer
Flow rate:	1.0 mL/min
Injection volume:	30 $\mu$ L
Detection:	UV @ 254 nm

---

## Results and Discussion

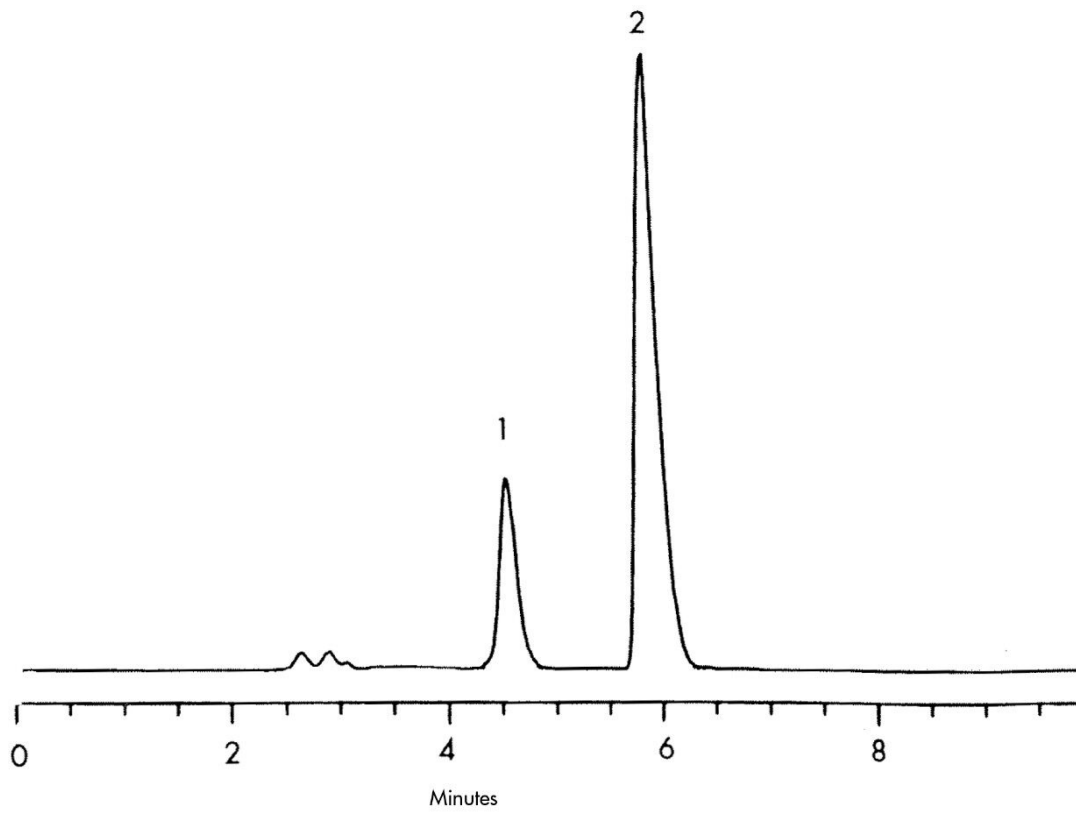

---

## Featured Products

WA31763.101, June 2003