

Nota applicativa

Paracetamol

Waters Corporation




This is an Application Brief and does not contain a detailed Experimental section.

Abstract

This application brief highlights the analysis of paracetamol using SymmetryShield RP₁₈.

Introduction

This application brief highlights the analysis of paracetamol using SymmetryShield RP₁₈.

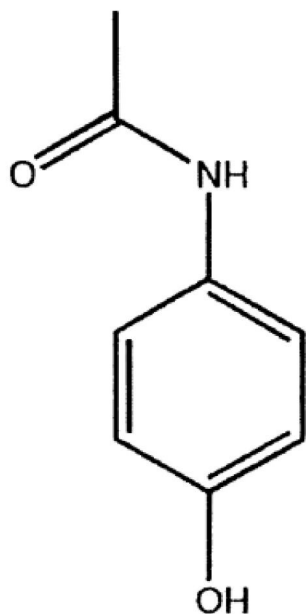

1. Paracetamol


Experimental

HPLC Method

Column:	SymmetryShield RP ₁₈ , 4.6 x 150 mm, 5 μm (p/n: 186000109)
Mobile phase:	20 mM phosphate buffer, pH 7.0/acetonitrile 85:15

Flow rate: 1.0 mL/min
Injection volume: 5 μ L
Detection: UV @ 254 nm
USP tailing factor: 1.106

Results and Discussion

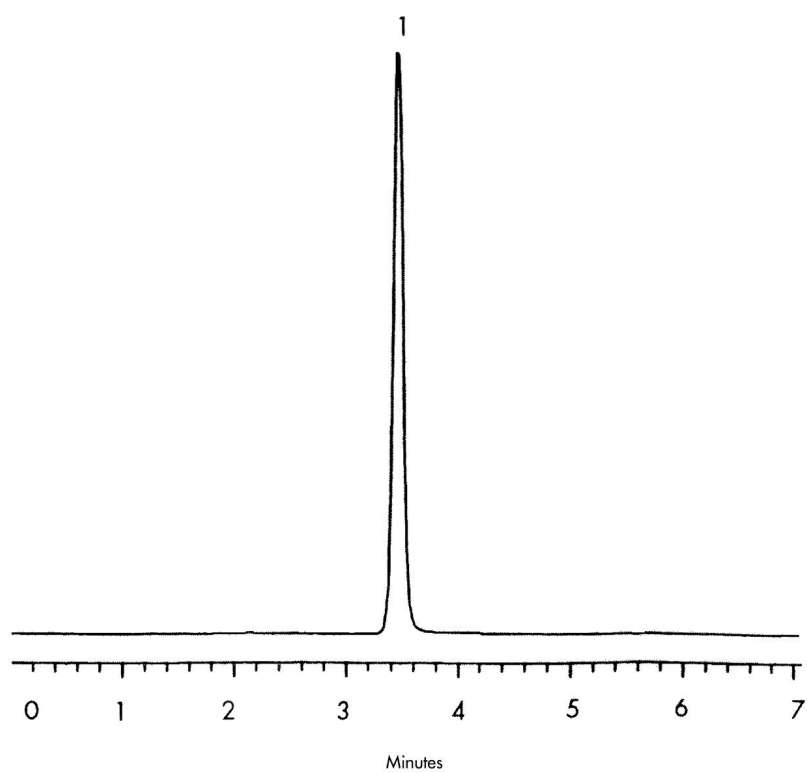

Featured Products

WA31763.119, June 2003

© 2021 Waters Corporation. All Rights Reserved.