

응용 자료

## Salicylic Acid and Related Compounds

---

Waters Corporation



This is an Application Brief and does not contain a detailed Experimental section.

---

### Abstract


This application brief demonstrates analysis of salicylic acid and related compounds.

---


## Introduction

The compounds used in this study are –


Compounds	USP Tailing Factors
Phthalic Acid	1.7
Salicylic Acid	1.1
Benzoic Acid	1.4

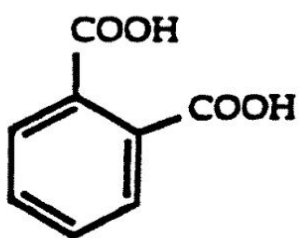

**1. Phthalic Acid**

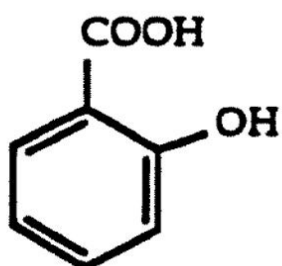

**2. Salicylic Acid**

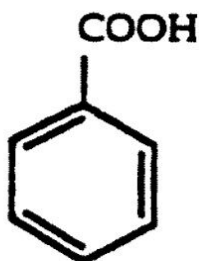

**3. Benzoic Acid**

---

## Experimental

### HPLC Method


Column:

Symmetry C<sub>18</sub>, 3.9 x 150 mm, 5 μm

Part numbers:	WAT046980
Mobile phase:	50 mM potassium phosphate pH 3.0/acetonitrile 70:30
Temperature:	40 °C
Flow rate:	0.7 mL/min
Injection volume:	6 µL of phthalic acid (0.5 mg/mL), benzoic and salicylic acids (1 mg/mL)
Detection:	UV @ 254 nm

---

## Results and Discussion

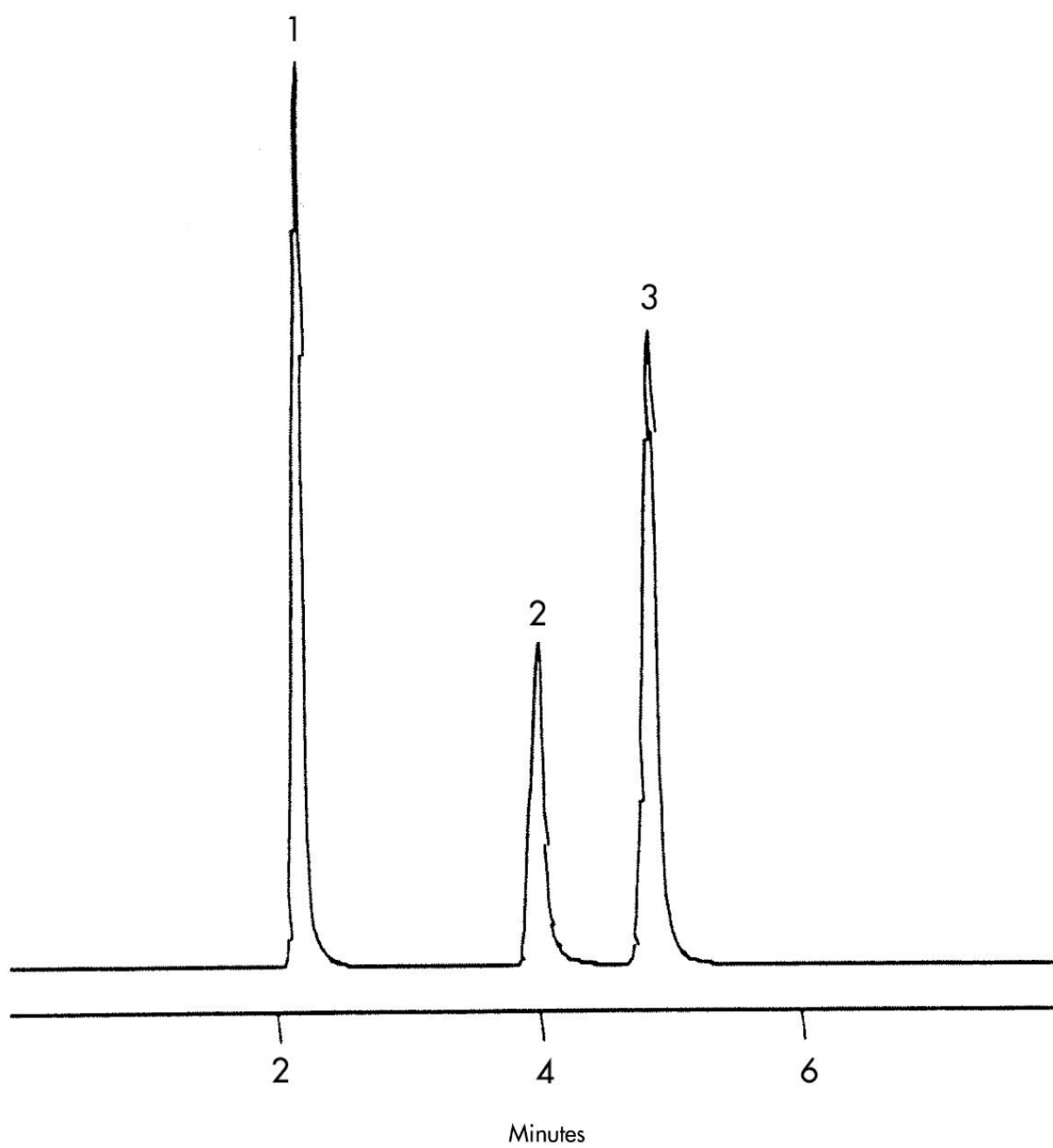

---

## Featured Products

WA31763.146, June 2003

© 2021 Waters Corporation. All Rights Reserved.