

Nota applicativa

## Steroids - Symmetry C<sub>8</sub> Column

---

Waters Corporation



This is an Application Brief and does not contain a detailed Experimental section.

---

### Abstract


This application brief demonstrates analysis of steroids using Symmetry C<sub>8</sub> Column.

---


## Introduction

The compounds used in this study are –


Compounds	USP Tailing Factors
Prednisolone	1.08
Betamethasone	1.18
Prednisolone Acetate	1.1
Betamethasone Valerate	1.31

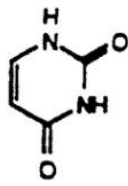


Vo: Uracil

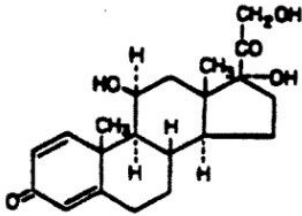


1. Prednisolone

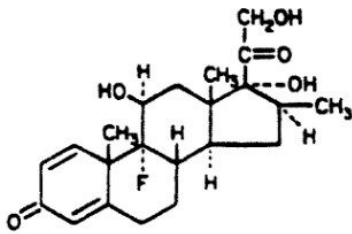

2. Betamethasone

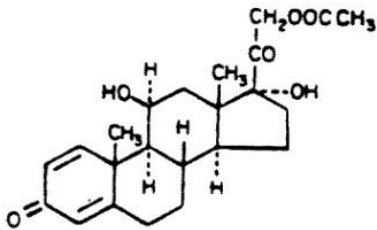

3. Prednisolone Acetate

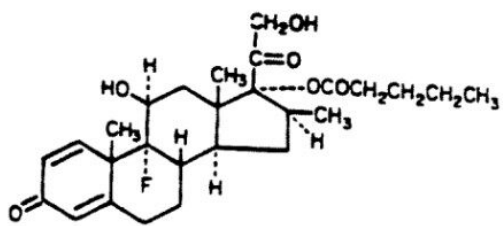

4. Betamethasone Valerate

---


Experimental

## HPLC Method

Column:	Symmetry C <sub>8</sub> ,3.9 x 150 mm, 5 μm
Guard column:	Sentry Guard Column 3.9 x 20 mm, 5μm
Part numbers:	Column - WAT046970, Guard - WAT054250
Mobile phase:	Acetonitrile/water 50:50
Flow rate:	1.0 mL/min
Injection volume:	5μL of 10 μg/mL sample
Detection:	UV @ 254 nm

---

## Results and Discussion

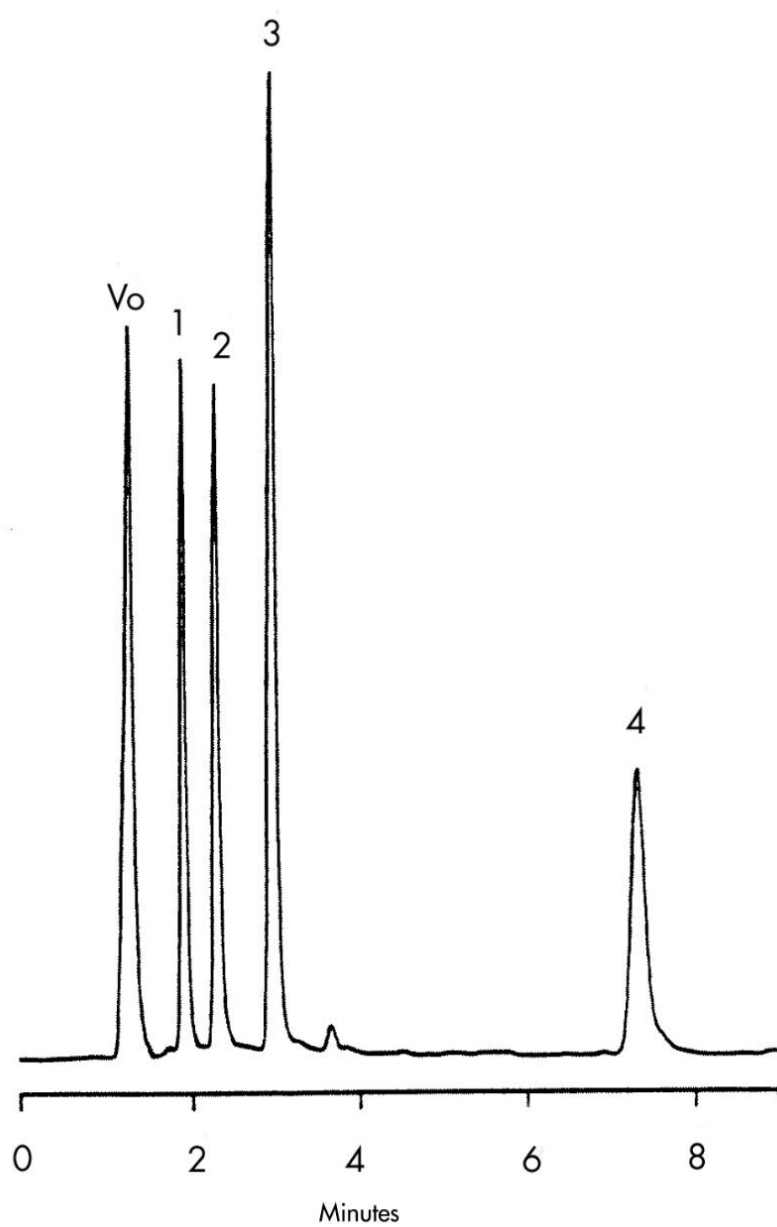

---

## Featured Products

WA31763.149, June 2003

© 2021 Waters Corporation. All Rights Reserved.