


Verapamil

Waters Corporation

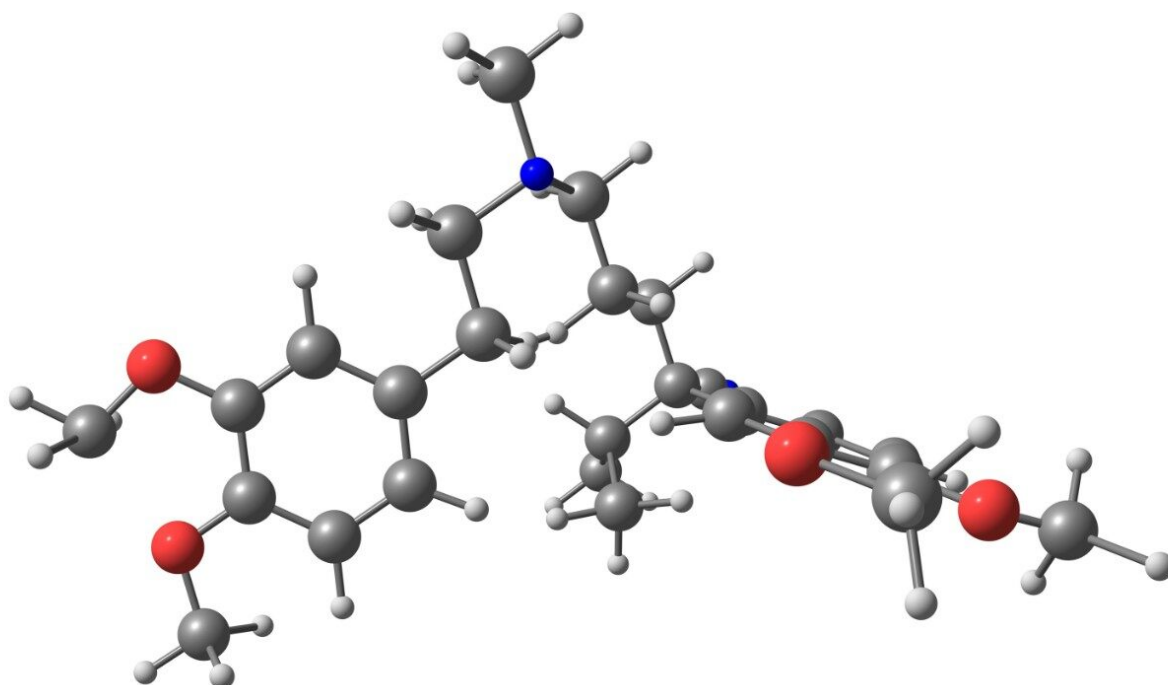


Introduction

The compound analyzed in this study is verapamil.

USP Tailing Factor

1. 1.27

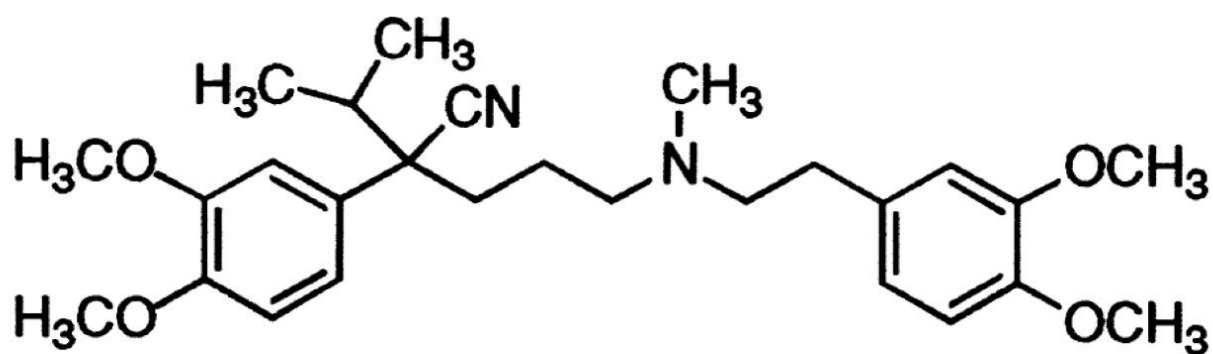

1. Verapamil


Experimental

HPLC Method

Column:	Symmetry C ₈ , 3.9 x 150 mm, 5 µm
Guard Column:	Sentry Guard Column 3.9 x 20mm, 5 µm
Part numbers:	Column - WAT046970, Guard - WAT054250
Mobile phase:	10 mM sodium acetate with 3.3% glacial acetic acid, pH 2.6/acetonitrile 70:30

Flow rate:	1.0 mL/min
Injection volume:	10 μ L of 100 μ g/mL sample
Detection:	UV @ 278 nm

Results and Discussion

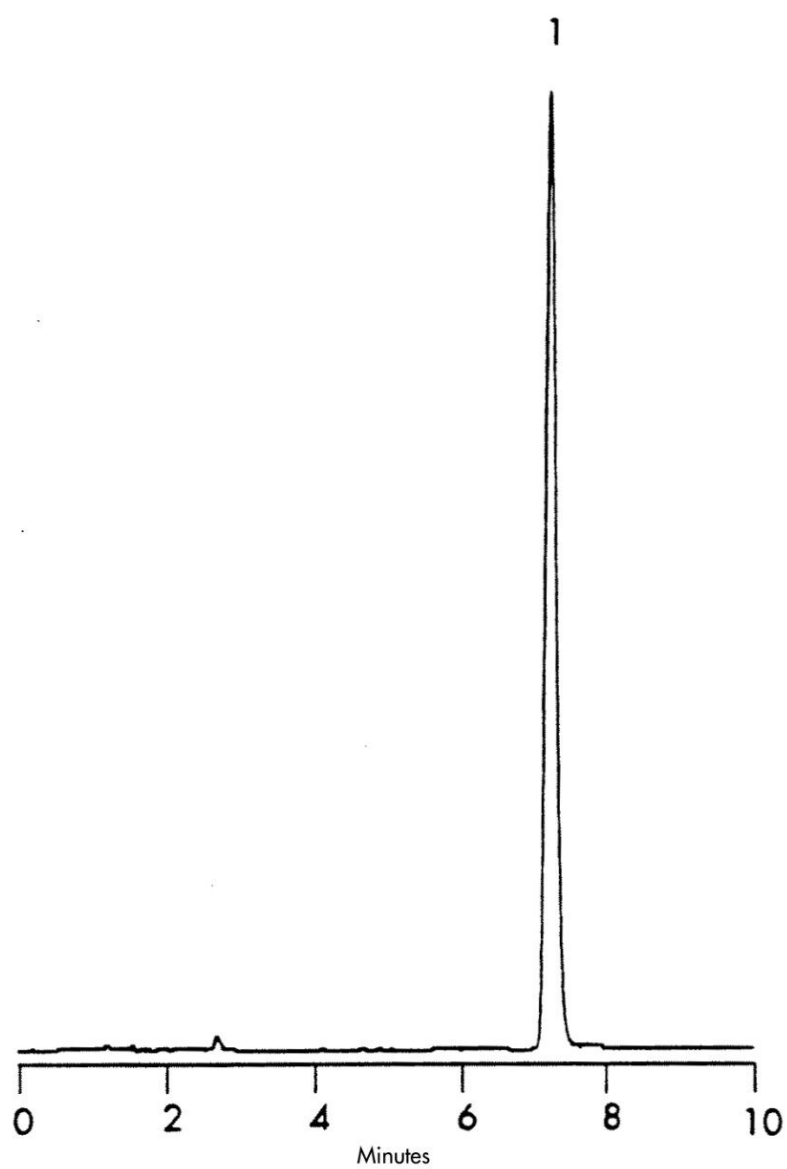

Featured Products

WA31763.177, June 2003

