

응용 자료

Beclomethasone

Waters Corporation



This is an Application Brief and does not contain a detailed Experimental section.

Abstract

This application brief demonstrates analysis of beclomethasone.

Introduction

The compounds used in this study are -

Compounds


USP Tailing Factors

Beclomethasone Dipropionate


1.1

Testosterone Propionate (I.S.)

1.0

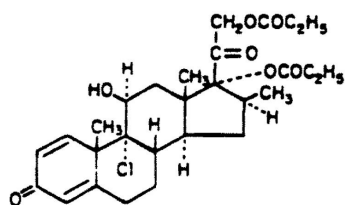

1. BeclomethasoneDipropionate

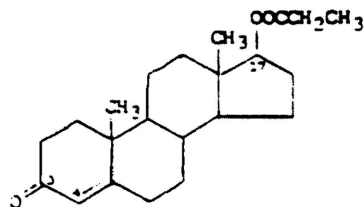

2. Testosterone Propionate (I.S.)

Experimental

HPLC Method

Column:


Symmetry C₁₈, 4.6 x 100 mm, 3.5 μm

Part numbers:

WAT066220

Mobile phase:	Acetonitrile/water 65:35
Flow rate:	1.4 mL/min
Injection volume:	15 μ L of 0.1 mg/mL of each beclomethasone-dipropionate and testosterone propionate (internal standard)
Detection:	UV @ 254 nm

Results and Discussion

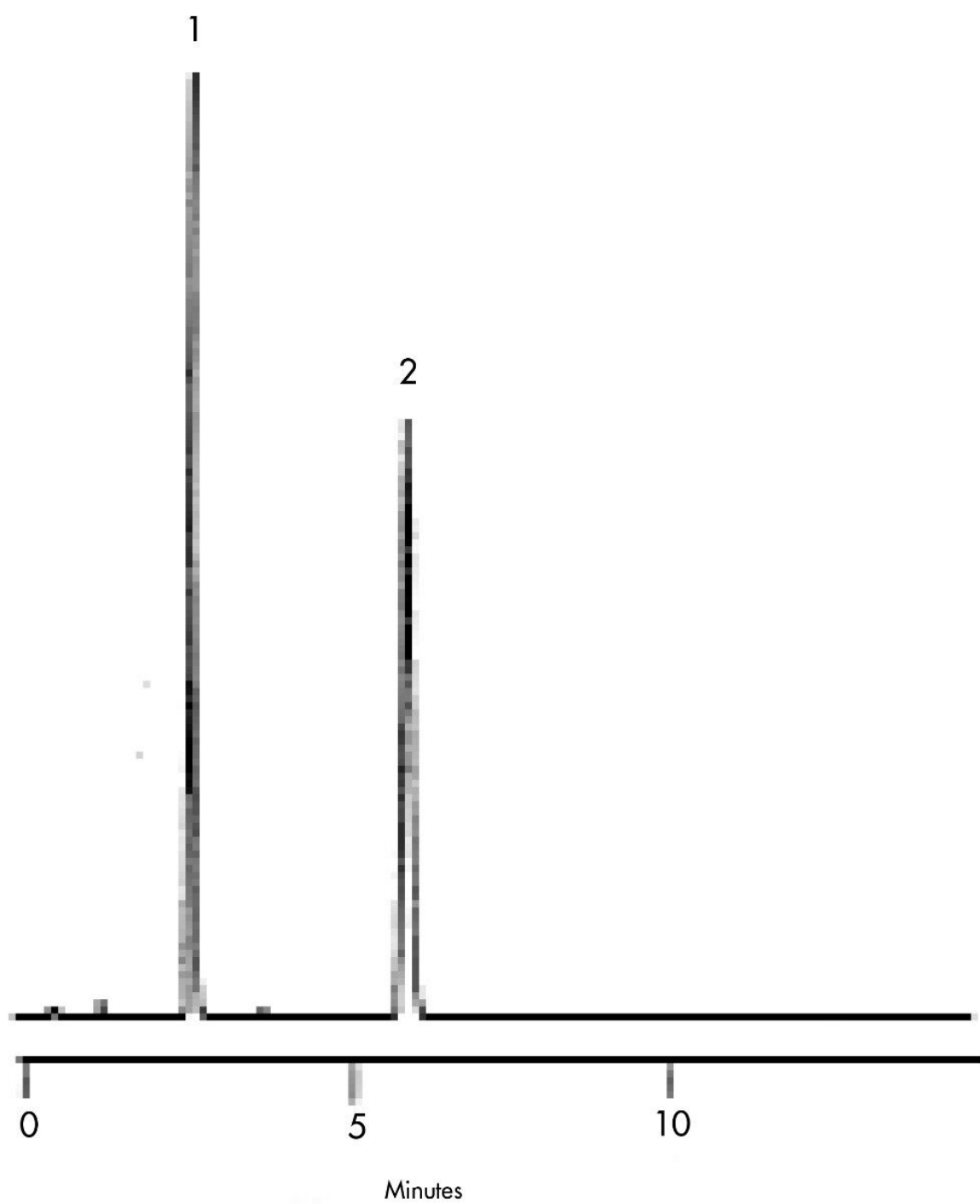

Featured Products

WA31763.35, June 2003

© 2021 Waters Corporation. All Rights Reserved.