



## Busiperone

---

Waters Corporation



This is an Application Brief and does not contain a detailed Experimental section.

---


## Abstract

This application brief demonstrates analysis busiperone.

---

## Introduction

The compound analyzed in this study is busiperone (USP Tailing Factors - 1.08).

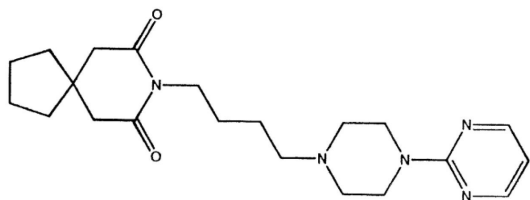

**1. Busiperone**

---


## Experimental

### HPLC Method

|                   |                                                       |
|-------------------|-------------------------------------------------------|
| Column:           | Symmetry Shield RP <sub>18</sub> , 4.6 x 150 mm, 5 µm |
| Part number:      | 186000109                                             |
| Mobile phase:     | 20 mM phosphate buffer, pH 7.0/acetonitrile<br>55:45  |
| Flow rate:        | 1.0 mL/min                                            |
| Injection volume: | 5 µL                                                  |
| Detection:        | UV @ 240 nm                                           |

---

## Results and Discussion

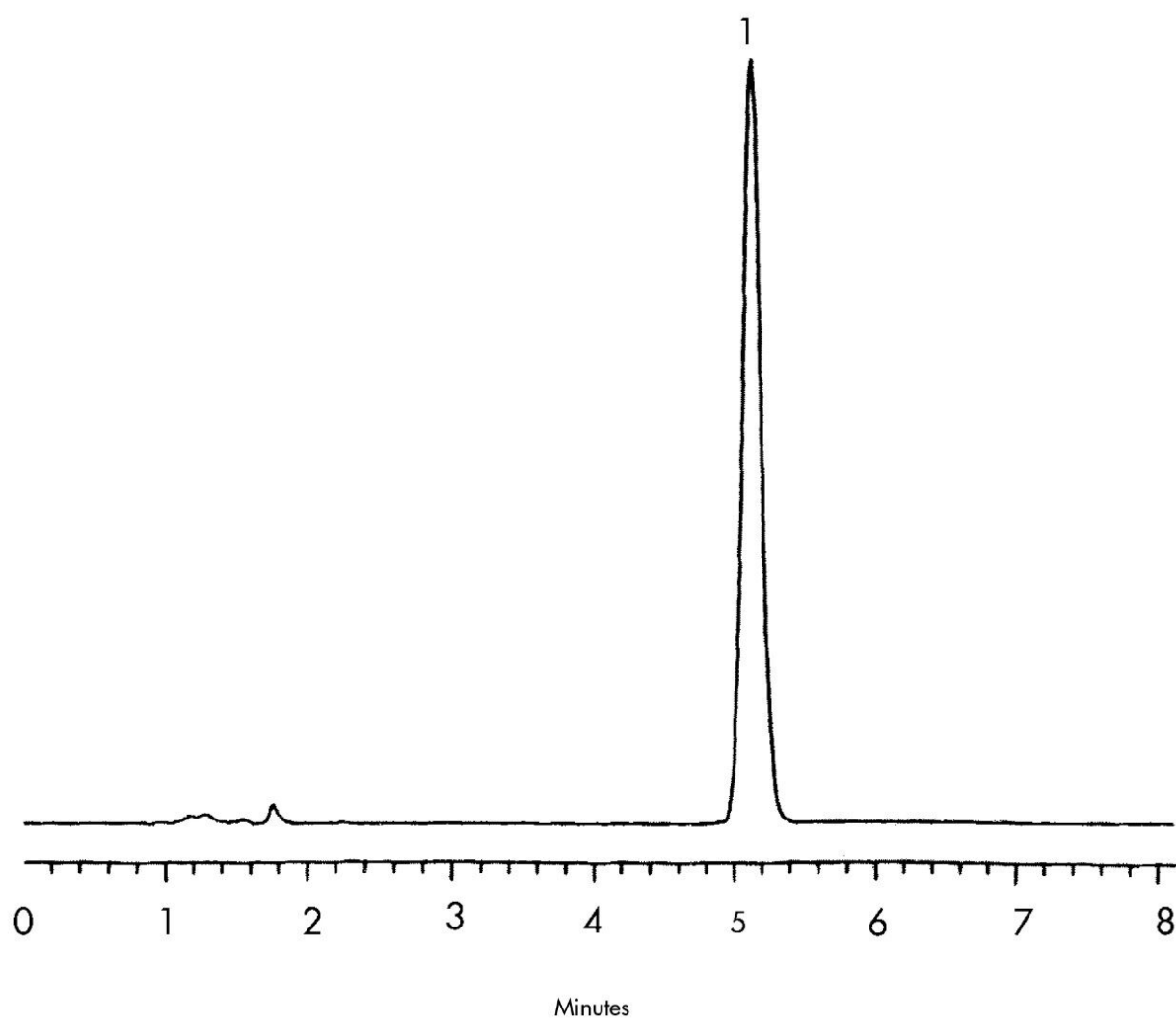

---

## Featured Products

WA31763.44, June 2003

