


Nota applicativa

## Carbadox

---

Waters Corporation

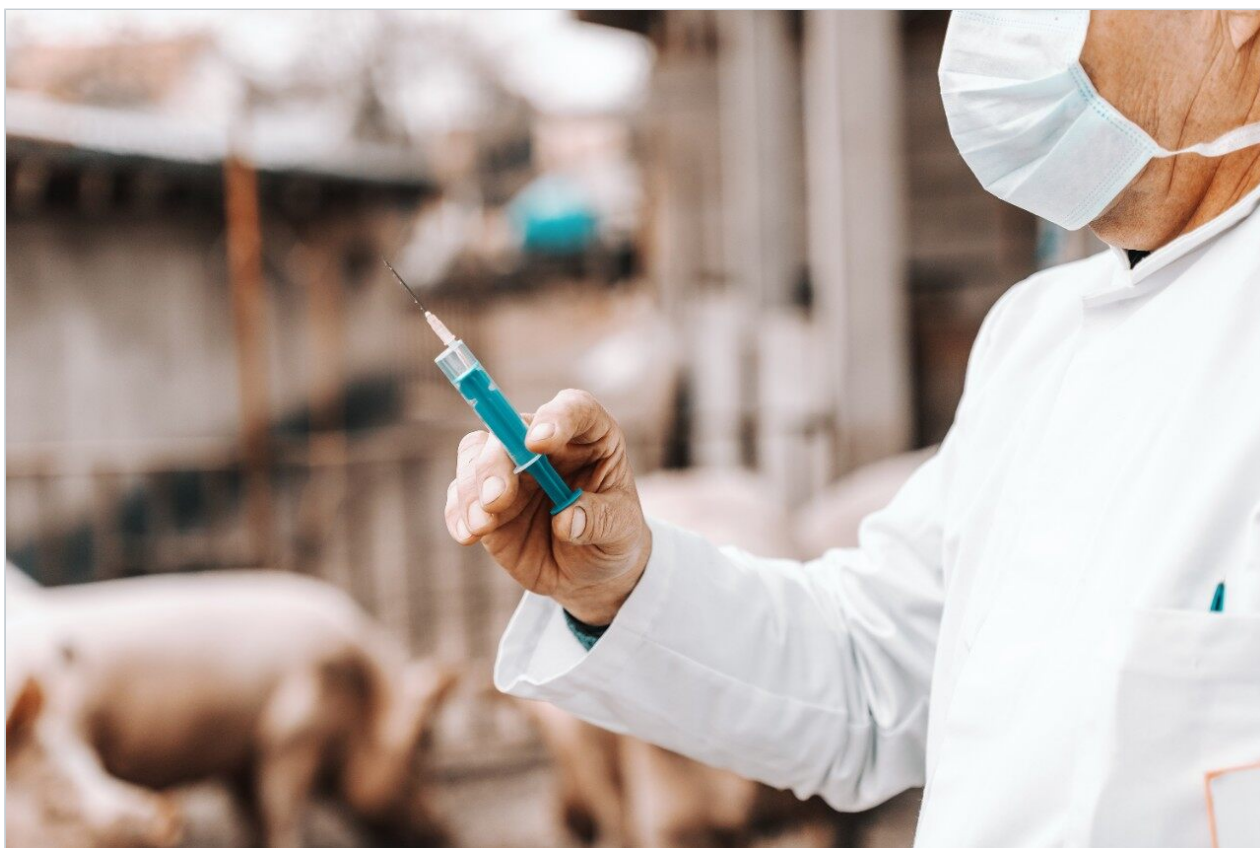

This is an Application Brief and does not contain a detailed Experimental section.

---


### Abstract

This application brief demonstrates analysis of carbadox.

---

## Introduction

The compound analyzed in this study is carbadox (USP Tailing Factor - 1.1).

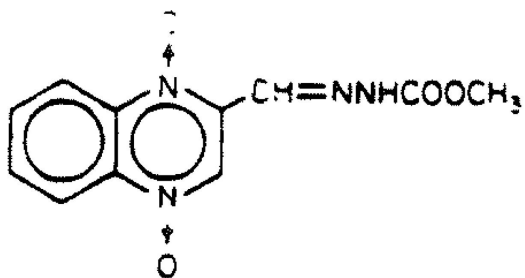

### 1. Carbadox

---

## Experimental

### HPLC Method


Column:	Symmetry C <sub>18</sub> , 4.6 x 100 mm, 3.5 μm
Part number:	WAT066220
Mobile phase:	20 mM ammonium acetate, pH 5.0/Acetonitrile 85:15
Flow rate:	1.4 mL/min
Injection volume:	15 μL of 0.025 mg/mL
Temperature:	30 °C

Detection:

UV @ 254 nm

---

## Results and Discussion

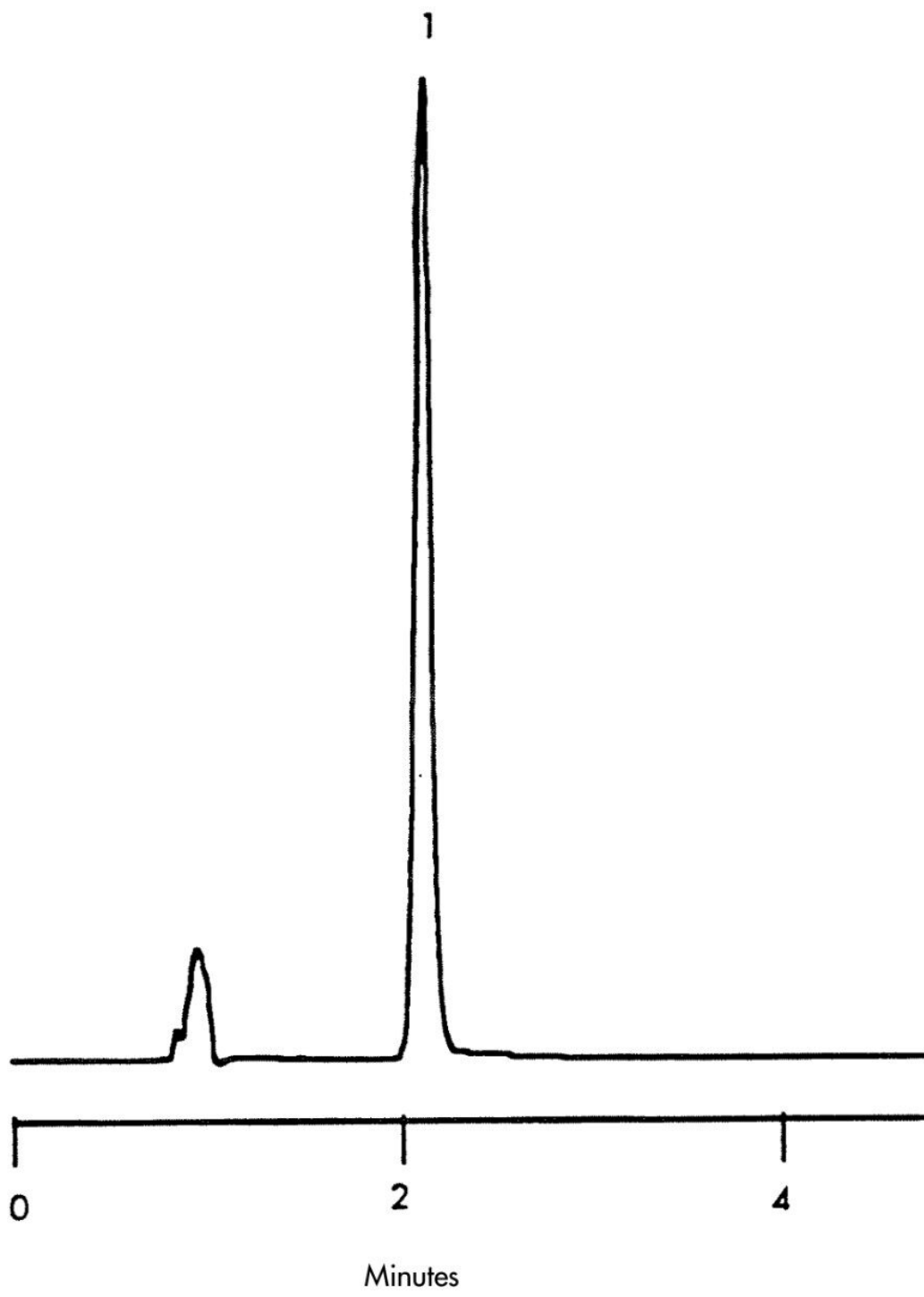

---

Featured Products

WA31763.47, June 2003

© 2021 Waters Corporation. All Rights Reserved.