

Cytochrome c Tryptic Digest - Recovery Study

Waters Corporation

This is an Application Brief and does not contain a detailed Experimental section.

Abstract

This application brief details on cytochrome c tryptic digest recovery study using Symmetry Columns.

Introduction

The compound used in this study is cytochrome c tryptic digest.

Experimental

HPLC Method

Column: Symmetry 300 C₁₈, 4.6 x 50 mm, 5 μm

Part number: WAT106209

Mobile phase A: 0.1% TFA in Water

Mobile phase B: 0.1% TFA in Acetonitrile

Flow rate: 1.0 mL/min

Injection volume: 20 μ L


Temperature: 35 $^{\circ}$ C

Detection: UV @ 214 nm

Gradient

Time (min)	Profile	
	%A	%B
0	100	0
45	70	30

Results and Discussion

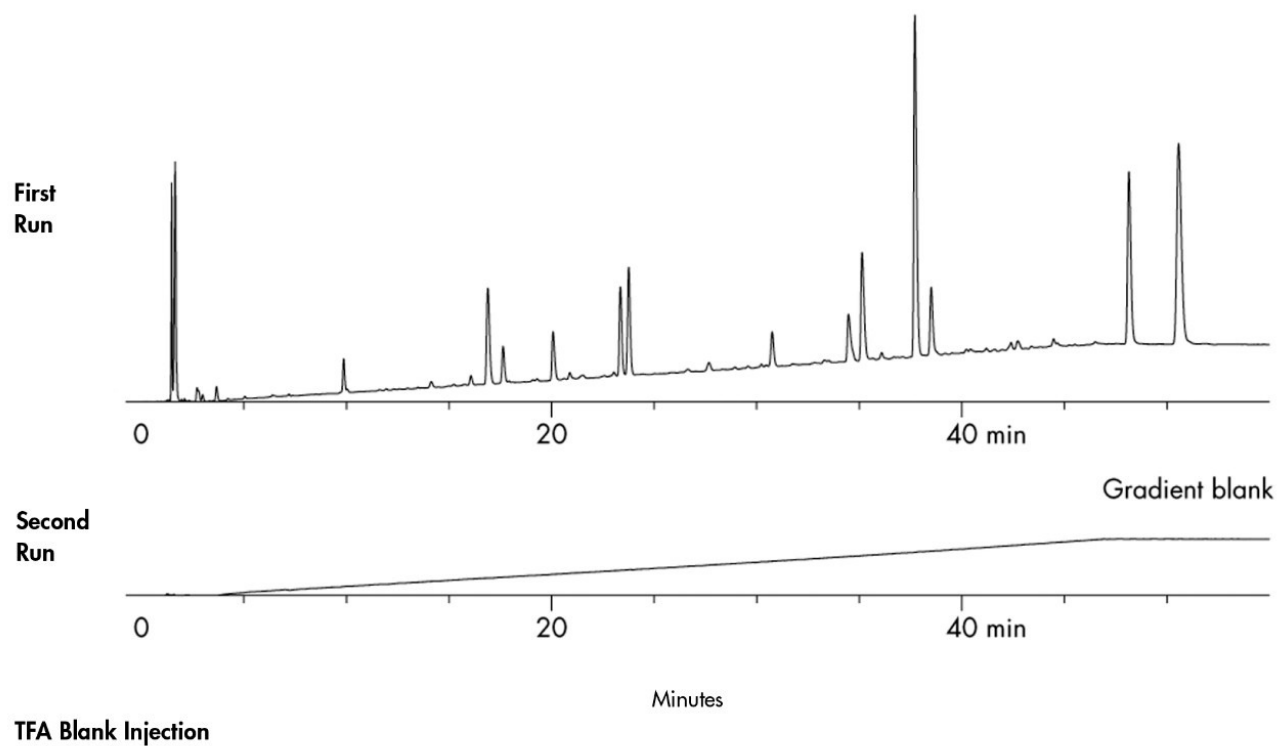

Featured Products

WA31763.68, June 2003

© 2022 Waters Corporation. All Rights Reserved.
