

Nota applicativa

Cytochrome c Tryptic Digest - Reproducibility Study

Waters Corporation

This is an Application Brief and does not contain a detailed Experimental section.

Abstract

This application brief details on cytochrome c tryptic digest reproducibility using Symmetry Columns.

Introduction

The compound used in this study is cytochrome c tryptic digest.

Experimental

HPLC Method

Column: Symmetry 300 C₁₈, 3.9 x 150 mm, 5 µm

Part number: WAT106154

Mobile phase A: 0.05% TFA in Water

Mobile phase B: 0.1% TFA in Acetonitrile

Flow rate: 1.0 mL/min

Injection volume: 20 µL


Temperature: 35 °C

Detection: UV @ 220 nm

Gradient

Time (min)	Profile	
	%A	%B
0	100	0
45	72	28

Results and Discussion

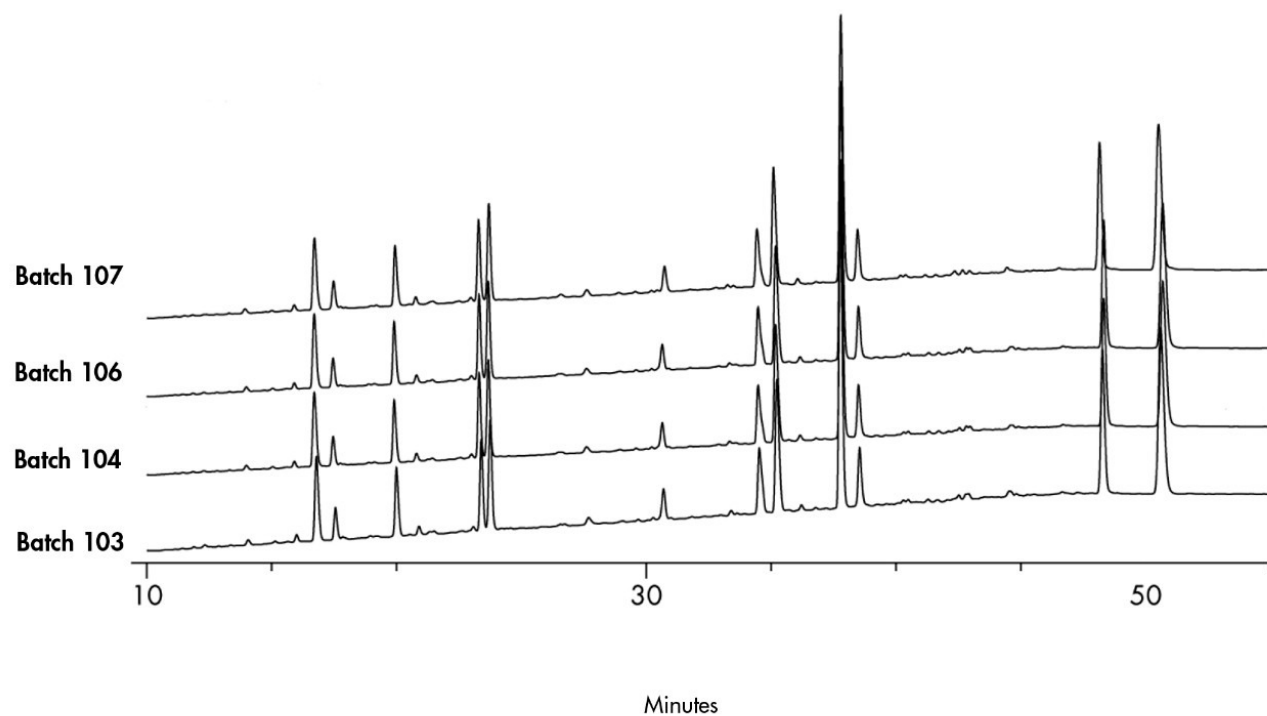

Featured Products

WA31763.69, June 2003

© 2022 Waters Corporation. All Rights Reserved.