

Nota applicativa

## Diphenoxylate HCL

---

Waters Corporation

This is an Application Brief and does not contain a detailed Experimental section.

---

### Abstract


This application brief demonstrates the analysis of diphenoxylate HCL using Symmetry Columns.

---

### Introduction

The compound analyzed in this study is Diphenoxylate HCL.

---

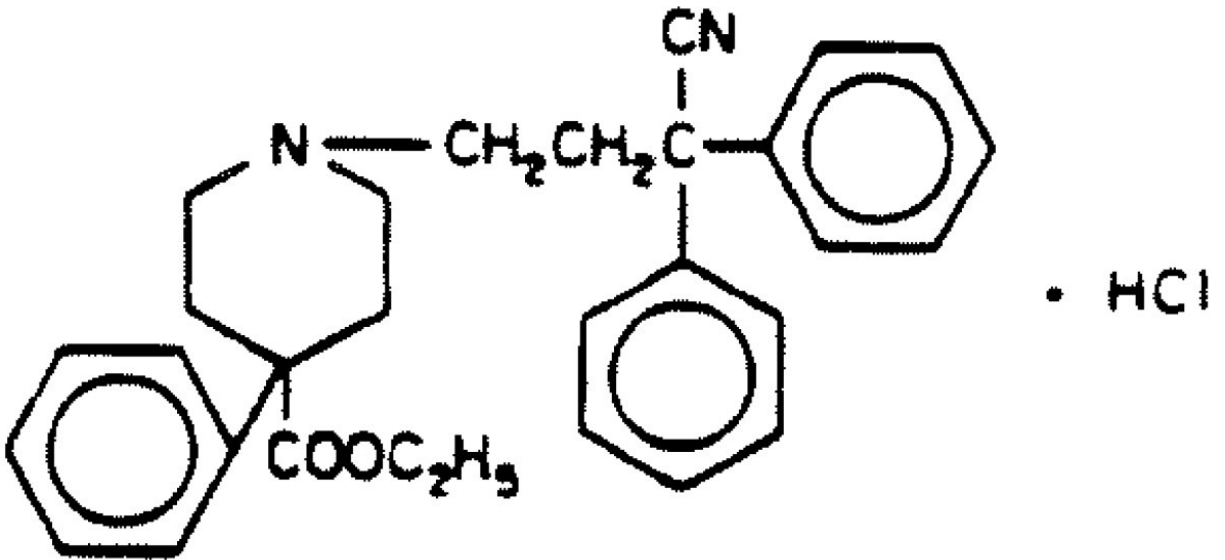

# 1. Diphenoxylate HCL

---

Experimental

HPLC Method

Column:	Symmetry C <sub>18</sub> , 3.9 x 150 mm, 5 μm
Part number:	WAT046980
Mobile phase:	50 mM potassium phosphate, pH 6.0/acetonitrile 30:70

Flow rate: 1.0 mL/min

Injection volume: 10 µL of 100 µg/mL Diphenoxylate hydrochloride


Detection: UV @ 210 nm

## USP Tailing Factors

1. 1.3

---

## Results and Discussion

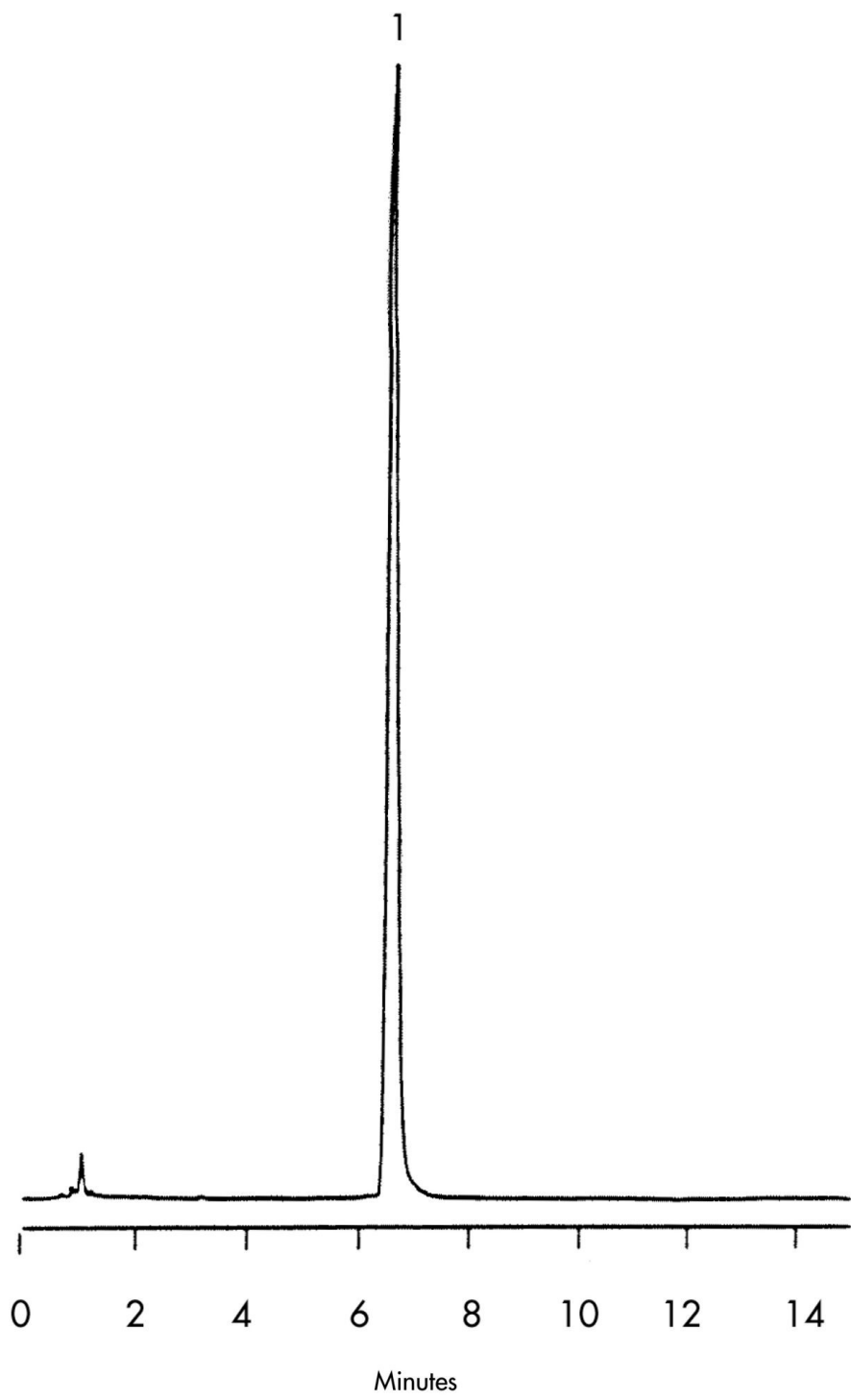

---

## Featured Products

WA31763.77, June 2003

© 2022 Waters Corporation. All Rights Reserved.