

응용 자료

Enalapril

Waters Corporation



This is an Application Brief and does not contain a detailed Experimental section.


Abstract

This application brief demonstrates analysis of enalapril.


Introduction

The compounds used in this study are –

1. Maleic Acid
2. Enalapril

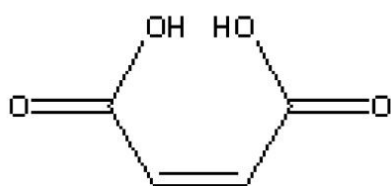

1. Maleic Acid

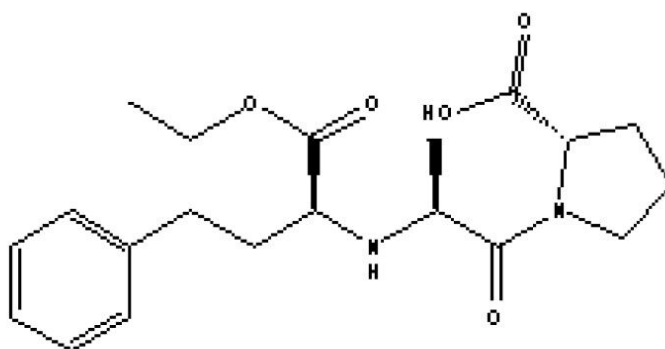

2. Enalapril


Experimental

HPLC Method

Column:	Symmetry C ₈ , 3.9 x 150 mm, 5 µm
Part Number:	WAT046970
Mobile phase:	10 mM potassium phosphate, pH 2/acetonitrile 75:25

Flow rate:	2 mL/min
Injection Volume:	2 μ L of 1.67 mg/mL
Temperature:	80 $^{\circ}$ C
Detection:	UV @ 215 nm

Results and Discussion

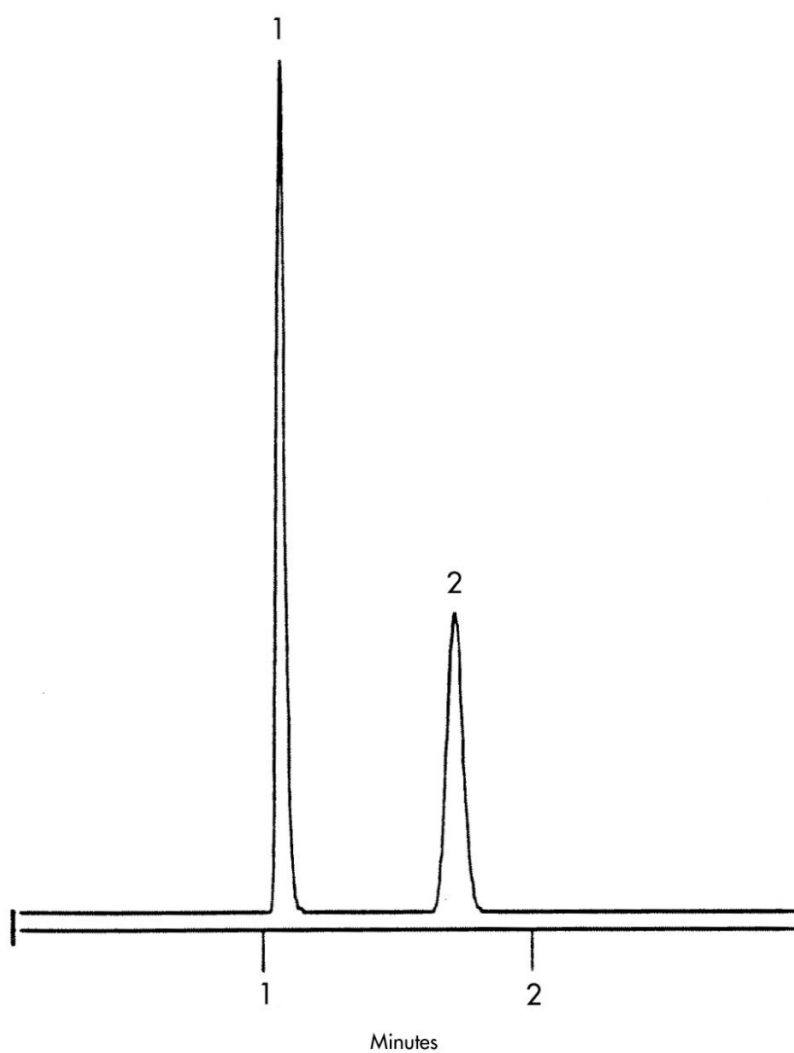

Featured Products

WA31763.80, June 2003

