

Nota de aplicación

Flunitrazepam

Waters Corporation




This is an Application Brief and does not contain a detailed Experimental section.

Abstract

This application brief demonstrates analysis of flunitrazepam.

Introduction

The compound analyzed in this study is flunitrazepam (USP Tailing Factors - 1.2).

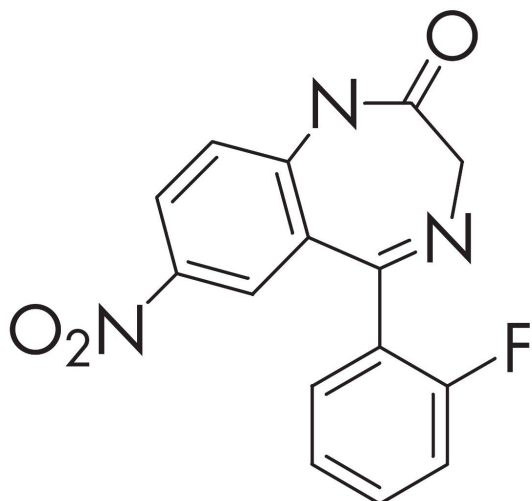

1. Flunitrazepam

Experimental


HPLC Method

| | |
|---------------|--|
| Column: | Symmetry C ₁₈ , 3.9 x 150 mm, 5 µm |
| Part number: | WAT046980 |
| Mobile phase: | 50 mM potassium phosphate pH 6.0/acetonitrile 60:40 |
| Flow rate: | 1.0 mL/min |

Injection volume: 10 μ L of 30 μ g/mL flunitrazepam

Detection: UV @ 219 nm

Results and Discussion

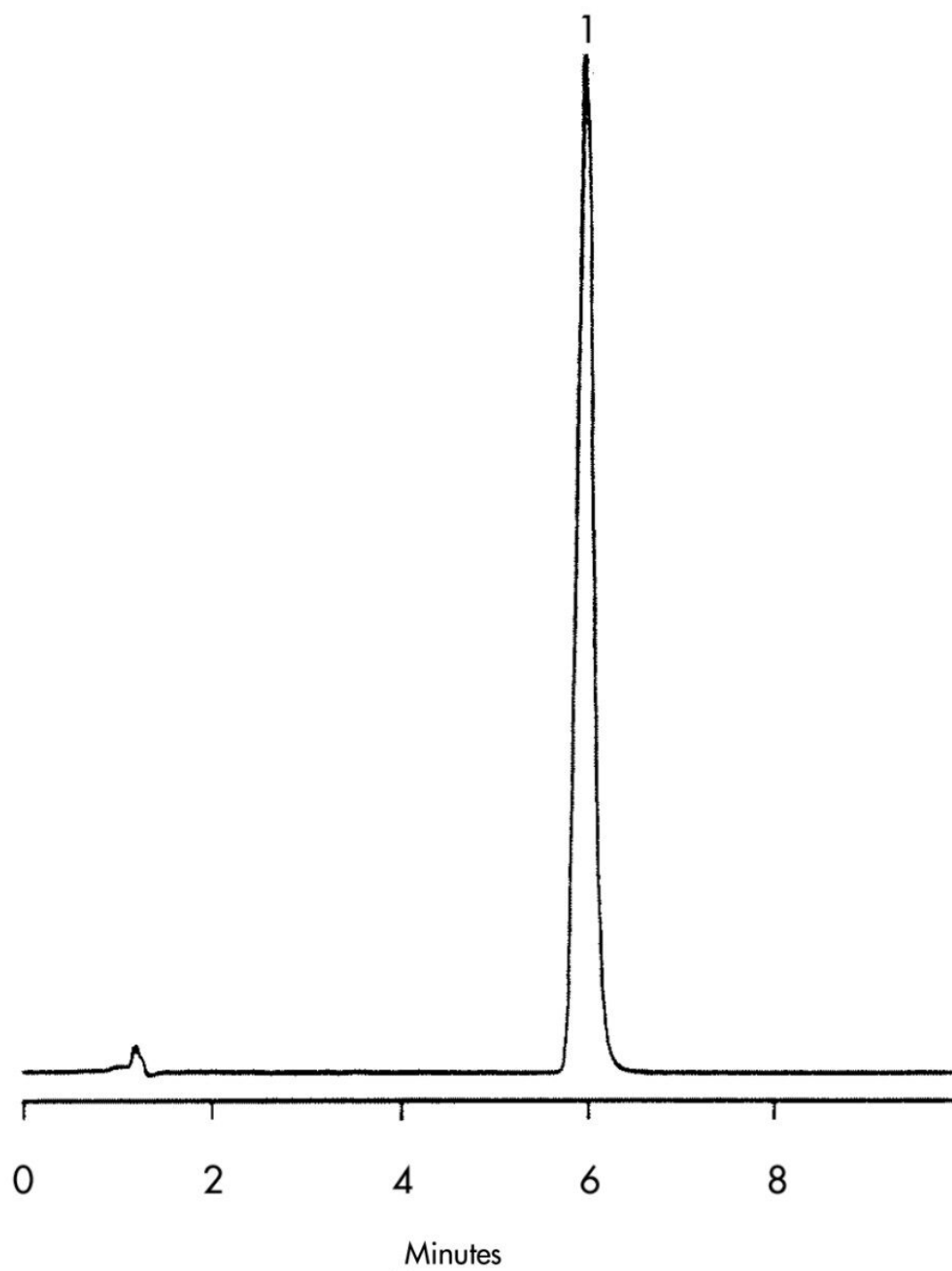

Featured Products

WA31763.85, June 2003

© 2021 Waters Corporation. All Rights Reserved.