

응용 자료

Four Drug Mixture- Warfarin, Dextromethorphan, Triprolidine, and Tetracaine

Waters Corporation



This is an Application Brief and does not contain a detailed Experimental section.

Abstract

This application brief demonstrates analysis of four drug mixture- warfarin, dextromethorphan, triprolidine, and tetracaine.

Introduction

The compounds used in this study are –

Compounds	USP Tailing Factors
Warfarin	1.3
Dextromethorphan HCL	1.6
Tripolidine HCL	1.1
Tetracaine	1.2



HPLC Method

Column:	Symmetry Shield RP ₁₈ , 3.9 x 150 mm, 5 μm
Part number:	186000108
Mobile phase A:	25 mM Potassium phosphate, pH 7
Mobile phase B:	Methanol
Flow rate:	1.0 mL/min


Injection volume: 10 µL of mixture with 0.3 mg/mL triprolidine, 0.4 mg/mL dextromethorphan, 0.15 mg/mL tetracaine, and 0.05 mg/mL of warfarin

Detection: UV @ 280 nm

Gradient

Time (min)	Profile	
	%A	%B
0	85	15
26	20	80

Results and Discussion

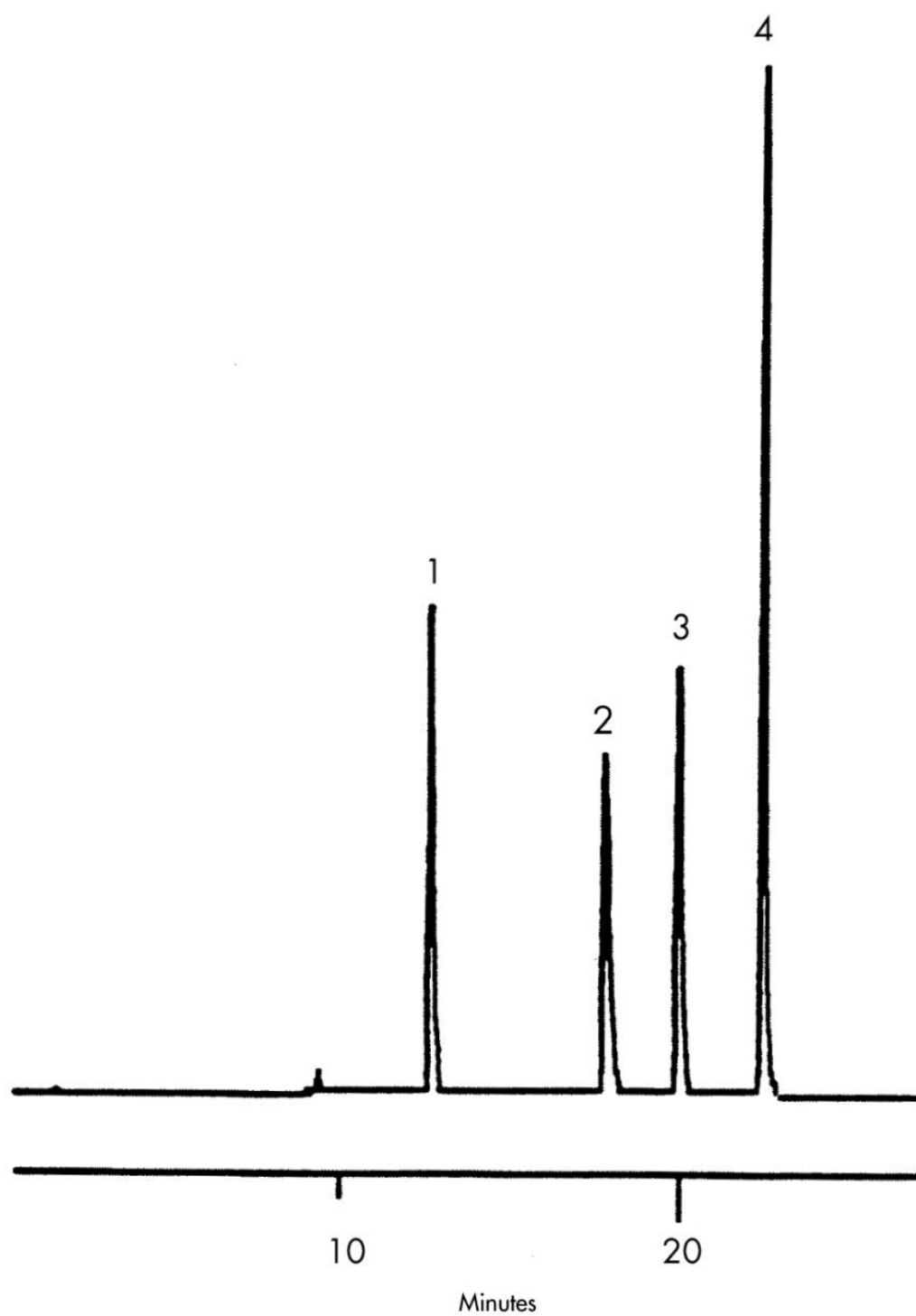

Featured Products

WA31763.87, June 2003

© 2021 Waters Corporation. All Rights Reserved.