

Pharmaceutical Residues in Environmental Samples – LC-MS And LC/PDA Results of LC/PDA Analysis, 2.5 ppb

Waters Corporation



This is an Application Brief and does not contain a detailed Experimental section.


Abstract

This application brief highlights the analysis of pharmaceutical residues in environmental samples using optimized Oasis MCX SPE method.


Introduction

Compounds used in this application brief are:


1. Acetaminophen
2. Phenylpropanolamine
3. Salicylic acid
4. Diphenhydramine
5. Clofibric acid
6. Ethynylestradiol
7. Tamoxifen
8. Ibuprofen

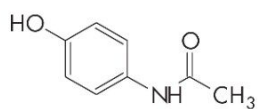


ACETAMINOPHEN

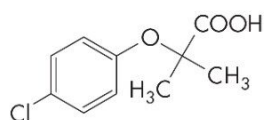


CLOFIBRIC ACID

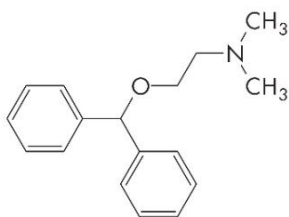


DIPHENHYDRAMINE

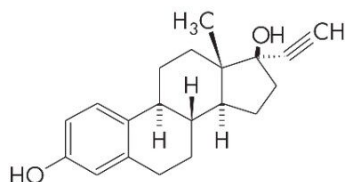


ETHYNYLESTRADIOL

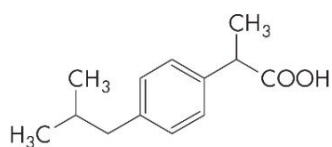


IBUPROFEN

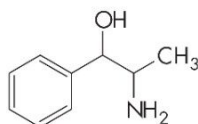

PHENYLPROPANOLAMINE

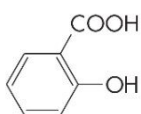

SALICYLIC ACID

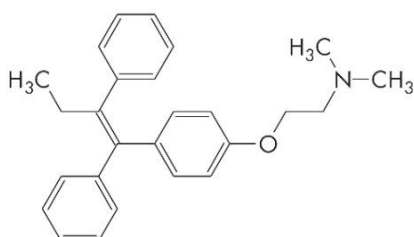

TAMOXIFEN

Experimental

HPLC Method

Column:

XTerra MS C₁₈ 4.6 x 100 mm, 3.5 μm (p/n:
186000436)

Mobile phase A:

15 mM NH₄COOH, pH 4.0

Mobile phase B:	MeOH
Flow rate:	1.0 mL/min
Injection volume:	40 µL
Detection:	MS ESI+
Instrument:	Alliance 2695, Micromass ZQ

Gradient

Time (min)	Profile	
	%A	%B
0.0	75	25
10.0	10	90

LC-MS Conditions Waters Micromass ZMD

Compound	MW	ESI+	ESI-
Acetaminophen	151.2	152(40) 110(40)	150(20)
Phenylpropanolamine	151.2	152(40) 134(40)	
Salicylic acid	138.1		137(17)
Diphenhydramine	255.3	256(25) 167(40)	
Clofibric	214.7		213(17)

Compound	MW	ESI+	ESI-
acid			127(17)
Ethinylestradiol	296.4		295(40) 145(80)
Tamoxifen	371.5	372(40)	
Ibuprofen	206.3		205(17) 159(30)

Optimized SPE Method For LC-MS Determination of Pharmaceutical Residues in Environmental Samples

Conditions for Oasis MCX 6 cc/150 mg (60 µm) Part Number 186000255


Oasis HLB Plus Part Number 186000132

OPTIMIZED SPE METHOD FOR LC/MS DETERMINATION OF PHARMACEUTICAL RESIDUES IN ENVIRONMENTAL SAMPLES


Conditions for

Oasis® MCX 6 cc/150 mg (60 µm) Part Number 186000255

Oasis® HLB Plus Part Number 186000132

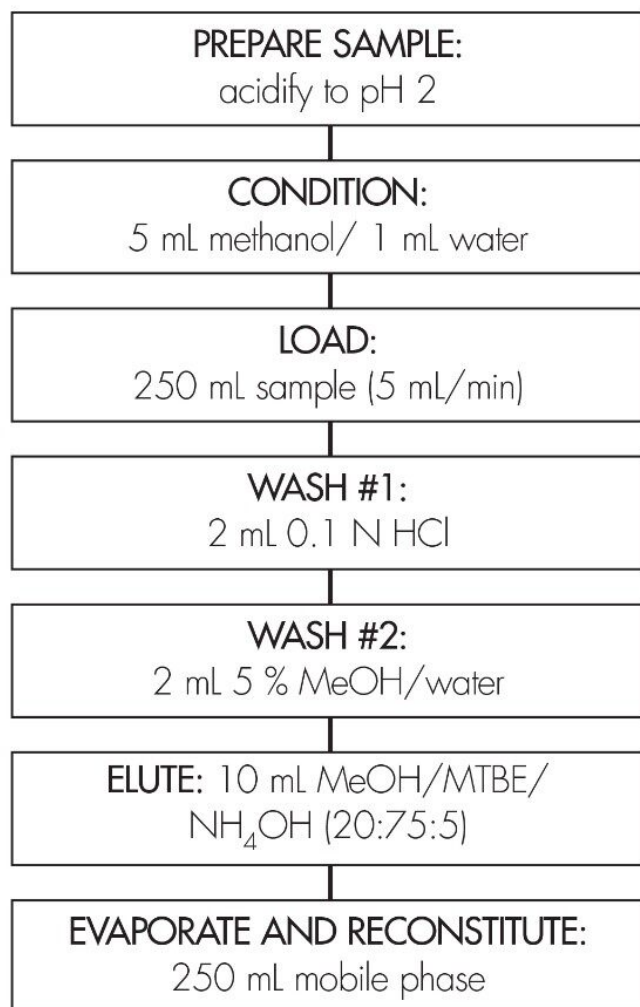

MeOH - methanol MTBE - methyl tbutyl ether
NH₄OH - 30% ammonium hydroxide

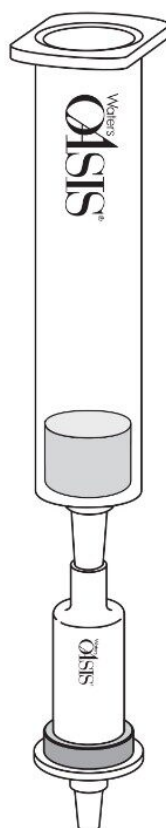

Oasis® MCX 6 cc,
150 mg (60 mm)
Part Number 186000255



Oasis® HLB Plus
Part Number 186000132

to increase retention
of acidic analytes.



Results and Discussion


COMPOUNDS:

1. Acetaminophen
2. Phenylpropanolamine
3. Salicylic acid
4. Diphenhydramine
5. Clofibric acid
6. Ethynylestradiol
7. Tamoxifen
8. Ibuprofen

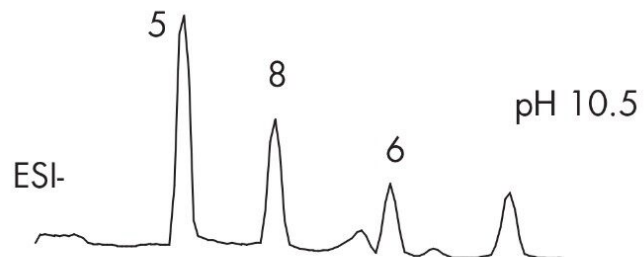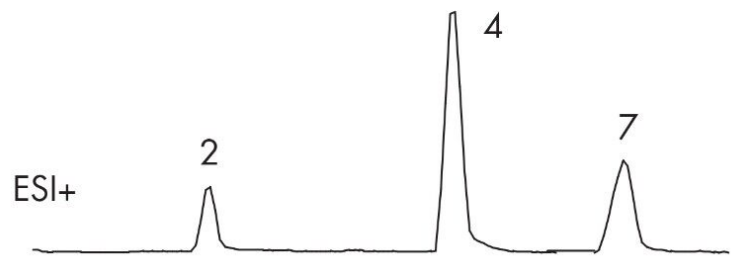

* note: ethynylestradiol
response only @ pH 10.5

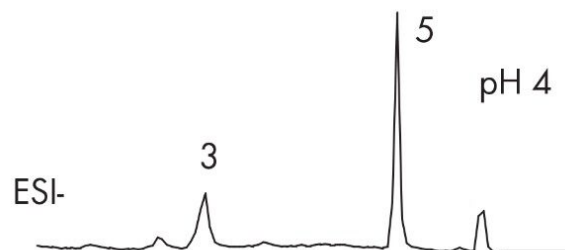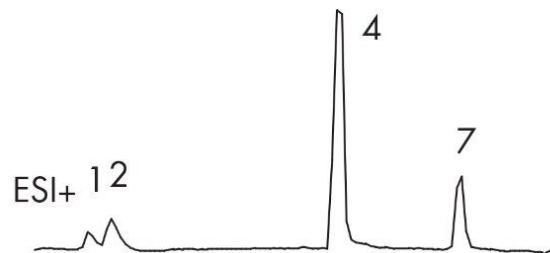

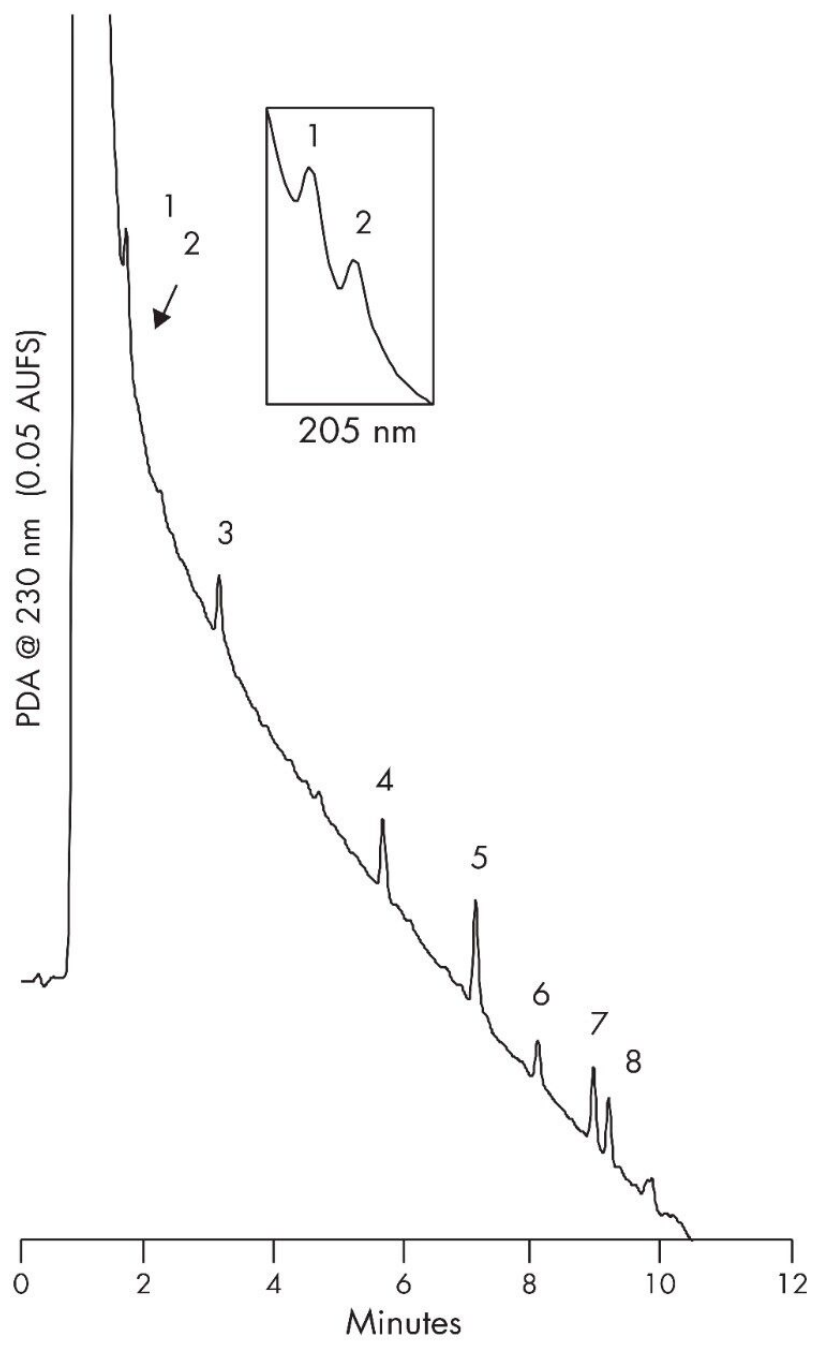

Compound	Response (area/ng/L)	R ²	%RSD/level*
1. acetaminophen	248.6	0.9999	11, 4.3, 5.1, 3.6
2. phenylpropanolamine	635	0.999	621, 22, 7.0, 3.6
3. salicylic acid (3 levels)	27.1	0.9999	nd, 15, 19, 9.0
4. diphenhydramine	3793	0.9985	35, 4.6, 19, 6.8
5. clofibric acid	62.35	0.9998	2.6, 2.3, 3.2, 3.9
6. ethynylestradiol (pH4)	no response		
6. ethynylestradiol (pH11)	104	0.9997	6.2, 5.1, 8.1, 11
7. tamoxifen	1092	0.9990	15, 11, 13, 7.9
8. ibuprofen	19.6	0.9999	20, 5.2, 6.4, 5.8

* %RSD observed for 5 replicates at 25, 100, 500, 2500 ng/L spike levels

Compound	Spike levels	% Recovery	%RSD
1. acetaminophen	500, 2500	int, 78	7.7, 10
2. phenylpropanolamine	500, 2500	nd, 61	nd, 5.7
3. salicylic acid	500, 2500	60, 61	20, 10
4. diphenhydramine	500, 2500	89, 86	3.7, 4.6
5. clofibric acid	500, 2500	101, 89	4.0, 4.2
6. ethynylestradiol	500, 2500	93, 86	7.3, 4.1
7. tamoxifen	500, 2500	76, 76	6.0, 2.1
8. ibuprofen	500, 2500	82, 77	7.7, 4.7

Featured Products

