

Pharmaceutical Compounds High Temperature

Waters Corporation



This is an Application Brief and does not contain a detailed Experimental section.


Abstract

This application brief demonstrates analysis of pharmaceutical compounds at high temperature.


Introduction

The compounds used in this study are –


1. Lidocaine
2. Prednisolone
3. Naproxen
4. Amitriptyline
5. Ibuprofen

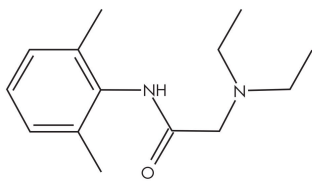


Lidocaine

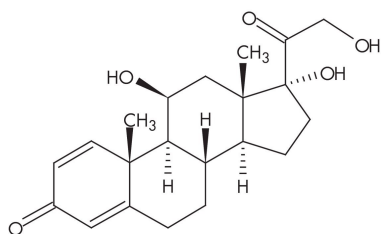


Prednisolone

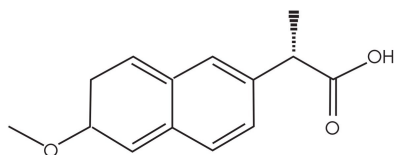

Naproxen

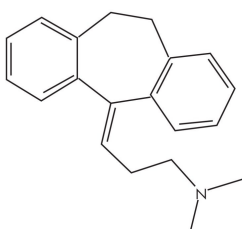

Amitriptyline

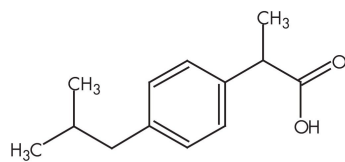

Ibuprofen

Experimental



LC Conditions

Column:	Xterra MS C ₁₈ , 2.1 x 20 mm <i>IS</i> , 2.5 µm, (p/n: 186001921)
Mobile phase A:	10 mM CH ₃ COOHN ₄ , pH 5.0/ACN (97/3)
Mobile phase B:	10 mM CH ₃ COOHN ₄ , pH 5.0/ACN (10/90)
Flow rate:	0.6 mL/min; 1.0 mL/min
Injection volume:	10 µL
Sample concentration:	20 µg/mL
Temperature:	30 °C; 60 °C
Detection:	UV @ 210 nm
Instrument:	Alliance 2795 with 996 PDA

Gradient

Time (min)	Profile	
	%A	%B
0.0	95	5
1.7	5	95

Results and Discussion

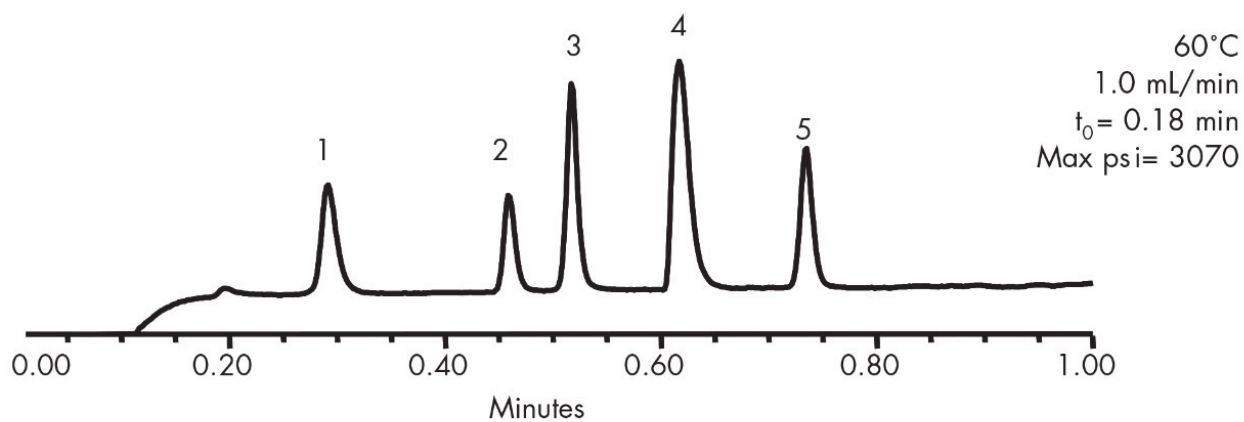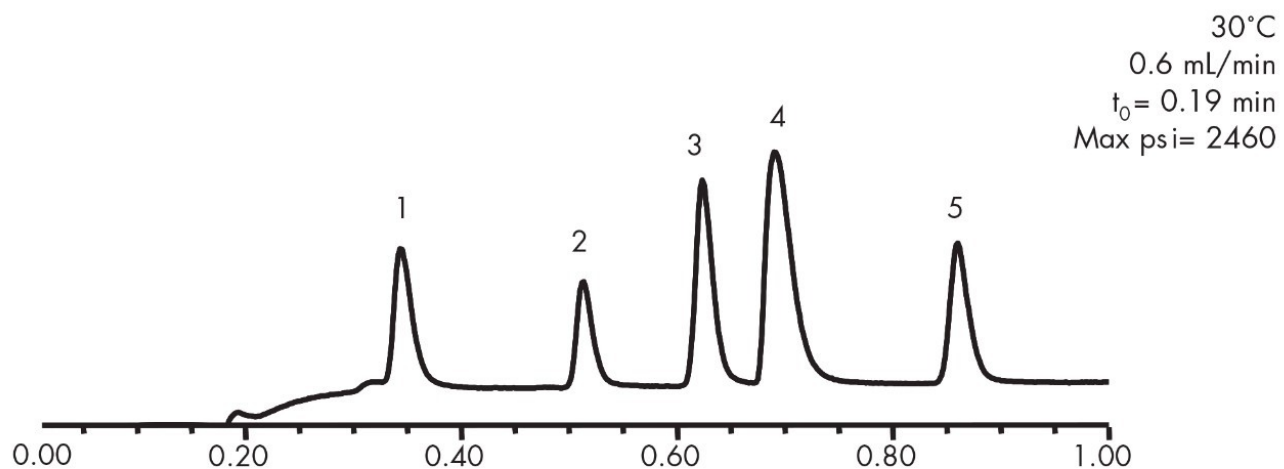

Featured Products

· [Alliance HPLC System <https://www.waters.com/534293>](https://www.waters.com/534293)

WA31787.20, June 2003



© 2021 Waters Corporation. All Rights Reserved.