

Note d'application

## Analysis of Tricyclic Antidepressants Using XBridge C<sub>18</sub>, pH 10.5

---

Waters Corporation



This is an Application Brief and does not contain a detailed Experimental section.

---


### Abstract

This application brief highlights the analysis of tricyclic antidepressants using XBridge C<sub>18</sub> Columns.


---

## Introduction


### Compounds

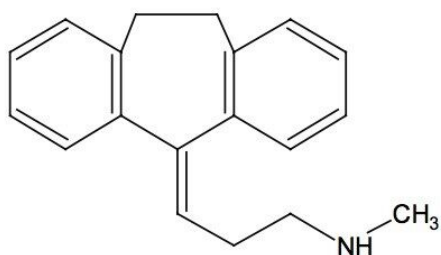


1. Nortriptyline

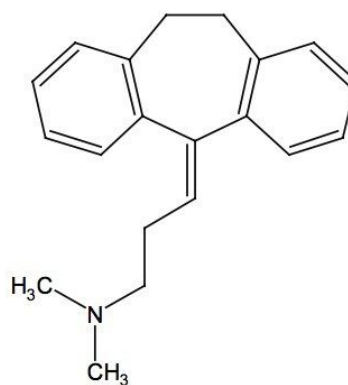

3. Amitriptyline

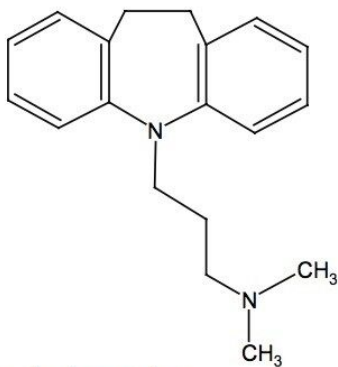

2. Imipramine

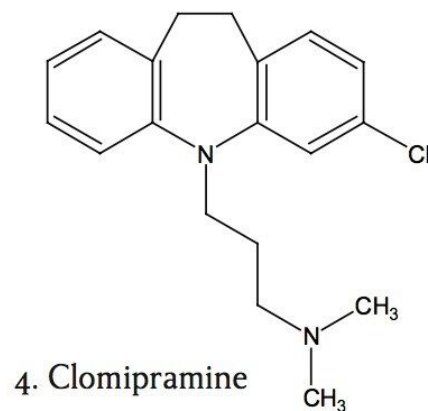

4. Clomipramine

---


## Experimental

### Test Conditions

Column:	XBridge C <sub>18</sub> , 5 µm, 4.6 x 150 mm p/n: 186003116
Isocratic Mobile Phase:	10 mM ammonium bicarbonate buffer, pH 10.5/acetonitrile (25:75)
Flow Rate:	1.5 mL/min
Injection Volume:	10 µL
Sample Concentration:	20 µg/mL in water
Temperature:	25 °C
Detection:	UV @ 254 nm
Instrument:	Alliance 2695 with 2996 PDA

---

## Results and Discussion

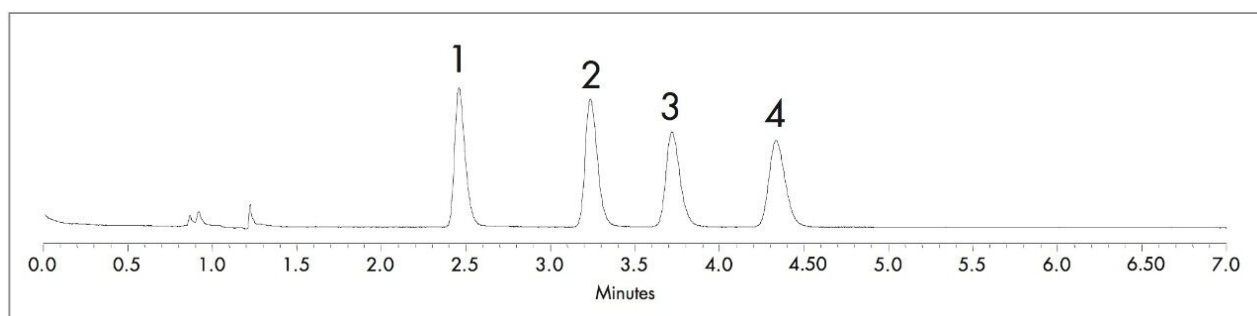


---

## Featured Products

Alliance HPLC <<https://www.waters.com/514248>>

WA60166, April 2007

© 2021 Waters Corporation. All Rights Reserved.