

## UPLC-MS Analysis of Food Sugars Using ACQUITY UPLC BEH Amide Columns with Acetone as Organic Modifier

---

Waters Corporation

This is an Application Brief and does not contain a detailed Experimental section.

---

### Abstract

This application brief highlights the UPLC-MS analysis of food sugars using ACQUITY UPLC BEH Amide Columns with acetone as organic modifier.

---






### Introduction

Compounds used for this study includes:

1. Fructose
  2. Glucose
  3. Sucrose
-

4. Maltose

5. Lactose

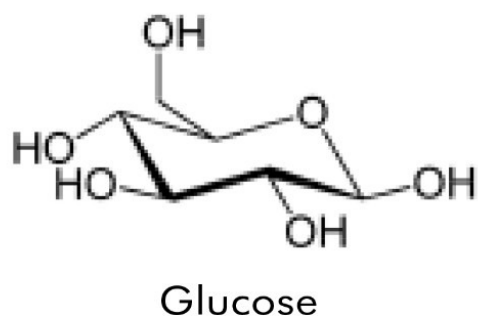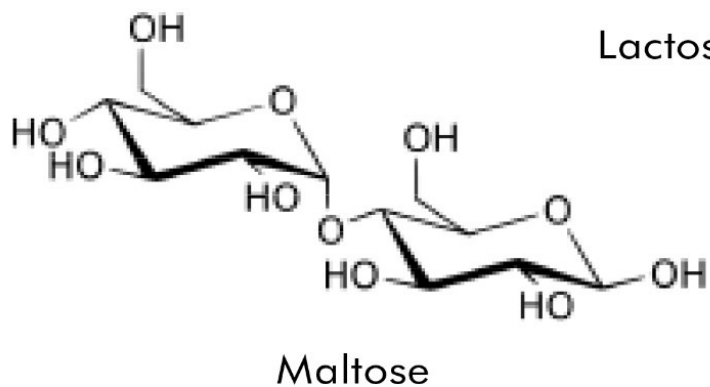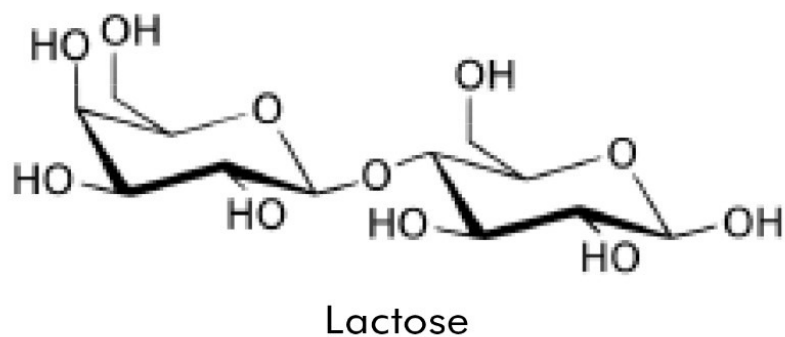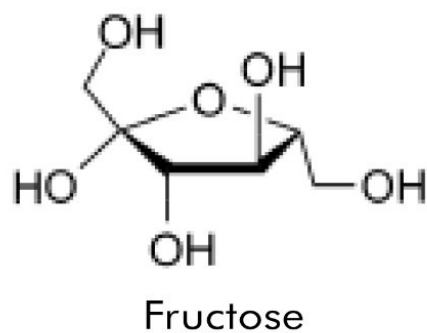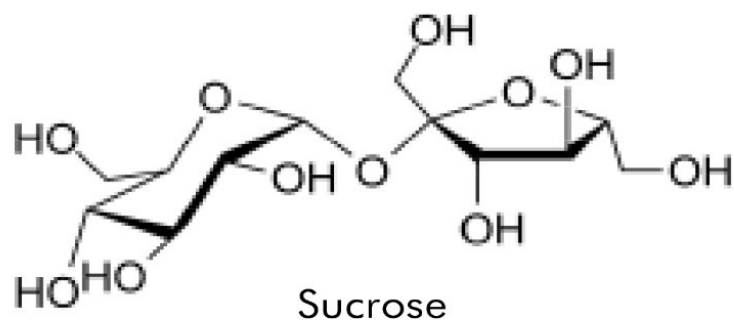

---

Experimental

Chromatographic Conditions

---


Column:	ACQUITY UPLC BEH Amide 2.1 x 50 mm, 1.7 μm
Part number:	186004800
Mobile phase A:	80/20 MeCN/H <sub>2</sub> O with 0.05% ammonium hydroxide [NH <sub>4</sub> OH]
Mobile phase B:	30/70 acetone/H <sub>2</sub> O with 0.05% ammonium hydroxide [NH <sub>4</sub> OH]
Flow rate:	0.13 mL/min
Flow profile:	94% A/6% B (77% acetone with 0.05% NH <sub>4</sub> OH)
Injection volume:	0.7 μL (PLNO)
Sample concentration:	10 μg/mL each
Sample diluent:	50/50 MeCN/H <sub>2</sub> O
Column temperature:	85 °C
Strong needle wash:	20/80 MeCN/H <sub>2</sub> O (800 μL)
Weak needle wash:	75/25 MeCN/H <sub>2</sub> O (500 μL)
Seal wash:	50/50 MeCN/H <sub>2</sub> O
Instrument:	Waters ACQUITY UPLC with ACQUITY TQD

## Mass Spectrometer Conditions

Ionization mode:	ES <sup>-</sup>
Capillary:	2.8 kV
Cone voltage:	25 V
Source temperature:	120 °C
Desolvation temperature:	350 °C
Desolvation gas flow:	500 L/Hr
Cone:	50 L/Hr
SIR ( <i>m/z</i> ):	179.2 (Fructose, Glucose); 341.3 (Sucrose, Maltose, Lactose)
Dwell time:	0.08 s

---

## Results and Discussion

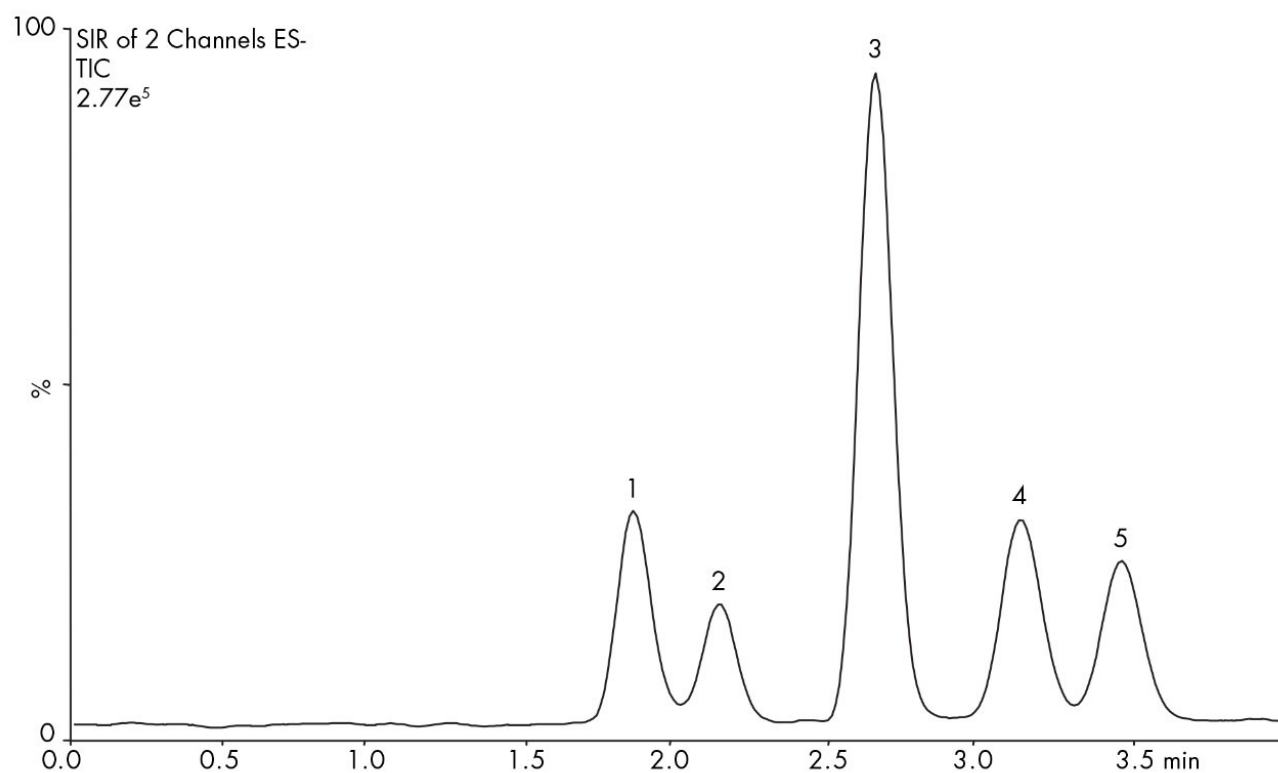

---

## Featured Products

WA60113, October 2009

© 2022 Waters Corporation. All Rights Reserved.

---